

UPDATED INVERSION MATRICES FOR ENGINE EXHAUST PARTICLE SIZER™ (EEPS™) SPECTROMETER MODEL 3090

APPLICATION NOTE EEPS-005 (A4)

Contents

Introduction.....	1
Background.....	2
Fundamental Operation of EEPS Spectrometer	2
Observed EEPS Spectrometer and SMPS Spectrometer Discrepancy.....	3
Cause of EEPS Spectrometer and SMPS Spectrometer Discrepancy	4
Matrix Development.....	6
IM-Compact.....	7
IM-Soot.....	8
Comparison of EEPS Spectrometer Inversion Matrices.....	9
Results Using New Matrices.....	10
Compact Inversion Matrix Results.....	10
Soot Inversion Matrix Results	12
Overall Sizing Comparison.....	14
Software and Firmware Updates	14
Using an Alternative Matrix in EEPS Spectrometer Software Version 3.2.5	15
Which Inversion Matrix Should I use?	15
References	16

Introduction

The TSI Engine Exhaust Particle Sizer™ (EEPS™) spectrometer is widely used to measure the sub-micrometer size distribution of fast changing aerosols. Using a set of 22 electrometers, the EEPS spectrometer records particle size distributions (PSD) at a rate of 10 Hz; allowing for the measurement and analysis of rapidly changing aerosols. While the EEPS spectrometer can record PSDs at a high rate, it is limited to a size resolution of 16 channels per decade, corresponding to 32 channels over the range of 5.6 nm to 560 nm. For slow transient or steady-state aerosols, higher resolution measurements can be made using the TSI Scanning Mobility Particle Sizer™ (SMPS™) spectrometer. While the SMPS spectrometer is broadly considered the standard in measuring sub-micrometer aerosols; measurements with an SMPS spectrometer take at least 5 seconds to complete for a size distribution of a high concentration aerosol (lower concentrations require longer scan times for reasonable results).



It has been reported that discrepancies exist between PSDs measured by EEPS spectrometer and those measured by SMPS spectrometer [1,2,3,4,5]. Above approximately 75 nm, particles measured by EEPS spectrometer are shifted to smaller sizes, resulting in a narrower distribution (having a lower geometric standard deviation (GSD)), and the concentrations of these larger particles are lower than those measured using SMPS spectrometer. It has also been shown that these discrepancies are more severe for aerosols composed of agglomerates (such as soot from engine exhaust) than compact (nearly spherical) aerosols.

To improve the agreement between EEPS spectrometer and SMPS spectrometer measurements, especially at larger particle sizes, TSI has developed two new matrices; one for compact aerosols (IM-compact) and one for soot-like agglomerate aerosols (IM-soot).

This application note details the observed inconsistencies between EEPS spectrometer and SMPS spectrometer when using the existing default matrix, the underlying cause for this discrepancy, the development and validation of the two new matrices, and a description of how one can use the new matrices with the most recent release of EEPS spectrometer software. The matrix development and validation was performed in conjunction with Dr. Xiaoliang Wang of Desert Research Institute (DRI) and the University of Minnesota Center for Diesel Research (CDR) and is discussed in greater detail in the companion papers [6,7].

Background

Fundamental Operation of EEPS Spectrometer

Figure 1 shows a simplified schematic of the EEPS spectrometer measurement region [8]. Sample flow of 10 L/min is drawn into the EEPS spectrometer and charged by a negative unipolar charger to bring the charge state of the aerosol to a uniform polarity. The aerosol is then charged by a positive unipolar charger and enters an annular region between an inner rod, and an outer column consisting of 22 electrometers. The inner rod is separated into three sections. Each section of the inner rod is independently charged to a high positive voltage. The aerosol is transported down the annular region by 39.4 L/min of clean, recirculating sheath air. As the aerosol is transported down the column, individual particles, now positively charged, are repelled from the center column and impact the electrometers that make up the outer column, transferring their charge. The specific electrometer on which a particle impacts depends on the particle's electrical mobility i.e. diameter and state of charge.

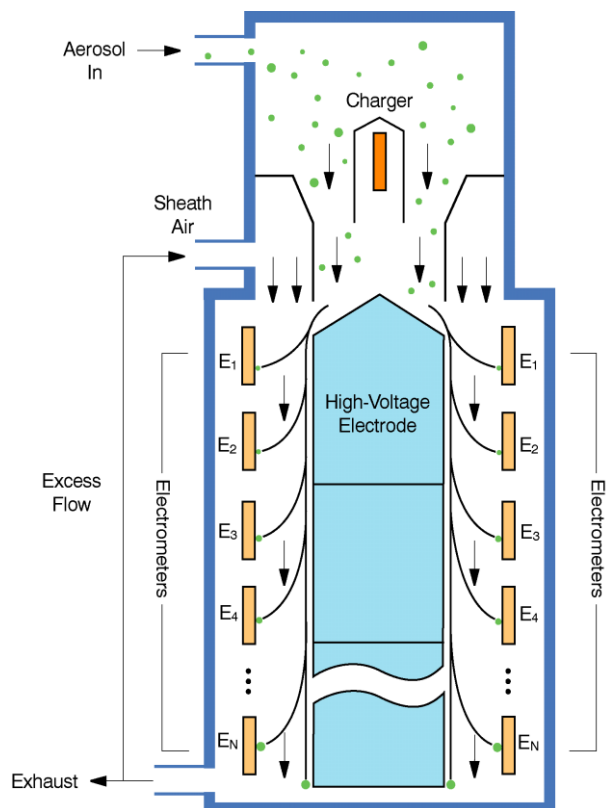


Figure 1: Model 3090 EEPS Spectrometer Schematic Diagram [8]

The electrical currents measured by each of the 22 electrometers are corrected for the image current induced on upper electrometers as particles pass them, destined for lower electrometers. Additionally, the electrometer signals are time aligned to correct for the delay between particles impacting the upper electrometers and larger particles traveling down the column and impacting the lower electrometers. This vector of time aligned and image current corrected signals is referred to as the instrument record, I_j , where j is the electrometer number ($j=1,2,\dots,22$).

The EEPS spectrometer determines the particle size distribution by solving for f , the particle size distribution in Equation 1 [6].

$$I_j = H_{j,i}f_i + u_j$$

Equation 1

Where I_j is the instrument record, $H_{j,i}$ is the EEPS spectrometer instrument matrix, f_i is the particle size distribution and u_j is a vector of noise on the electrometers. $i=1,2,\dots,17$ and corresponds to the boundaries of each of the primary size channels (17 channels equally spaced on a log basis over the 5.6 to 560nm size range of the EEPS spectrometer). The EEPS spectrometer instrument matrix (IM) is what the EEPS spectrometer uses to convert the electrometer signals to particle size distributions. The default EEPS spectrometer matrix, IM-default, was developed in 2004 by Dr. Aadu Mirme at the University of Tartu and is based upon a combination of theoretical modeling of and experimental characterization with prototype EEPS spectrometer units [9].

The same instrument geometry and data inversion algorithm is used with the TSI Fast Mobility Particle Sizer™ (FMPS™); however the FMPS spectrometer is limited to a sample frequency of 1Hz and lacks the analog outputs of the EEPS spectrometer.

Observed EEPS Spectrometer and SMPS Spectrometer Discrepancy

Discrepancies between EEPS spectrometer and SMPS spectrometer PSDs have been reported by several sources. It has been reported that for small particles ($< \sim 75$ nm) the size and concentration of particles as measured by EEPS spectrometer is consistent with SMPS spectrometer measurements; however, above 75 nm, the EEPS spectrometer begins to undersize particles, resulting in narrower particle distributions. It has also been reported that the discrepancy in sizing is more severe when sampling agglomerate particles such as engine soot.

Jeong and Evans [2], and Zimmerman et al. [10,11], examined the discrepancies between EEPS spectrometer and SMPS spectrometer and showed that above ~ 75 nm EEPS spectrometer measured particle diameters were roughly 60% of those measured by SMPS spectrometer, and that above ~ 75 nm, PSDs measured by EEPS spectrometer were narrower than with SMPS spectrometer. The discrepancy was found to be greater for agglomerates such as engine exhaust, than ambient aerosols. An empirical post-processing method was developed by Jeong and Evans, and expanded upon by Zimmerman et al. [10,11] that showed significantly improved agreement between EEPS spectrometer and SMPS spectrometer data.

Similar findings were reported by Asbach et al. [1], who used diesel soot and NaCl to evaluate four SMPS systems in comparison to an FMPS spectrometer. The concentrations measured between the TSI SMPS spectrometer and FMPS spectrometer were in reasonable agreement; however, the FMPS spectrometer resulted in broader distributions with mode diameters that were 9 to 22% below those measured by the SMPS systems. Leskinen et al. [5] showed that EEPS spectrometer and SMPS spectrometer measurements of compact ammonium sulfate particles were consistent in concentration and size for smaller particles but undersized larger particles. For TiO₂ agglomerates, the EEPS spectrometer measured particle size was much smaller than SMPS spectrometer, and that this discrepancy increased at larger particle sizes. Similarly, Kaminski et al. [3] determined that for small compact aerosols (< 40 nm NaCl), sizing and concentration between FMPS spectrometer and SMPS spectrometer agreed well but for larger particles and agglomerates, FMPS spectrometer undersized particle diameter by up to 40%.

Lee et al. [4] compared SMPS spectrometer and FMPS spectrometer measurements of ambient aerosols, as well as lab generated ammonium sulfate, ammonium nitrate, and ammonium chloride aerosols and determined that the mode diameters as measured by FMPS spectrometer were roughly 40% smaller than those measured by SMPS spectrometer. Furthermore, applying a linear correction to the data significantly improved the sizing agreement.

Awasthi et al. [12] examined the difference between EEPS spectrometer and SMPS spectrometer diameters in agglomerates compared to compact aerosols by measuring the PSD of silver nanoparticles before sintering (agglomerates) and after sintering (spherical particles). The undersizing for agglomerates was much more pronounced than the difference when measuring compact particles.

Cause of EEPS Spectrometer and SMPS Spectrometer Discrepancy

A particle's electrical mobility depends on both the size of the particle and the number of charges on the particle. Therefore, to determine particle size using electrical mobility, the charge of the particle must be known. SMPS spectrometer and EEPS spectrometer both impart repeatable, steady, charge states on the aerosol prior to mobility selection. SMPS spectrometer uses a bipolar neutralizer, whereas EEPS spectrometer uses a unipolar corona charger. Charging in a bipolar charger is dominated by diffusion, and the steady state charge distribution was determined by Fuchs [13] and approximated by Wiedensohler [14]. Figure 2 shows the bipolar charge distribution for a soft x-ray neutralizer, covering the size range of the EEPS spectrometer.

In this size range, bipolar charging is roughly independent of particle morphology [15], and the majority of particles exiting a bipolar charger have zero or one charge. As particle size increases the fraction of particles carrying multiple charges also increases.

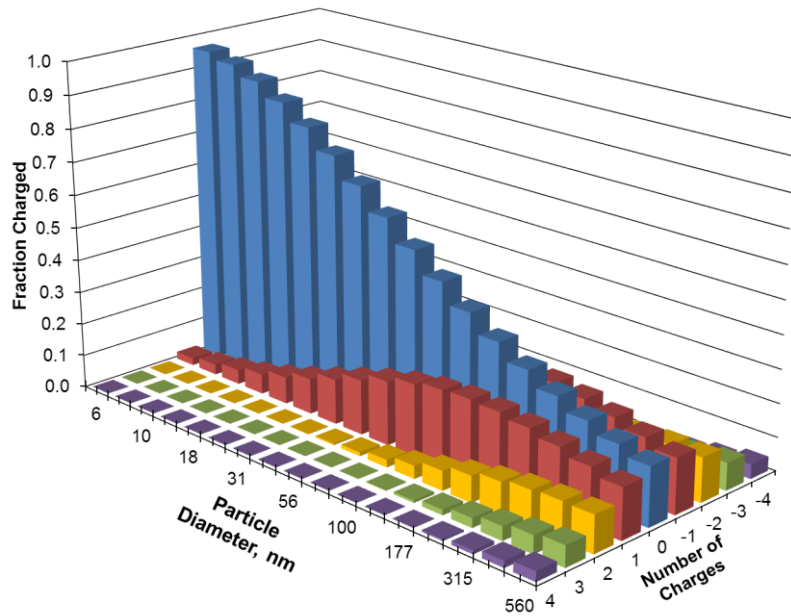


Figure 2: Bipolar charge distribution over EEPS spectrometer size range calculated using Wiedensohler approximation [14].

When using a unipolar charger, the fraction of particles that are charged is much higher, and the increase in frequency of multiple charges with particle size is much more rapid. Additionally, particle shape affects unipolar charging due to the influence of electrostatic forces and increased capacitance [12]. Oh et al. [16] showed that unipolar diffusion charging theory for spheres underestimates the charge level of agglomerates by 30%. Awasthi [12] also showed in Figure 3 that the mean charge per particle for agglomerates is 20 to 50% higher than for spheres, depending on particle size.

Therefore, for particles of the same size as determined by SMPS spectrometer, on average, agglomerates carry more charges than spherical particles. As a result, when using a unipolar charger, agglomerates will appear smaller than spheres. This effect of particle morphology on charging, and therefore particle sizing in an EEPs spectrometer and SMPS spectrometer is shown by Awasthi et al [12] in Figure 4.

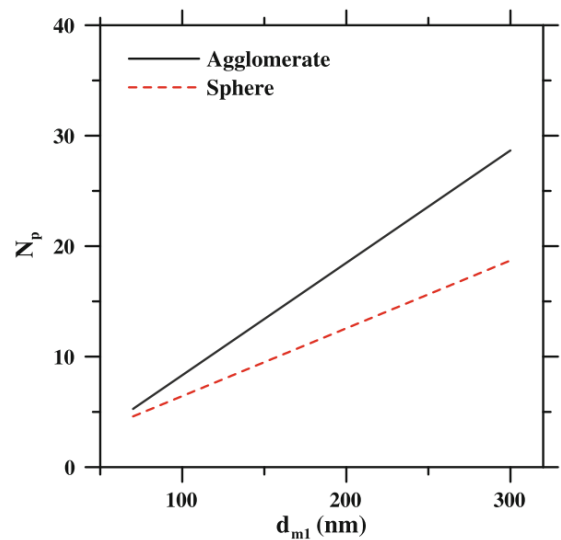


Figure 3: Mean charge per particle for agglomerates and spherical particles as a function of mobility diameter [12].

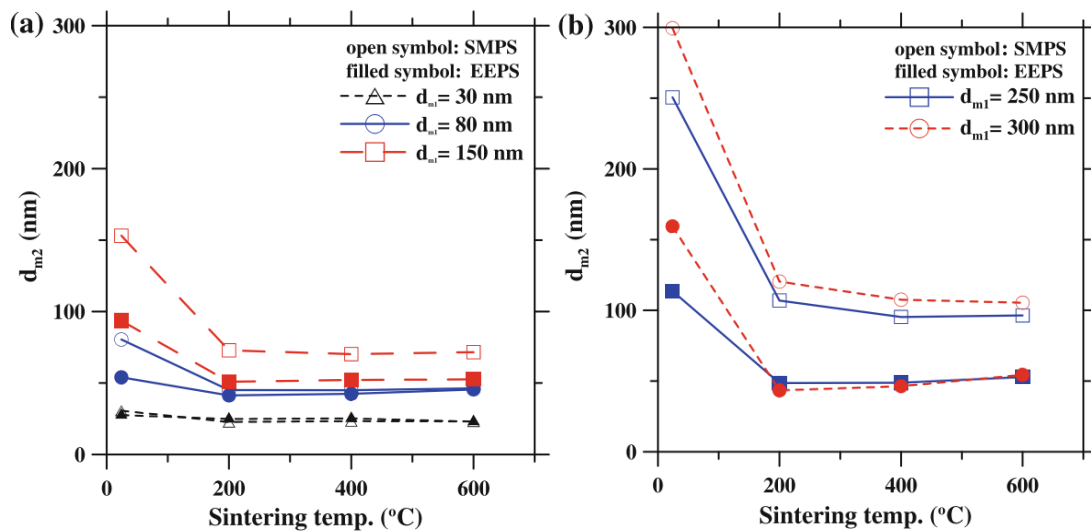


Figure 4: Effect of sintering temperature on particle size as measured by EEPS spectrometer and SMPS spectrometer. As the Sintering temperature increases, the particulates become more spherical (and therefore less agglomerate), and the discrepancy in particle size decreases. [12].

TiO₂ particles were measured with EEPS spectrometer and SMPS spectrometer after a sintering furnace at different temperatures. As the temperature of the sintering furnace increases, and the initially agglomerate particle becomes more spherical, and the discrepancy between the size measured by EEPS spectrometer and by SMPS spectrometer decreases.

Therefore, to accurately size agglomerates and spheres in a device with a unipolar charger, morphology-specific calculations must be used. To address this issue and improve agreement between EEPS spectrometer and SMPS spectrometer, two new inversion matrices were developed.

Matrix Development

The matrix development and validation was performed in conjunction with the help of Dr. Xiaoliang Wang of Desert Research Institute and is discussed in greater detail in the companion papers [6,7]. Whereas the original EEPS spectrometer inversion matrix IM-default was developed using a combination of theoretical calculations and experimental measurements; the IM-compact and IM-soot matrices were developed by empirically calibrating the EEPS spectrometer to SMPS spectrometer measurements of compact, nearly spherical aerosols, and agglomerate aerosols, respectively.

The experimental setups for the development of both matrices were similar; a polydisperse aerosol was generated and then classified using one or two differential mobility analyzers (DMA) to generate a monodisperse aerosol. The monodisperse aerosol was then delivered to an SMPS spectrometer, EEPS spectrometer, and condensation particle counter (CPC) as shown in Figure 5. The EEPS spectrometer, SMPS spectrometer, and CPC sampled simultaneously and the size distributions measured by the SMPS spectrometer, the total concentration measured by SMPS spectrometer and CPC, and the instrument record measured by the EEPS spectrometer were used to generate the EEPS spectrometer inversion matrix, $H_{j,i}$ using Equation 2 [6].

$$H_{j,i} = I_j / N \times \text{dlog}Dp_i$$

Equation 2

Where $H_{j,i}$ is the instrument matrix, I_j is the instrument record measured by the EEPS spectrometer, N is the total number concentration measured by SMPS spectrometer and CPC, and Dp_i is the mode particle diameter of the monodisperse aerosol as measured by the SMPS spectrometer.

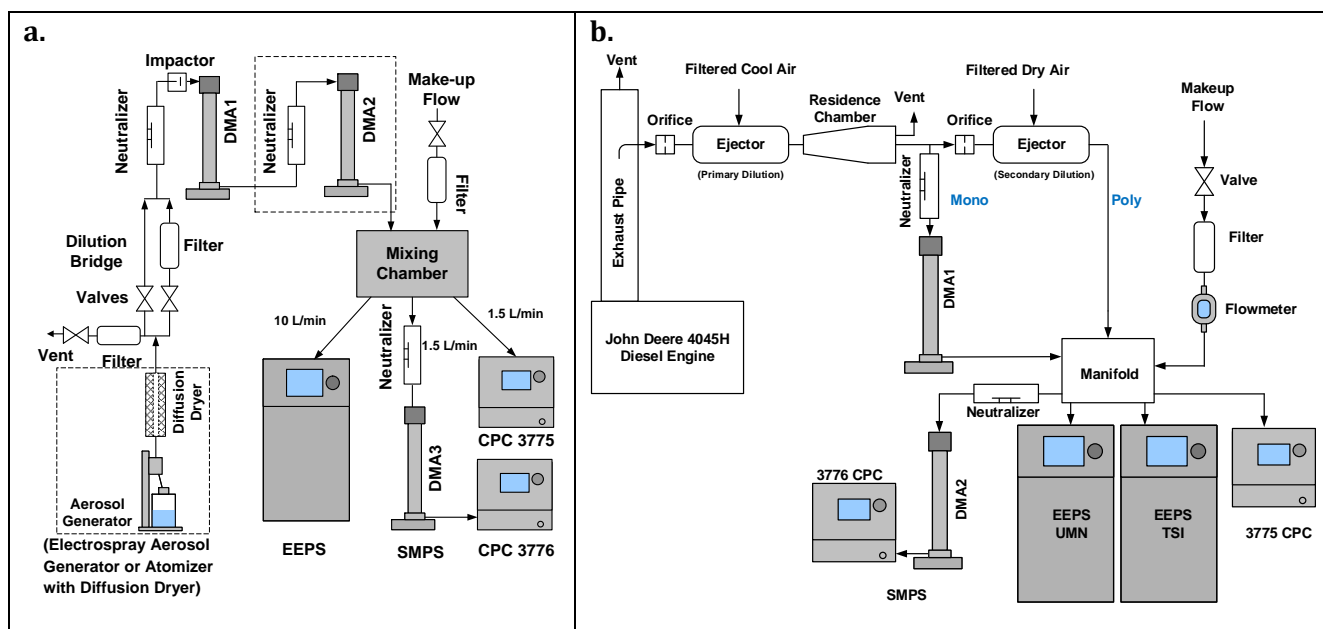


Figure 5: a. Equipment configuration for the development of IM-compact [6] b. Equipment configuration for the development of IM-soot [7]. The sample line lengths from the mixing chamber/manifold to the instruments were proportional to the sample flow rates of each instrument to ensure consistent diffusion losses across all instruments.

IM-Compact

Monodisperse compact aerosols were generated at each of the EEPS spectrometer 17 primary channel diameters (i) by classifying a polydisperse aerosols using a TSI model 3082 electrostatic classifier. The polydisperse aerosols were generated using two aerosol generators.

For particles smaller than 30 nm an electrospray aerosol generator (EAG, TSI model 3480) was used to generate solid sucrose and liquid poly-alpha-olefin (PAO) oil particles. For particles between 30 and 560 nm, a Collison-type atomizer (TSI model 3076) was used to generate solid NaCl and liquid PAO particles.

To reduce the number of multiply-charged particles being delivered to the instruments, the classified sizes were always chosen to be larger than the mode of the polydisperse distribution. Additionally, when the initial aerosol concentration was high enough, two classifiers were used in series to limit the multiply-charged fraction. For the largest sizes, a single classifier and impactors were used to remove larger, multiply-charged particles.

The resulting monodisperse distributions were measured using EEPS spectrometer, a 3775 CPC, and a 3938 SMPS spectrometer consisting of a model 3082 classifier and a 3776 CPC. For sizes below 40 nm model 3085 nano-DMAs were used; for sizes above 40 nm model 3081 long DMAs were used.

For the sizes in which the multiply-charged particles persisted, the instrument matrix column was corrected for the contribution of doubly-charged particles by determining the fraction of doubly charged particles in the sample using SMPS spectrometer data, and determining the instrument kernel function of the doubly charged particle size by interpolating the kernel functions of the two neighboring columns. The fraction of the multiply charged particles was then multiplied by the interpolated kernel function to determine the contribution to the uncorrected instrument matrix column. The corrected matrix column was then determined by subtracting this contribution from uncorrected instrument matrix column. This procedure is examined in detail in [6].

The resulting instrument matrix is plotted in Figure 6 where each curve corresponds to a column of the instrument matrix, H_{ji} from Equation 2.

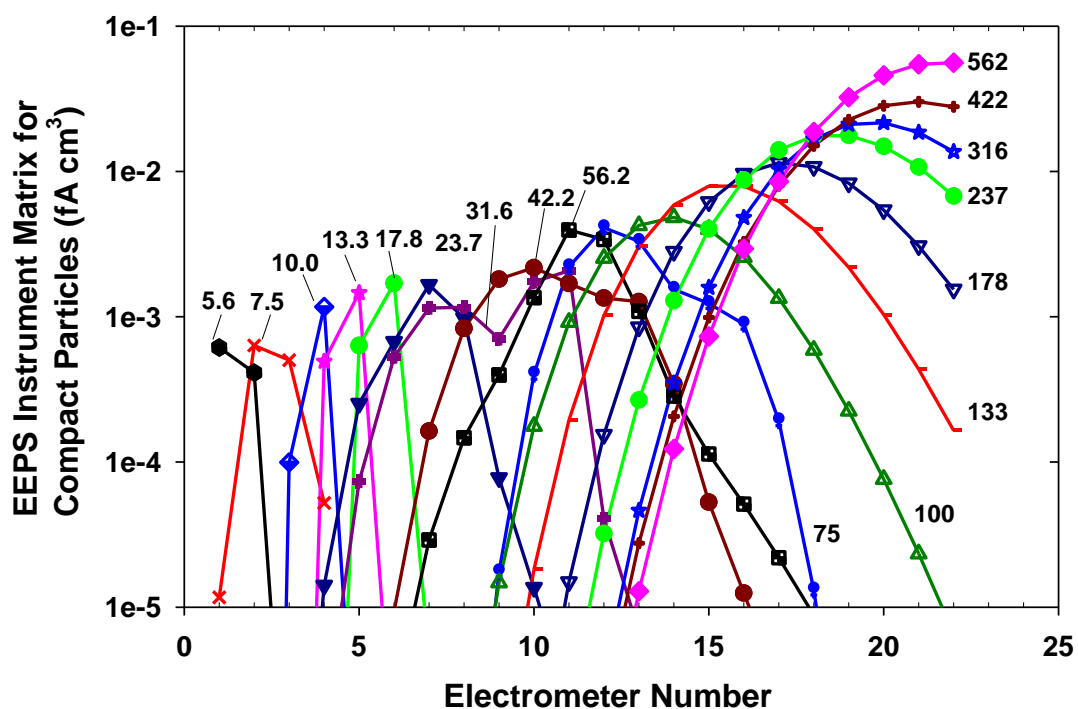


Figure 6: EEPS spectrometer "compact" inversion matrix (IM-compact) resulting from calibrating EEPS spectrometer to SMPS spectrometer measurements of compact aerosols. Each curve corresponds to a column of the instrument matrix. [6].

IM-Soot

As in the development of IM-compact, IM-soot was developed by calibrating the EEPS spectrometer to SMPS spectrometer measurements using monodisperse aerosols classified from a polydisperse aerosol. However, the initial polydisperse aerosol used for IM-soot was diluted diesel engine exhaust. Calibration was performed in conjunction with the University of Minnesota CDR using a John Deere 4045H (turbocharged, 4 cylinder, 4.5L, 75 kW, Tier 3) diesel engine fueled with BP6 diesel fuel with a sulfur content of 6 ppm.

As shown in Figure 5, the exhaust was sampled and diluted with cooled filtered air using an ejector diluter at a dilution ratio of approximately 4:1. The diluted exhaust was aged for approximately 1 second in a dilution tunnel before being classified by a TSI model 3080 electrostatic classifier. The monodisperse aerosol was then mixed with filtered air in a manifold before being simultaneously sampled by two EEPS spectrometers (one belonging to TSI and one to the University of Minnesota) as well as a 3775 CPC and 3936 SMPS spectrometer consisting of a 3080 classifier and a 3776 CPC. For sizes below 40 nm model 3085 nano-DMAs were used; for sizes above 40 nm model 3081 long DMAs were used. The engine was operated at a variety of speeds and loads to generate particle distributions across the full operating range of the engine and the full size range of the EEPS spectrometer.

Additionally, the exhaust sample was secondarily diluted at a ratio of approximately 12:1 and the resulting polydisperse distributions were sampled with the suite of instruments to evaluate the matrix performance.

To reduce the multiply-charged particles passing through the classifier, monodisperse sizes were consistently chosen to be larger than the mode particle diameter of the polydisperse distribution. Additionally, the contribution of multiply charged particles was subtracted from the resulting matrix as described in the above IM-Compact section and further detailed in [6].

The polydisperse distributions were controlled by the engine operating condition, and the monodisperse classified sizes were chosen to be just larger than the mode particle diameter; as a result, the sizes in the resulting matrix columns did not necessarily correspond with EEPS spectrometer primary channel boundary diameters ($i=1,2,\dots,17$ from Equation 1) as required by the EEPS spectrometer instrument matrix.

To develop instrument matrix columns at each of the primary channel boundary diameters, the following process was used and is further detailed in [7]. Lognormal fits [17] were applied to each of the measured matrix columns. The controlling parameters for each lognormal fit were plotted and best fit functions fit to the plotted parameters. These best fit functions were used to determine lognormal fit parameters at the EEPS spectrometer primary channel boundary diameters. These fit parameters were then used to develop an instrument matrix at the corresponding EEPS spectrometer primary channel boundary diameters. The resulting IM-soot matrix is plotted in Figure 7 where each curve corresponds to a column of the instrument matrix, $H_{j,i}$ from Equation 2.

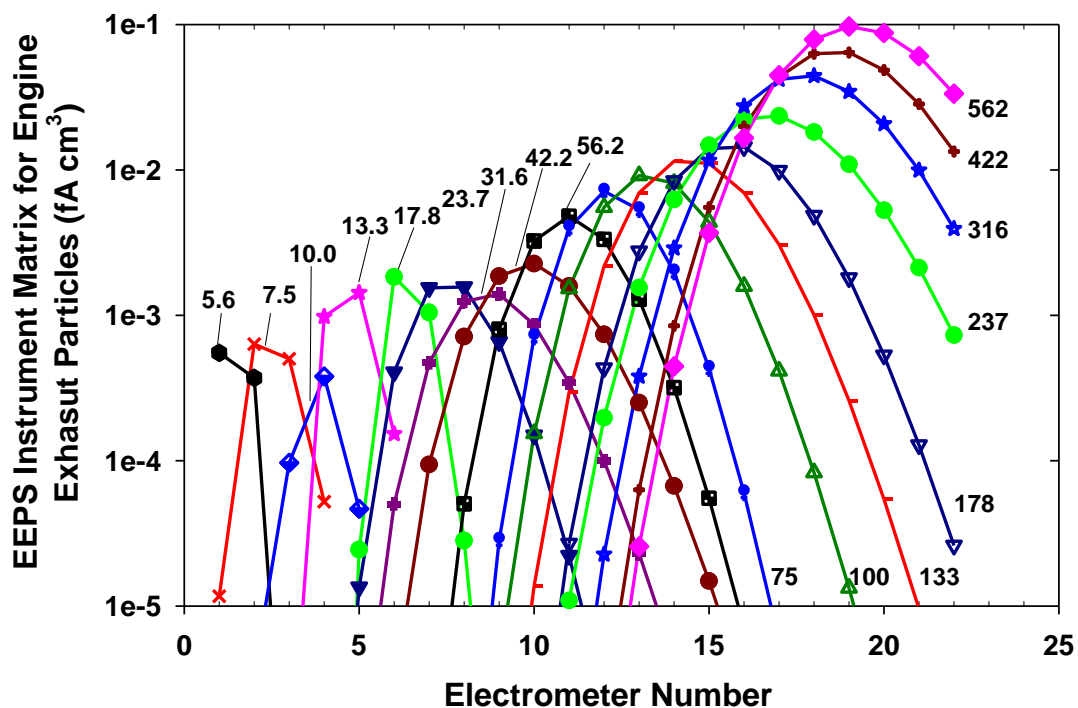


Figure 7: EEPS spectrometer soot inversion matrix (IM-soot) resulting from calibrating EEPS spectrometer to SMPS spectrometer measurements of agglomerate aerosols, specifically diesel engine exhaust. Each curve corresponds to a column of the instrument matrix. [7].

Comparison of EEPS Spectrometer Inversion Matrices

In Figure 6 and Figure 7 each curve represents a column of the instrument matrix i.e. the concentration normalized electrometer signal recorded for each of the monodisperse aerosols listed. As the mean particle size increases, the electrometer number on which the particles are detected also increases. Additionally, since the particle charge distribution broadens as size increases, the larger particles are detected on a wider range of electrometers, broadening the instrument matrix columns for larger sizes. It should also be noted that the IM curves are often overlapping, indicating that a single electrometer can detect particles of multiple sizes. This overlapping limits the EEPS spectrometer sizing resolution of the EEPS spectrometer to 16 channels per decade. Additionally, the y-axes of Figure 6 and Figure 7 are plotted on log scales to show the wide range of signal strengths on each electrometer for each size.

Figure 8 shows the 42 nm and 420 nm instrument matrix columns for the three matrices, IM-default, IM-compact, and IM-soot. Since with a unipolar charger, agglomerates obtain a higher state of charge than spherical particles, the agglomerates have a higher electrical mobility than their compact counterparts of equivalent size. Thus, the 420 nm column in IM-soot corresponds to a lower number electrometer than the 420 nm column in IM-compact or IM-default. Additionally, it can be seen that for the smaller 42 nm columns, since the agglomerates are more compact at this size, the difference between matrices is much less. The differences in the plotted columns for different aerosols are consistent with the discrepancies reported by others in the Background section, as agglomerates would show up as smaller particles using IM-default.

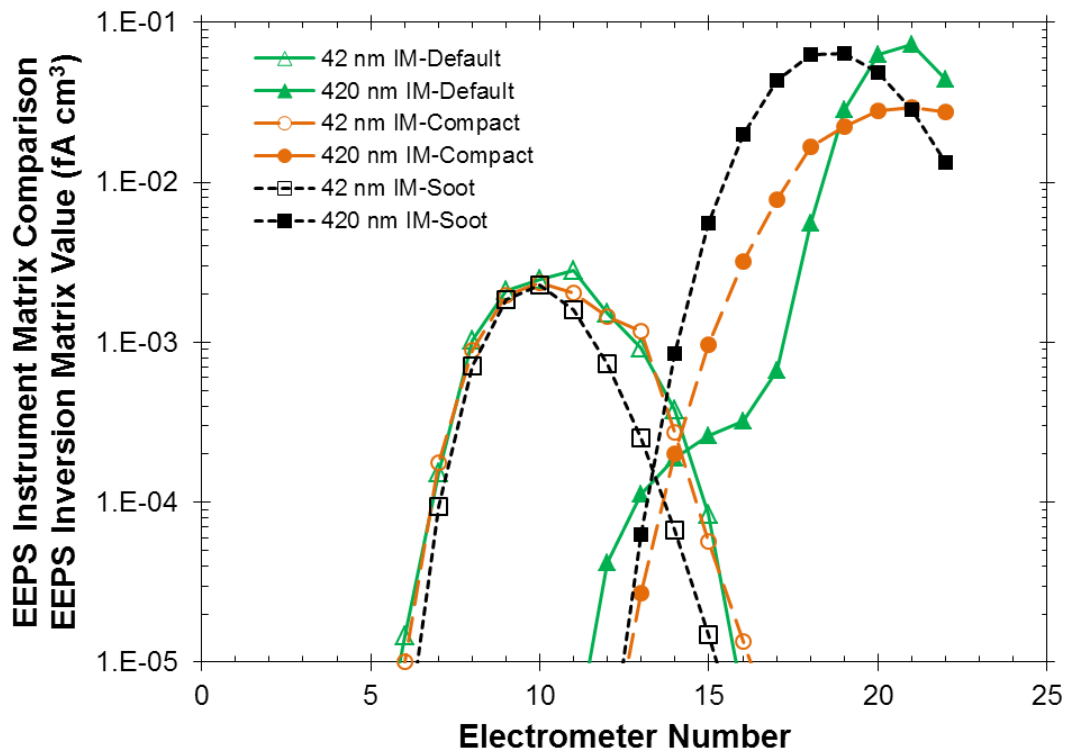


Figure 8:
Comparison of 420 nm and 42 nm EEPS spectrometer inversion matrix columns for IM-default, IM-compact, and IM-soot matrices

Results Using New Matrices

The IM-compact and IM-soot matrices were evaluated by comparing SMPS spectrometer measurements to EEPS spectrometer measurements processed with the compact and soot matrices for compact aerosols and engine exhaust, respectively. Multiple EEPS units were tested on a variety of aerosol types to determine the bulk effect of the alternative matrices.

Compact Inversion Matrix Results

A comparison of the EEPS spectrometer compact and default matrices measuring a variety of compact aerosols are presented in Figure 9. As shown in a, c, and e, the use of the EEPS spectrometer compact matrix results in much better agreement with the SMPS spectrometer size distributions. As the mode particle diameter increases, from 40 nm in a, to 160 nm in c, and 280 nm in e, the undersizing of the default matrix becomes apparent. However, if the compact matrix is used, the resulting EEPS spectrometer distributions much more closely match the SMPS spectrometer distribution in both sizing and concentration.

In addition to the undersizing issue shown with the default matrix, the entire distribution as measured by EEPS spectrometer is on average much narrower than the distribution as measured by SMPS spectrometer. As a result, when the data are plotted on a volume weighted scale, the default matrix greatly underestimates particle volume. When using the compact matrix, the discrepancy in volume is significantly reduced. Further details can be found in [6].

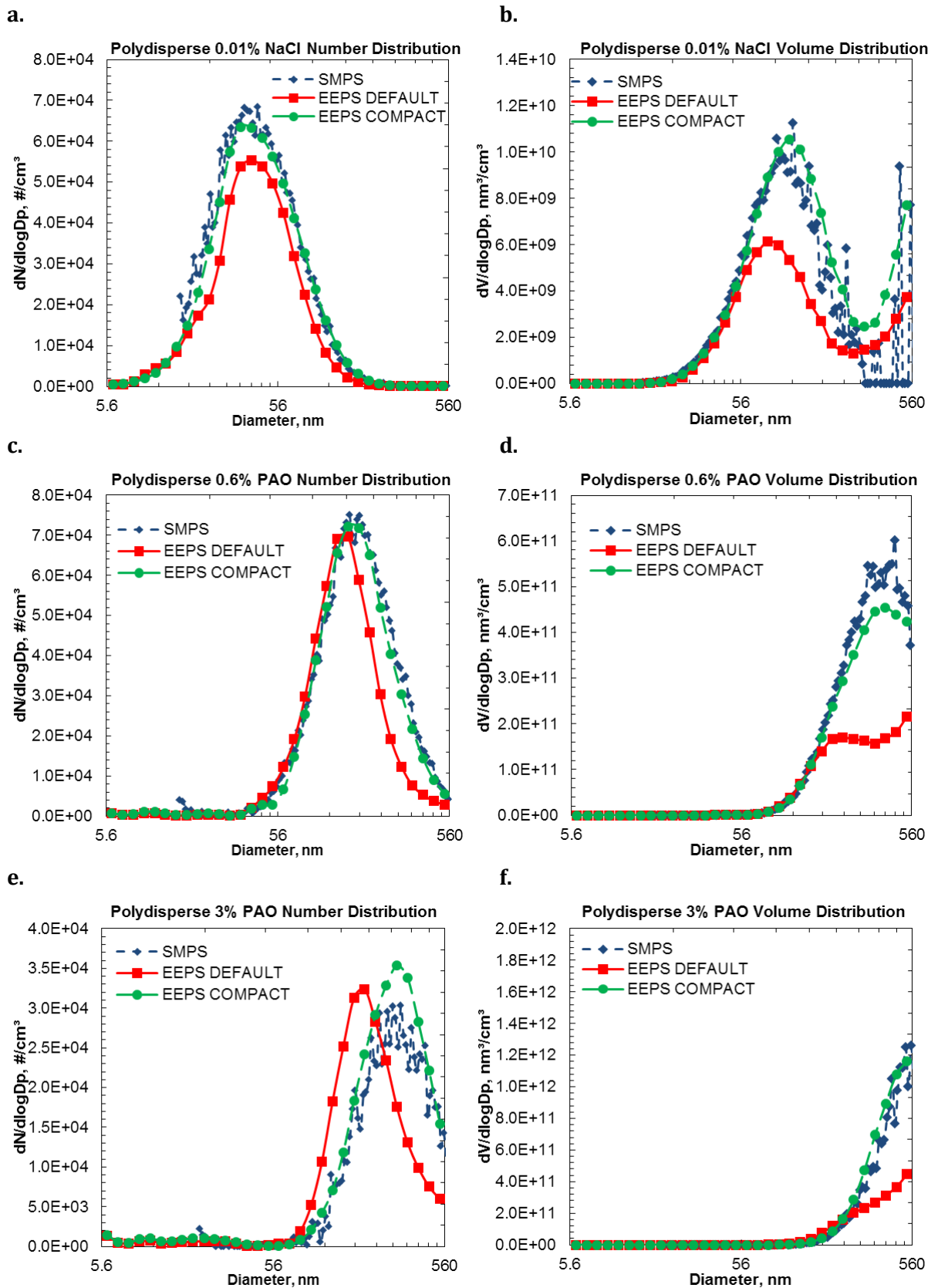


Figure 9: EEPS spectrometer Default and compact instrument matrix number and volume weighted particle size distribution compared to SMPS spectrometer for atomized NaCl (a, b) and atomized PAO (c,d,e,f) of various sizes.

Soot Inversion Matrix Results

Data were collected by EEPS spectrometer and SMPS spectrometer on multiple engine types running at several different conditions and compared to evaluate the performance of the EEPS spectrometer soot matrix. Figure 10 shows example measurements taken from a John Deere 4045h heavy duty engine running at two loads [7].

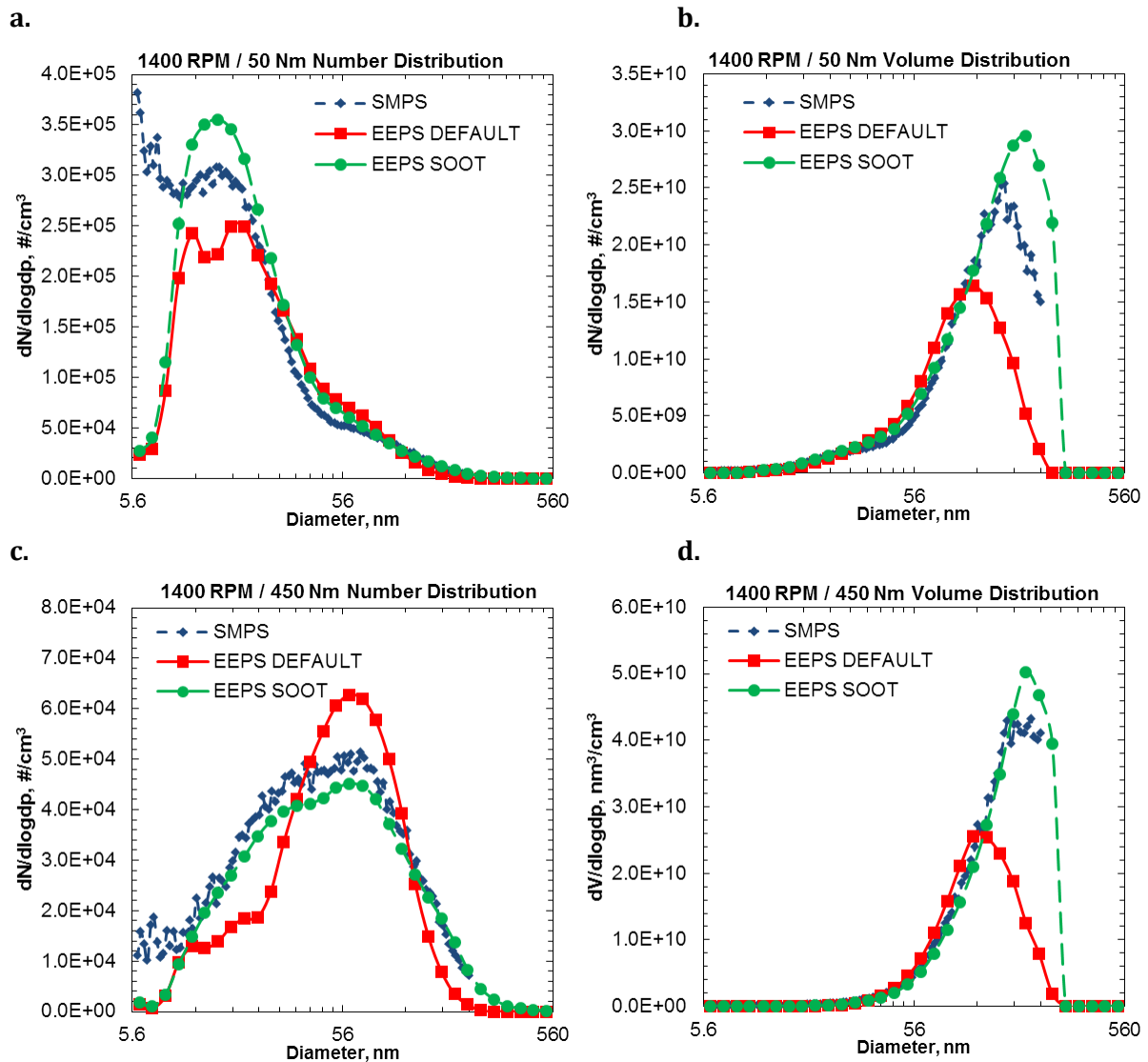


Figure 10:

EEPS spectrometer Default and soot instrument matrix processed number and volume weighted particle size distribution compared to SMPS spectrometer for a heavy duty diesel engine running at 1400 RPM, 50 Nm (a, b) and 1400 rpm, 450 Nm (c, d) [7].

For both the low load (a) and high load (b) engine conditions, the EEPS spectrometer data processed with the soot matrix more closely match SMPS spectrometer data than EEPS spectrometer data processed using the default matrix. Additionally because the sizing of particles larger than 75 nm is drastically improved, when weighted for volume, the EEPS spectrometer soot matrix again agrees with SMPS spectrometer much more closely. Further details can be found in [7].

Figure 11 shows example results from tests completed in conjunction with the University of Minnesota CDR on a GM A20DTH 2.0L light duty turbo charged diesel engine. The dilution and sampling system was similar to the one described in Figure 5b. When plotted on a number weighted scale, the mode diameter (as measured by SMPS spectrometer) is 18 nm. While the mode diameter of the EEPS spectrometer default matrix and EEPS spectrometer soot matrix are both within 5% (19 nm, and 17 nm respectively), the concentration of the default matrix is much lower than that of the SMPS spectrometer, while still reporting a lower concentration than the SMPS spectrometer, the soot matrix results in a distribution that more closely matches that of SMPS spectrometer. The same is true when plotted on a volume distribution.

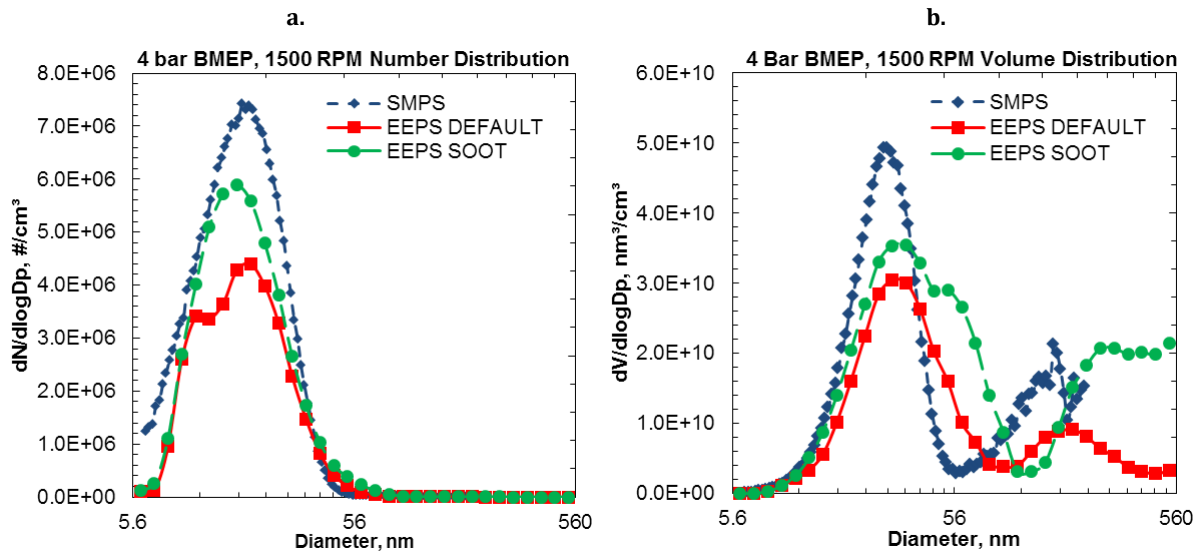


Figure 11: EEPS spectrometer Default and soot instrument matrix processed number (a) and volume (b) weighted particle size distribution compared to SMPS spectrometer for a light-duty diesel engine running at 1500 rpm and 4 bar BMEP.

Data were also collected from a Mazda 3 on a chassis dynamometer and are shown in Figure 12. As with the heavy duty and light duty diesel engines, the soot matrix more closely matches SMPS spectrometer measurements when number weighted as well as when volume weighted.

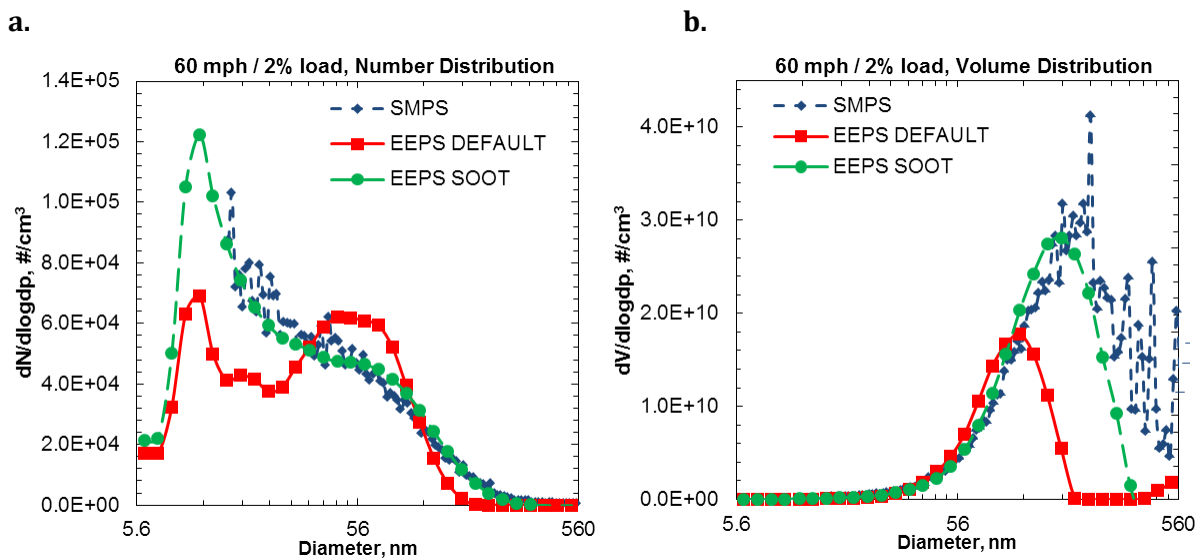


Figure 12: EEPS spectrometer Default and soot instrument matrix processed number (a) and volume (b) weighted particle size distributions compared to SMPS spectrometer for a passenger car running at 60 mph and 2% load [7].

As a result of improved sizing performance, the agreement between SMPS spectrometer measured volume distributions and EEPS spectrometer measured volume distributions is significantly improved. Together with an assumption of soot density, the use of the soot inversion matrix significantly improves the ability of EEPS spectrometer to measure soot mass.

Solid particle data were collected from the GM A20DTH at the University of Minnesota CDR a two stage partial flow dilution system and catalytic stripper. When compared over the same size range (i.e. greater than 23 nm) EEPS data processed with the soot inversion matrix agreed with TSI model 3790A engine exhaust CPC concentrations to within approximately $\pm 10\%$.

Overall Sizing Comparison

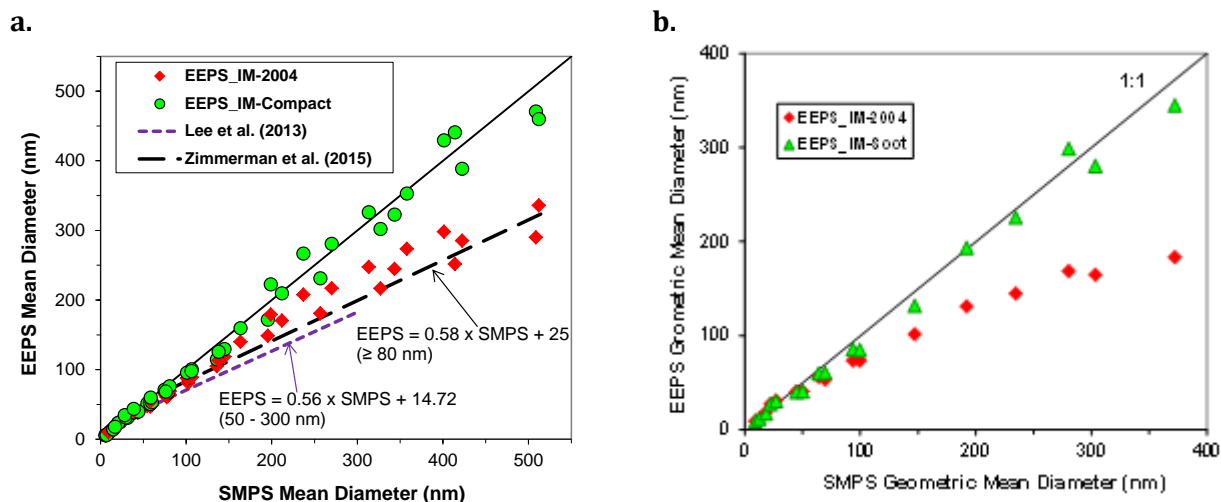


Figure 13: Mean diameter of number weighted particle size distributions as measured by EEPS spectrometer compared to those as measured by SMPS spectrometer for compact aerosols (a) and engine exhaust (b) processed with the default matrix (represented as IM-2004) and the compact matrix (a) and soot matrix (b) [6,7].

To examine the bulk sizing performance of the soot and compact EEPS spectrometer matrices, the mean diameter of several compact aerosols as measured by EEPS spectrometer (compact and default) were plotted against the mean diameter as measured by SMPS spectrometer in Figure 13a. Using the default matrix, the mean diameter as measured by EEPS spectrometer is roughly 60% of that as measured by SMPS spectrometer. When the compact matrix is applied, the data follow a 1:1 trend across the entire size range. Likewise, for several different engine models and operating conditions, the geometric mean diameters as measured by EEPS spectrometer using the soot matrix follow a 1:1 trend with the geometric mean diameters as measured by SMPS spectrometer as shown in Figure 13b. However, when the default matrix is used, the EEPS spectrometer geometric mean diameters are up to 50% smaller than those measured by SMPS spectrometer.

Software and Firmware Updates

TSI has released MCU firmware version 3.13 and software version 3.2.5. Firmware version 3.13 allows an inversion matrix other than the default to be uploaded onto the EEPS instruments itself, displaying data on the front panel processed with the uploaded matrix.

With EEPS spectrometer software version 3.2.5, data can be collected using IM-default, IM-compact, IM-soot, or any custom matrix. Additionally, software version 3.2.5 can be used to post-process previously recorded EEPS spectrometer data files provided that the files contain the raw instrument record. The raw instrument record is saved in any EEPS spectrometer data file recorded with software version 3.1.0 or later *and* MCU firmware version 3.05 or later. Data files collected on instruments with an MCU firmware version older than 3.05 and or software version older than 3.1.0 do not contain the raw instrument record and thus cannot be post processed using the EEPS spectrometer software.

The EEPS spectrometer MCU firmware version is displayed on the EEPS spectrometer front panel screen in the moments after turning on the EEPS spectrometer. The EEPS spectrometer software version can be found by opening the EEPS spectrometer software and navigating to Help>About EEPS spectrometer Software. A simple way to tell if a file contains the raw instrument record is to open the file using EEPS spectrometer software and navigate to the Export Data Options window by clicking **File>Export...** If "Raw Data (Instrument Record)" is displayed as an option in the "Data Types" section, as shown in Figure 14a, the file contains the raw instrument record and can be post-processed with an alternative inversion matrix in software version 3.2.5. If the option is missing as in Figure 14b., the raw instrument record was not recorded with the file and the data cannot be post-processed using an alternative matrix.

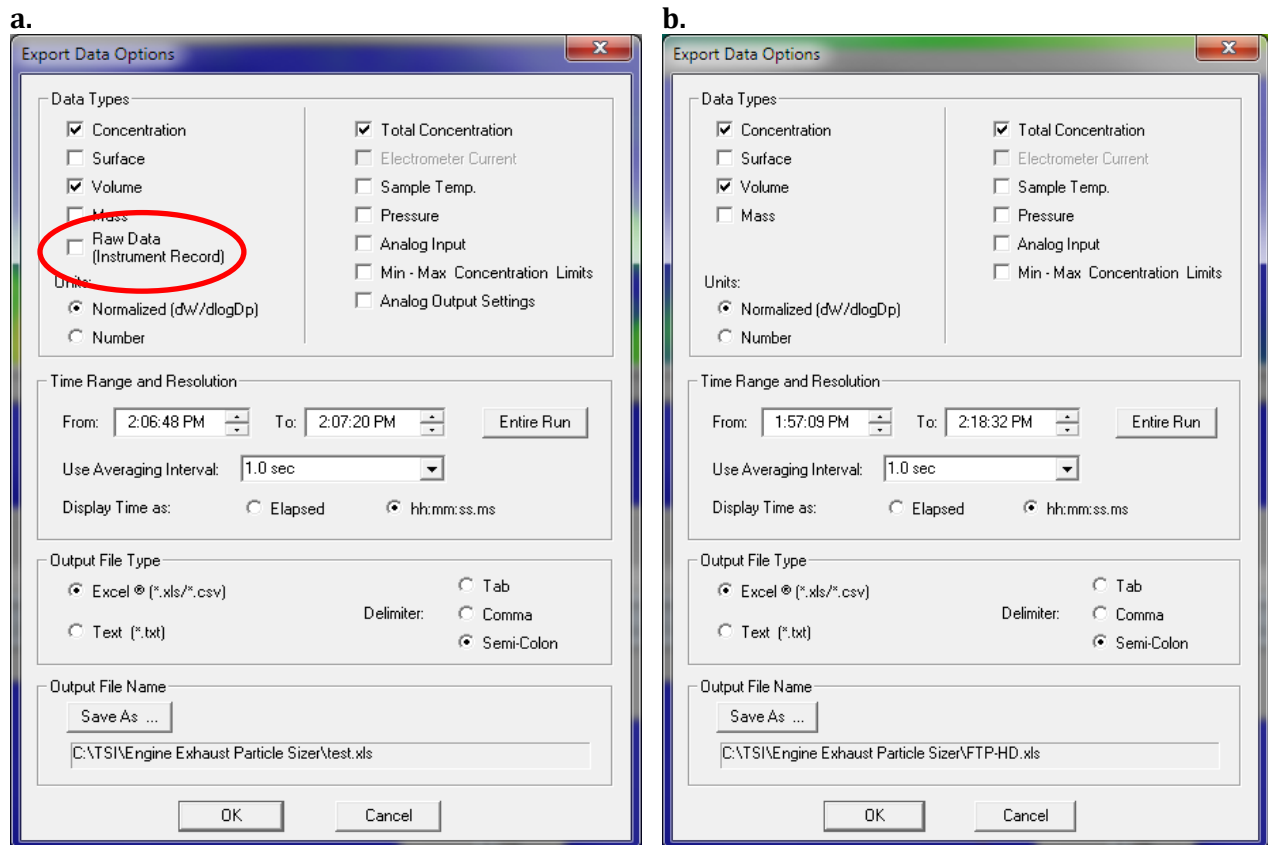


Figure 14: Export data screen for an EEPS spectrometer data file containing the raw instrument record (a.) and an EEPS spectrometer data file without the raw instrument record (b.).

Software version 3.2.5 and firmware version 3.13 are available as free upgrades. Contact TSI support at particle@tsi.com, with the serial number and firmware version of the instrument to be upgraded. To upgrade to firmware version 3.13, the existing EEPS spectrometer firmware version must be version 3.0 or greater. EEPS instruments with an MCU firmware version older than 3.0 cannot be upgraded to version 3.13 without being sent back to TSI for service due to hardware changes. Contact your TSI sales representative for information on upgrading a unit with firmware older than 3.0.

Using an Alternative Matrix in EEPS Spectrometer Software Version 3.2.5

Information on using an alternative EEPS spectrometer matrix can be found in the EEPS spectrometer instruction manual located in the EEPS spectrometer software installation directory. Briefly, to reprocess data from an existing file, open the data file and navigate to the properties window by clicking **File>Properties**. In the "Instrument Matrix" section choose from "DEFAULT", "COMPACT", "SOOT", or a custom matrix from the drop down menu to select one of the alternate matrices. Click **Apply** to reprocess the data file. The file can then be saved under a different file name or can be saved over the existing file.

Which Inversion Matrix Should I use?

EEPS spectrometer Software version 3.2.5 allows the user to collect, and process data using any of the three included matrices (default, compact, soot) in addition to any custom matrices that may be created. In general, the soot matrix should be used for engine exhaust and the compact matrix reserved for studies in which it is known that the aerosol being sampled is near-spherical in nature.

References

- 1 Asbach C, Kaminski H, Fissan H, Monz C, Dahmann D, Mulhopt S, Paur HR, Kiesling HJ, Herrmann F, Voetz M, et al. Comparison of four mobility particle sizers with different time resolution for stationary exposure measurements. *Journal of Nanoparticle Research*. 2009;**11**:1593-1609.
- 2 Jeong CH, Evans GJ. Inter-Comparison of a Fast Mobility Particle Sizer and a Scanning Mobility Particle Sizer Incorporating an Ultrafine Water-Based Condensation Particle Counter. *Aerosol Science and Technology*. 2009;**43**(4):364-373.
- 3 Kaminski H, Kuhlbusch TAJ, Rath S, Götz U, Sprenger M, Wels D, Polloczek J, Bachmann V, Dziurawitz N, Kiesling HJ, et al. Comparability of mobility particle sizers and diffusion chargers. *Journal of Aerosol Science*. 2013;**57**:156-178.
- 4 Lee BP, Li YJ, Flagan RC, Lo C, Chan CK. Sizing Characterization of the Fast-Mobility Particle Sizer (FMPS) Against SMPS and HR-ToF-AMS. *Aerosol Science and Technology*. 2013;**47**:1030-1037.
- 5 Leskinen J, Joutsensaari J, Lyyränen J, Koivisto J, Ruusunen J, Järvelä M, Tuomi T, Hämeri K, Auvinen A, Jokiniemi J. Comparison of nanoparticle measurement instruments for occupational health applications. *Journal of Nanoparticle Research*. 2012;**14**.
- 6 Wang XL, Grose MA, Avenido AS, Stolzenburg MR, Caldow R, Osmondson BL, Chow JC, Watson JG, Mirme A. Improvement of Engine Exhaust Particle Sizer (EEPS) Size Distribution Measurement - I. Algorithm and Applications to Compact Aerosols. *Journal of Aerosol Science*. 2015a Under Preparation.
- 7 Wang XL, Grose MA, Avenido AS, Caldow R, Osmondson BL, Swanson JJ, Chow JC, Watson JG, Kittelson DB, Li Y, et al. Improvement of Engine Exhaust Particle Sizer (EEPS) Size Distribution Measurement - II. Engine Exhaust Aerosols. *Journal of Aerosol Science*. 2015b Under Preparation.
- 8 TSI. Engine Exhaust Particle Sizer (EEPS) Spectrometer Model 3090/3090AK. Instruction Manual. Shoreview, MN: TSI; 2014 November.
- 9 Johnson T, Caldow R, Pöcher A, Mirme A. A New Electrical Mobility Particle Sizer Spectrometer for Engine Exhaust Particle Measurements. *SAE Technical Paper*. 2004.
- 10 Zimmerman N, Pollit KJG, Jeong CH, Wang JM, Jung T, Cooper JM, Wallace JS, Evans GJ. Comparison of three nanoparticle sizing instruments: The influence of particle morphology. *Atmospheric Environment*. 2014;**86**:140-147.
- 11 Zimmerman N, Jeong CH, Wang JM, Ramos M, Wallace JS, Evans GJ. A source-independent empirical correction procedure for the fast mobility and engine exhaust particle sizers. *Atmospheric Environment*. 2015;**100**:178-184.
- 12 Awasthi A, Wu BS, Liu CN, Chen CW, Uang SN, Tsai CJ. The Effect of Nanoparticle Morphology on the Measurement Accuracy of Mobility Particle Sizers. *Mapan - Journal of Metrology Society of India*. 2013;**28**:205-215.
- 13 Fuchs NA. On the Stationary Charge Distribution on Aerosol Particles in a Bipolar Ionic Atmosphere. *Pure and Applied Geophysics*. 1963;**56**(3):185-193.
- 14 Wiedensohler A. An Approximation of the Bipolar Charge Distribution for Particles in the Submicron Size Range. *Journal of Aerosol Science*. 1988;**19**(1):387-389.
- 15 Rogak SN, Flagan RC. Bipolar diffusion charging of spheres and agglomerate aerosol particles. *Journal of Aerosol Science*. 1992;**23**(7):693-703, 705-710.
- 16 Oh H, Park H, Kim S. Effects of Particle Shape on the Unipolar Diffusion Charging of Nonspherical Particles. *Aerosol Science and Technology*. 2004;**38**:1045-1053.
- 17 Hinds WC. *Aerosol Technology*. New York: Wiley-Interscience; 1998.
- 18 Xue J, Li Y, Wang X, Durbin TD, Johnson KC, Karavalakis G, Asa-Awuku A, Villela M, Quiros D, Hu S, et al. Comparison of Vehicle Exhaust Particle Size Distributions Measured by SMPS and EEPS During Steady-State Conditions. *Aerosol Science and Technology*. **49**:984-996, 2015.
- 19 Xue J, Li Y, Wang X, Durbin TD, Johnson KC, Karavalakis G, Quiros D, Hu S, Huai T, Ayala A, et al. Assessing particle mass emission from light-duty vehicles over a transient cycle using integrated particle size distribution (IPSD) method. In Preparation.

TSI and TSI logo are registered trademarks of TSI Incorporated.
Engine Exhaust Particle Sizer, EEPS, Fast Mobility Particle Sizer, FMPS, Scanning Mobility Particle Sizer, and SMPS are trademarks of TSI Incorporated.



UNDERSTANDING, ACCELERATED

TSI Incorporated – Visit our website www.tsi.com for more information.

USA	Tel: +1 800 874 2811	India	Tel: +91 80 67877200
UK	Tel: +44 149 4 459200	China	Tel: +86 10 8219 7688
France	Tel: +33 1 41 19 21 99	Singapore	Tel: +65 6595 6388
Germany	Tel: +49 241 523030		



Distributed by:
Kenelec Scientific Pty Ltd
1300 73 22 33
sales@kenelec.com.au
www.kenelec.com.au