

VAISALA

USER'S GUIDE

vLog 4.5

Vaisala Veriteq vLog Software



M211334EN-A

Copyright

Copyright ©2011 Vaisala Canada Inc. All rights reserved.

Trademarks

Vaisala Canada Inc, Vaisala Veriteq, Vaisala Veriteq vLog, and Vaisala Veriteq Spectrum are trademarks of Vaisala Canada Inc.

Windows is a registered trademark of Microsoft Corporation.

All other trademarks referred to are the property of their respective owners.

Disclaimer and Limitation of Liability

Vaisala Canada Inc and its subsidiaries assume no responsibility for any damage or loss resulting from the use of this manual.

Vaisala Canada Inc and its subsidiaries assume no responsibility for any loss or claims by third parties which may arise through the use of this software. Vaisala Canada Inc and its subsidiaries assume no responsibility for any damage or loss caused by deletion of data as a result of hardware malfunction. Be sure to make backup copies of all important data to protect against data loss.

Vaisala Canada Inc makes no warranties, either expressed or implied, regarding Vaisala Veriteq data loggers or software, its merchantability, or its fitness for any particular purpose. The exclusion of implied warranties is not permitted by some states. As such, the exclusion may not apply to you.

Vaisala Veriteq vLog has been tested using a variety of network systems. However, the large number of possible hardware and network configurations makes testing under every circumstance impossible. If you have trouble using Vaisala Veriteq vLog software, contact Vaisala Canada Inc.

Technical Support

For technical support in North America, please call 1-866-861-3388, or for customers outside North America, see "Getting Help" on page 9.

Email: veriteqsupport@vaisala.com

Document part number: M211334EN

Release date: July 2011

No copying, distribution, publication, modification, or incorporation of this document, in whole or part, is permitted for commercial purposes without the express written permission of Vaisala Canada Inc.

Vaisala Canada Inc
13775 Commerce Parkway
Richmond, BC V6V 2V4
Canada

Contents

Chapter 1: Introducing vLog	1
Overview	2
vLog SP or vLog VL?	3
Preparing to use Vaisala Veriteq vLog Software	3
Understanding the vLog Window	5
Using Menus & Keyboard Shortcuts	6
Using the Tool Bar	7
Conventions Used in this Document	9
Getting Help	9
Chapter 2: Getting Started	11
System Requirements	12
Installing vLog Software	12
Setting Up and Using vLog—Overview	14
Ensuring Loggers are Calibrated	15
Validating Loggers (vLog VL)	15
Connecting Loggers to PCs	16
Logging on to vLog SP	17
Logging on to vLog VL	17
Setting Audit Trail Options (vLog VL)	21
Setting Security Levels (vLog VL)	22
Working with User Accounts (vLog VL)	23
Setting up PC Communication Ports	25
Configuring Loggers	26
Linking Loggers to an Audit Trail File (vLog VL)	31
Configuring Logger Description	31
Setting Logger Sample Timing	33
Setting the Sample Interval	37
Enabling and Disabling Logger Channels	39
Configuring Channels to Work with Transducers	39
Configuring Warm Up Time	42
Configuring Logger Channels to Work with Thermocouples	43
Clearing Logger Samples	44

Configuring Loggers: Batch Setup (vLog VL)	46
Mounting Loggers	48
Using External Temperature Probes	49
Caring for Relative Humidity Sensors	50
Chapter 3: Transferring Logger Files	51
Transferring Data	52
Setting Transfer Preferences	52
Configuring the Default File Name Format	53
Transferring Logger Data to the PC	54
Logger File Creation Report (vLog VL)	60
Logger File ID Number (vLog VL)	63
Setting System Date and Time	64
Batch Transfer of Logger Data to the PC (vLog VL)	64
Chapter 4: Working with Graphs	67
Displaying and Creating Graph Files	68
Inserting Logger Files	68
Inserting Multiple Logger Files	70
Removing Channels from a Graph	70
Scrolling Through Graphs	71
Pinpointing Graph Values and Times	71
Changing the Time Zone	72
Adding Graph Titles	72
Selecting a Graph Font	73
Zooming In	73
Zooming Out	76
Formatting Measurement Units	77
Specifying the Y-axis Scale	77
Changing the Decimal Place Precision	78
Changing the Temperature Scale	78
Setting Sample Smoothing	79
Adding Horizontal Lines to Graphs	80
Verifying a Graph (vLog VL)	82

Copying and Pasting Graphs	83
Exporting Graph Data	83
Saving Graphs	84
Printing Graphs	85
Printing Line Numbers on Graphs	85
Printing Footers on Graphs	86
Print Preview	86
Printing Graphs	86
Chapter 5: Viewing Reports and Statistics	87
Historical Data Reports	88
Adding Statistics to Historical Data Reports	89
Customizing Headers, Footers and Font Sizes in Historical Data Reports	91
Audit Trail Reports (vLog VL)	94
Viewing and Filtering Audit Trail Reports	94
Configuring Date and Time Columns in Audit Trail Reports	97
Customizing Headers, Footers and Font Sizes in Printed Audit Trail Reports	97
Batch Setup Reports (vLog VL)	98
Customizing Headers and Footers in Printed Batch Setup Reports	99
Graph Security Problems Reports (vLog VL)	100
Copying and Pasting Report Data	101
Saving Reports	102
Printing Reports	102
Appendix: FAQs and Troubleshooting	103
Software Error Messages	103
Frequently Asked Questions	104
Glossary	109
Index	117

Chapter 1: Introducing vLog

The vLog 4.5 User's Guide includes the information you need to install, configure and operate vLog. You can also review the vLog Quick Guide for an overview of vLog installation steps.

This chapter includes:

- overview of Vaisala Veriteq vLog VL and vLog SP software features
- hardware and software requirements

Overview

Vaisala Veriteq vLog 4.5 software provides high-integrity recordings of temperature, humidity and other conditions. It can also provide you with the level of file security required by the pharmaceutical industry, depending on the security key you use to install vLog. Your security key provides you with access to either vLog SP features or additional vLog VL features.

Note: vLog SP is designed to be used with Vaisala Veriteq SP-Series Data Loggers, while vLog VL can be used with both SP- and VL-Series loggers.

With Vaisala Veriteq vLog 4.5 you can:

- Create electronic records of temperature, humidity and other conditions. Using vLog VL, you can also generate secure files and Audit Trail Reports that ensure logger files have not been altered or modified in any way.
- Create graphs and reports from logger files. Additionally, if you use vLog VL, you can ensure your graphs and reports are secure and traceable. Traceability is achieved using Administrative Options (security levels), the Audit Trail Report, and Logger File ID Numbers. These features help ensure that you do not accidentally use the wrong logger files in your graphs and reports.
- Produce detailed graph reports and, in vLog VL, document the secure, verified and calibration status of both the logger and the data file.
- Overlay logger files within the same graph to check relationships between logger files.
- Zoom-in and zoom-out for a customized view of the graph.
- Add manual and automatic maximum, minimum, and average horizontal lines to graphs.
- Use the on-screen graph cursor to pinpoint exact values, times and dates.
- Open multiple graph windows at one time.
- Create statistical summaries.
- View, customize and print tabular-format reports with dates, times, logged values, interval statistics, headers and footers.

- Export graph data and reports to external databases, word processors or spreadsheets.
- Manually or automatically scale graphs.
- Use multiple Y-axes scaling and ranges.
- Enable and disable input channels.
- Set logger start and stop times, sample intervals and other logger settings.

vLog SP or vLog VL?

While the core vLog features are the same, additional features are available if you install vLog using a vLog VL security key:

- Compatible with both Vaisala Veriteq SP- and VL-Series Data Loggers.
- Create secure graph and logger files - the security status of your data is identified on every report. You know whether files are secure or not, and if files have security problems or tampering issues.
- Data and user security - vLog uses Windows password authentication—domain level or local account authentication. Your IT staff can implement what they already use to protect other company assets.
- Audit Trail - A built-in audit trail report generator creates reports for all logger data or for a time period you specify. No changes to the system go unrecorded.
- Batch Setup and Transfer - Save time using the settings from one logger in additional loggers—start and stop time, sample interval, logger and channel descriptions, and channel flags. You can also configure the Batch Setup feature to support a range of COM ports.

Preparing to use Vaisala Veriteq vLog Software

Vaisala Veriteq vLog software supports one or more Vaisala Veriteq SP-Series Data Loggers, a Vaisala Veriteq cable (PC-IC cable if connecting by Serial port, INT-USB-DL cable if connecting by USB, or Digi and Ethernet cables if connecting via Digi or Vaisala Veriteq vNet device through the network), and a PC with an attached printer.

If you use vLog VL, you can also use Vaisala Veriteq VL-Series Validatable Data Loggers.

Once you connect loggers to a PC, you can use vLog software to configure logger settings such as the sampling interval, then place the loggers in the area where you want to measure conditions and record data. Your loggers collect the information and store it in their internal memory. Once the desired information is collected, you transfer it to a PC running vLog and import (“insert”) the data into a vLog graph. Use the graph to display, analyze, and print graphical or tabular reports, and, if needed, export the information to use in an analysis or presentation application.

VL-Series Validatable Data Loggers

Each Vaisala Veriteq VL-Series Validatable Data Logger is factory-calibrated against National Institute of Standards and Technology (NIST) traceable standards. All calibration information is stored digitally in the logger’s memory. For logger specifications, see www.vaisala.com/lifescience-hitech or contact Vaisala Canada Inc.

There are many models of VL-Series loggers to choose from, including the VL-1000, VL-1200, VL-1400, VL-1700, VL-2000, VL-4000 and VLT series. Each logger has one or more channels and can monitor and record data for FDA-regulated applications, including ambient temperature and relative humidity. Many loggers can also have optional plug-in external probes attached.

Understanding the vLog Window

The following illustration shows the main sections of the vLog window, with items described in Table 1:

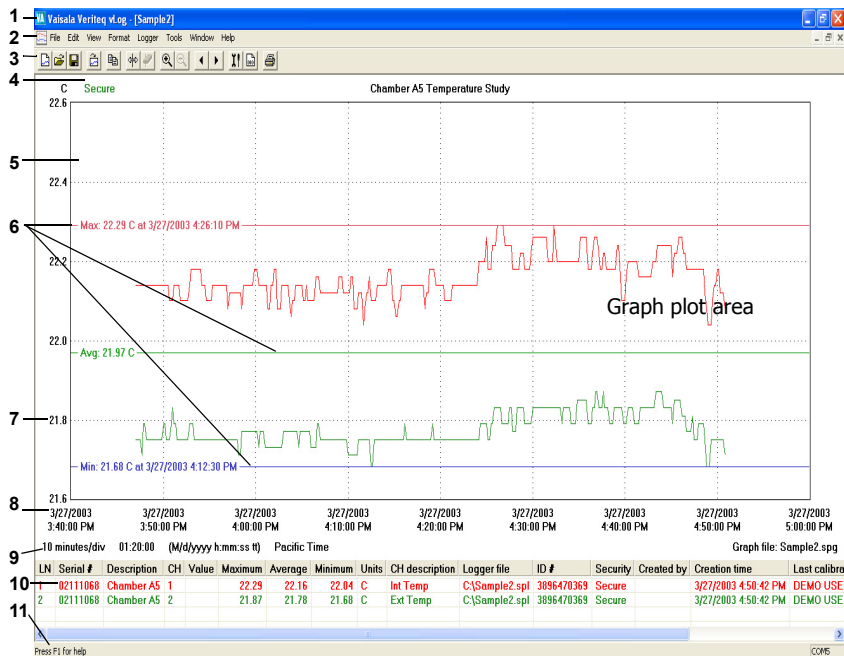


Table 1: Items in the vLog Window

1 - Title bar	Shows the name of the active Graph file
2 - Menu bar	Allows access to the various software functions
3 - Toolbar	Allows quick access to the most used software functions
4 - File security status	vLog VL: Shows the security level of the Graph file (note: If the checksum in a logger file or graph file is invalid, vLog will not open the file). Tampered means at least one logger is a VL-Series logger but one of the loggers has an invalid checksum. Secure means each logger is a VL-Series logger and all data and graphs have valid checksums. Not Secure means there are security problems, such as one SP-Series logger shown in the graph, or one of the graphs or logger files does not have a checksum.

Table 1: Items in the vLog Window

File verification status	(Not shown - vLog VL) Shows the verification level—only appears if a transfer password was used when creating the currently-displayed logger file (see page 106 for more information).
5 - Graph plot area	The Graph file is displayed here.
6 - Horizontal lines	Automatic: Lines showing maximum, minimum and average readings. Manual: Lines showing measurements set by you.
7 - Left-side Y-axis scale	Shows the scale for the data displayed in the graph (this example displays temperature).
Right-side Y-axis scale	(Not shown) This scale is used when two or more scales are displayed on the graph.
8 - X-Axis time scale	Shows the time scale.
9 - Graphing grid details	Shows the time per division, graph duration, date and time format and the time zone. The time zone is based on the time zone setting of the PC when the Graph file was created.
10 - Channel information	Shows the name(s) of the logger files currently being displayed, along with other identifying information including the logger serial number, logger description, logger channel number, measurement units, measurement parameter, and Logger File Name and ID Number. Also shows channel statistics such as the maximum, average and minimum values of data currently shown in the graph. The color shown for each channel corresponds to the color of the graph lines.
11 - Status bar	Displays information about the menu items and the toolbar. If the software is idle, the message Press F1 for Help is displayed. You can alternatively display or hide the status bar by selecting View Status Bar .

Using Menus & Keyboard Shortcuts

There are multiple ways to initiate tasks in vLog - using keyboard shortcuts, function keys, menu options, or toolbar icons. The instructions in the vLog User's Guide provide you with the menu option and additional options are indicated where applicable.

You may wish to take a few minutes to become familiar with the menu items, or right-click on various parts of the window to discover additional ways to access menu items.

Here are some common keyboard shortcuts and function keys:

- Create a new blank Graph File: Ctrl+N (**File | New**)
- Open an existing Graph file: Ctrl+O (**File | Open**)
- Save a Graph file: Ctrl+S (**File | Save**)
- Insert a Logger file: Ctrl+I (**File | Insert**)
- Print a graph file: Ctrl+P (**File | Print**)
- Produce an historical report associated with the active graph: F6
- Scroll Backward: Alt+Left arrow
- Scroll Forward: Alt+Right arrow
- Open the Logger Setup window: F11
- Transfer information from the logger to your PC: F12
- Help: F1 (opens an online, printable version of the User's Guide)

Using the Tool Bar

The vLog Toolbar allows quick access to the most commonly used features of the program. You can hide or display the toolbar by selecting **View | Toolbar**. The following table describes the function of each of the toolbar icons.

Table 2: Toolbar Functions

icon	Function
	Creates a new graph you can insert Logger files into. Same as File New.
	Opens an existing graph file (.spg)
	Saves the currently-active graph using the same file name.

Table 2: Toolbar Functions

	Inserts a Logger file into the currently-active graph.
	Copies graph to the Clipboard.
	Positions a vertical cursor line on the active graph to determine exact graph values and times.
	Removes the vertical cursor line from the active graph.
	Zooms-in for a closer look at graph data.
	Zooms-out.
	Scrolls the graph to the left. Time base is shifted backward by one graph division.
	Scrolls the graph to the right. Time base is shifted forward by one graph division.
	Opens the Logger Setup window.
	Transfers logger data to the computer and creates a Logger file.
	Prints the currently-active graph.

Conventions Used in this Document

This document uses the following conventions:

- Menu selections, items you select, and the names of buttons are shown in **bold**.
- A sequence of actions is indicated by a list separated by a vertical line. For example:
 'In vLog, choose **Tools | Options**'
- Click: To select a screen button with a click of the mouse.

Getting Help

If you need help, technical support is available.

North America

Contact Vaisala Canada Inc, 8am-4pm PST Monday to Friday, at 1-866-861-3388 (or 604-273-6850) or email veriteqsupport@vaisala.com. See also www.vaisala.com/lifescience-hitech.

For sales, pricing, quotations, or general information, please call 1-800-273-8374 (or 604-273-6850).

Vaisala Headquarters (Finland)

Vanha Nurmijärventie 21
01670 Vantaa
FINLAND
Industrial Instruments
Phone: +358 9 8949 2658
Fax: +358 9 8949 2295

Japan Service Center

42 Kagurazaka 6-Chome
Shinjuku-ku
Tokyo 162-0825
JAPAN
Phone: +81 3 3266 9611
Fax: +81 3 3266 9610

China Service Center

Floor 2, EAS Building

No. 21, Xiao Yun Road, Dongsanhuan Beilu

Chaoyang District

Beijing 100027

CHINA

Phone: +86 10 5827 4100

Fax: +86 10 8526 1155

Chapter 2: Getting Started

This chapter will assist you by providing information about:

- system requirements
- installing and logging on to vLog
- configuring loggers

System Requirements

vLog software requires a PC with the following minimum configuration:

- 32-bit or x64 version of: Microsoft Windows® XP, Windows Server® 2003, Windows Vista®, Windows Server® 2008, and Windows® 7
- one available Serial or USB communication port (for transferring logger data files; you can view graphs and reports for previously transferred logger files without this port).


Installing vLog Software

To install vLog software:


Note: vLog installation automatically uninstalls vLog versions 4.3 and 4.4; earlier versions of vLog are not uninstalled automatically. Also note that the uninstall program will not remove logger files (.spl), Graph files (.spg), Batch Setup files (.bsf) or Audit trail files (.vat).

- 1 Ensure you have administrator privileges.
- 2 Start Windows, log in as the administrator, and close all running applications. If you will be using the USB Logger Cable, do not connect the cable until installation is complete.
- 3 Insert the Vaisala Veriteq vLog CD into the CD-ROM drive and wait for it to start automatically. If it does not start automatically, from the **Start** menu, select **Run**. Type d:\setup and press Enter. If d:\ is not your CD-ROM drive, type the appropriate drive letter.

Select your install language, then click **Next**. The vLog Setup window opens.

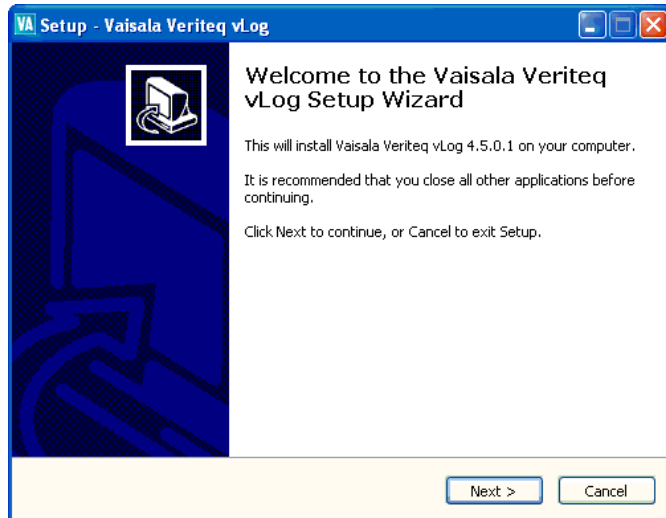

- 4 Click **Next**.
- 5 Review and accept the License Agreement, then click **Next**.
- 6 Enter your security key, including spacers (for example, xxxx-xxxx-xxxx-xxxx-xxxx).
- 7 Accept the default Destination Folder, or click **Browse** to select a new Destination Folder.

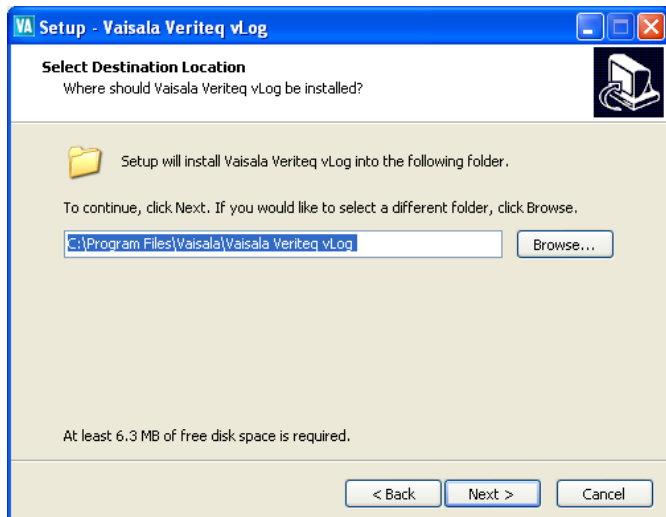

- 8 Click **Next**.
 - 9 Accept the default Start Menu Folder or click **Browse** to select a different folder. Click **Next**.
 - 10 If you want to add a shortcut to your desktop, select this option then click Next.
 - 11 Click **Install**.
 - 12 Click **Finish**.
- vLog is now installed.

Setting Up and Using vLog—Overview

Once you have installed vLog, you may find the following checklist a useful overview of steps for setting up and using vLog:

- 1 Ensure that loggers are calibrated (see page 15).
- 2 Validate loggers (see page 15).
- 3 Connect loggers to the PC (see page 16).
- 4 Log on to vLog software (see page 17).
- 5 Set Administrator Options (see page 21).
- 6 Select a communications port (see page 25).
- 7 Configure loggers (see page 26).
- 8 Place loggers in the area you want to validate and collect samples for the desired time (see page 48).
- 9 Transfer raw data in the form of a logger file (.spl) to the PC (see “Transferring Data” on page 52).
- 10 Insert logger files (.spl) into Graph files (.spg) (see page 68).
- 11 View, copy, customize and print Graphs (see **Chapter 4: Working with Graphs.**)
- 12 View, copy, customize and print Reports (see **Chapter 5: Viewing Reports and Statistics.**)
- 13 Have the printed Graph or Report approved and signed.

Ensuring Loggers are Calibrated

Vaisala Veriteq loggers are calibrated at Vaisala Canada Inc and should be calibrated on a suitable frequency as determined by the accuracy requirements of the application in which the product is used.

If you use vLog VL, the Calibration frequency should be specified when the loggers are purchased, and the next calibration date pre-set in the logger. Vaisala Canada Inc recommends a 90-day interval for the most critical applications, while a one-year interval may be adequate for less demanding situations.

To calibrate your loggers, contact Vaisala Canada Inc for pricing, lead times, and delivery arrangements.

Note: Treat the calibration process very carefully. The consequences of poor calibration can be very expensive. Vaisala Veriteq loggers are high-accuracy instruments. In order to maintain the specified accuracy, it is essential to have trained calibration personnel, strict calibration procedures and proper test and calibration facilities. If you have questions about the calibration process, contact veriteqsupport@vaisala.com.

Validating Loggers (vLog VL)

For companies that are required to validate the vLog system (including VL-Series Validatable Data Loggers and vLog VL software) prior to use, Vaisala Canada Inc provides a comprehensive protocol document.

Each package includes installation qualification (IQ) and operational qualification (OQ) procedures, and forms. The information and forms are provided as hard copies in a 3-ring binder and in electronic form (Microsoft Word documents) on a CD-ROM. The electronic files provide an easy way for you to customize the forms for the application.

For more information, contact Vaisala Canada Inc.

Connecting Loggers to PCs

To connect your logger to a PC using USB or Serial cable:

- 1 Connect the logger to a Vaisala Veriteq USB or Serial cable.



- 2 Connect the other end of this cable to a USB or Serial port on the PC.



To connect your logger to a PC using a Digi or vNet device:

- 1 Connect Digi or vNet device (will be referred to as device, following) to power supply and Ethernet outlet.
- 2 Connect logger to device (Digi requires Vaisala Veriteq cable).
- 3 Obtain a static IP address for device from your IT department. If your networking policy requires you to reserve IP addresses using DHCP, see www.vaisala.com/lifescience-hitech for technical support.
- 4 Insert device driver CD into PC. The Device Setup Wizard launches automatically. Click **Next**.
- 5 Select the device that matches the MAC address printed on the bottom of your device. Click **Next**.


Note: If your device is on a different subnet from the PC, you will need to know the IP address of the device and launch the RealPort Setup.exe installer from the RealPort folder on the device drivers CD.

- 6 In the Configure Network Settings window, enter an IP address (provided by your IT department). Click **Next** two times.
- 7 In the Configure RealPort Settings screen, select **Install (Digi) RealPort on this computer**. Click **Next**.
- 8 Click **Next** again. The settings are saved.
- 9 Click **Finish**. Drivers required to connect to your logger through the device have been installed.

Logging on to vLog SP

To log on to vLog SP:

- Double-click the vLog desktop icon, , or select **Start | Programs | Vaisala | Vaisala Veriteq vLog | Vaisala Veriteq vLog**. An empty vLog graph opens:

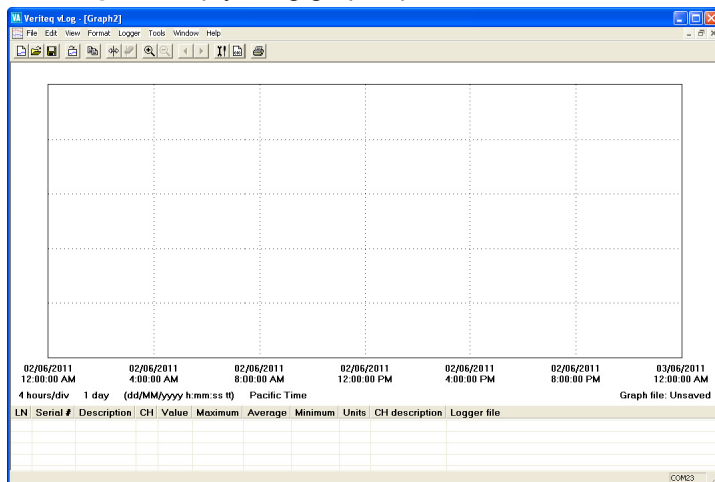


Logging on to vLog VL


The first time you run vLog VL, you need to run the vLog security wizard, which asks you to enter a security key and confirm that you have a valid Windows Administrator user name and password. You will also select where to save the audit trail file and what level of security to apply to users of this software.

To log on to vLog VL for the first time:

- 1 From the Start menu, choose **Vaisala | Vaisala Veriteq vLog | Vaisala Veriteq vLog**.
- 2 The vLog security wizard confirmation window opens, prompting you to create a new vLog security file. Click **Yes**.

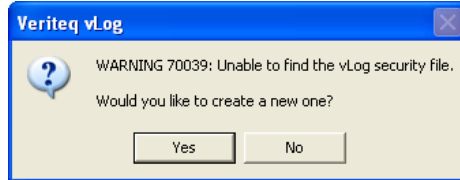

- 3 The vLog security wizard begins. Follow the instructions carefully. Click **Next**.
- 4 Enter the vLog Security Key that you received in your installation package (you do not need to enter the dashes).

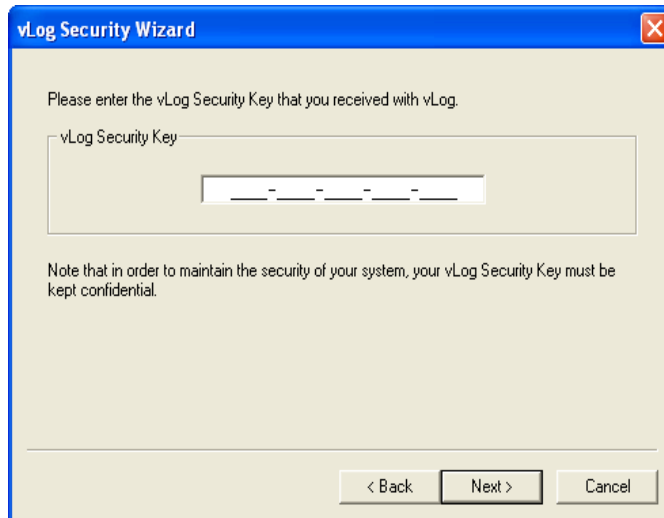

- 5 Click **Next**.
- 6 The next screen prompts you to create a vLog administrator account.
 - Enter an existing Windows Administrator User name.
 - Enter the password (this field is case-sensitive).
 - In the **Log on to** drop-down list, select the computer or domain you want to log on to.



The next step is to create a vLog administrator account. Note that the specified account must also be a valid Windows account.

Administrator Account

User name: (e.g. jsmith)

Password:

Log on to: MYDESKTOP (this computer) ▼

Full name: (e.g. John Smith)

< Back Next > Cancel

- Enter your full name.
- 7 Click **Next**. The Security Wizard shows you the path it intends to use for your audit trail folder.



The next step is to specify the audit trail folder. Note that all vLog users require write access to this folder.

Audit Trail Folder:


...

< Back Next > Cancel

Specify the folder where you want to keep the audit trail file, either on the local computer or stored in a network location. You can modify this option later using **Tools | Administrator Options**. For more information, see “Setting Audit Trail Options (vLog VL)” on page 21.

The audit trail file is a log of operations performed in vLog by various users. Multiple vLog installations on different PCs can use the same audit trail file, and it is recommended that all VL-Series loggers be linked to an audit trail file. If multiple PCs use the same loggers, they should point to the same audit trail file.

- 8 Click **Next**.
- 9 Select the security level you want to apply to all users with access to this software. Choose one of the following levels:
 - **Security off** - All users will have full access
 - **Normal security** - Users will be prompted to enter their user name and password when vLog opens
 - **High security** - Similar to Normal security, with the additional requirement that users enter their password whenever they change a logger's configuration.

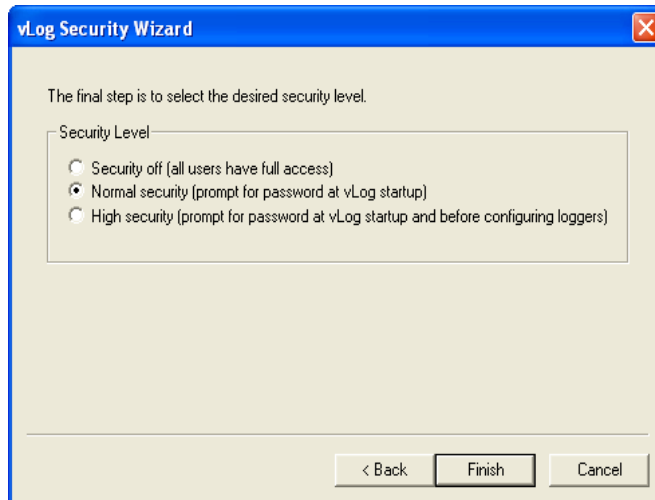


You can modify this option later (see “Setting Security Levels (vLog VL)” on page 22.

- 10 Click **Finish**.

To log on to vLog VL (after successfully logging on the first time):

- 1 From the Start menu, choose **Vaisala | Vaisala Veriteq vLog | Vaisala Veriteq vLog** (or use the desktop icon if you have one). The Log On to vLog window opens.

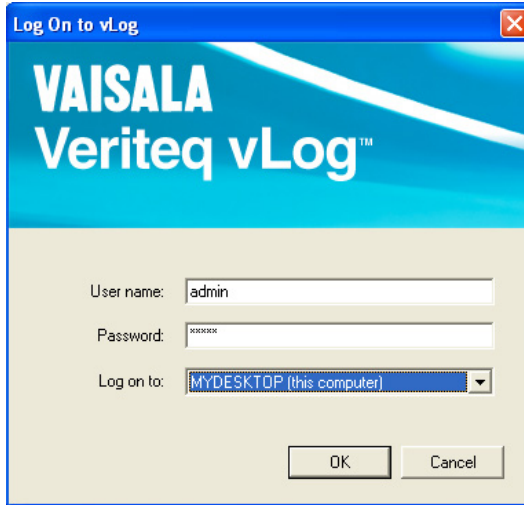

- 2 Confirm your user name or enter a different name if you are switching between users.
- 3 Enter your password.

Note: Because vLog VL log on security is tied to Windows security, you may be locked out of vLog temporarily if you exceed the maximum number of incorrect log on attempts (4). Try to use the correct user name and password when logging on.

- 4 Select the computer or domain you are logging onto.
- 5 Click **OK**.

vLog VL opens, displaying an empty graph window. Next, set your audit trail path, security options and create user accounts.

Setting Audit Trail Options (vLog VL)

If you have Administrator privileges in vLog VL, you can use Administrator Options to change vLog security levels, add

and change user accounts and roles, and change where the audit trail file is saved.

To change the path to the audit trail file:

- 1 In vLog VL, choose **Tools | Administrator Options**.
- 2 From the Audit Trail tab, change the folder where the audit trail file is saved.
- 3 Click **OK** to save, or **Cancel** to close without saving.

Setting Security Levels (vLog VL)

vLog VL has three security levels:


- **Security off** - all users will have full access.
- **Normal security** - users will be prompted to enter their user name and password when vLog opens.
- **High security** - similar to Normal security, with the additional requirement that users enter their password whenever they change a logger's configuration.

User accounts also have roles which define their task permissions: Administrators can change security levels and other administrative options, while Users and Guests cannot. For more on user accounts and roles, see "Working with User Accounts (vLog VL)" on page 23.

To change user security levels:

- 1 In vLog VL, choose **Tools | Administrator Options**.

- 2 From the Security tab, select either **Off**, **Normal**, or **High** security. This setting was originally set when you ran the vLog Security Wizard.

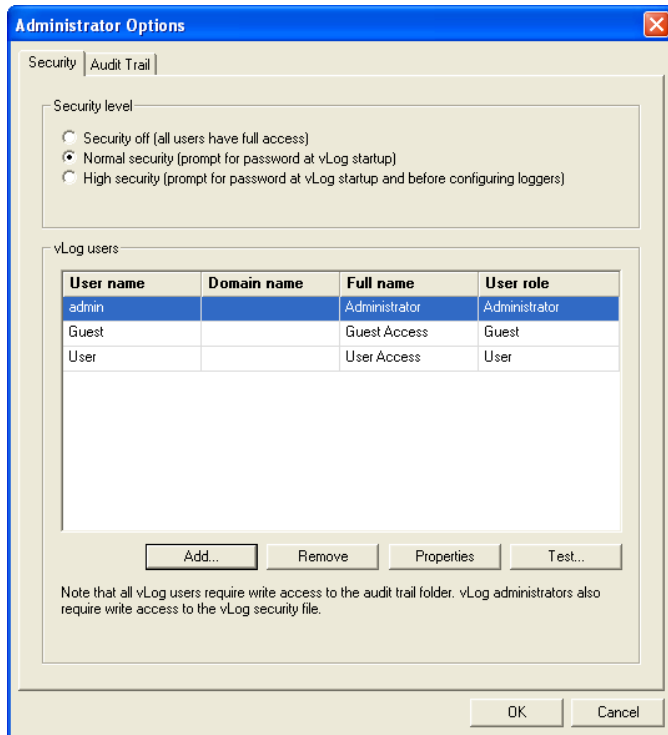

- 3 Click **OK** to save, or **Cancel** to close without saving.


Working with User Accounts (vLog VL)

To log in to vLog VL, users must have an account in vLog VL and on the PC where vLog VL is installed. User Accounts also have defined roles, including Administrator (full use of vLog; can change administrator options), User (full use of vLog; cannot change administrator options) and Guest (read-only access to vLog). Create user accounts for people who will be working with vLog VL.

Security levels are different from user account roles. For more on security levels, see “Setting Security Levels (vLog VL)” on page 22.

To add a user account in vLog:

- 1 In vLog, choose **Tools | Administrator Options**.
- 2 Click the **Security** tab. Here, you can view, add, and change users currently authorized to access vLog.
- 3 Click **Add**. The Add User window opens.

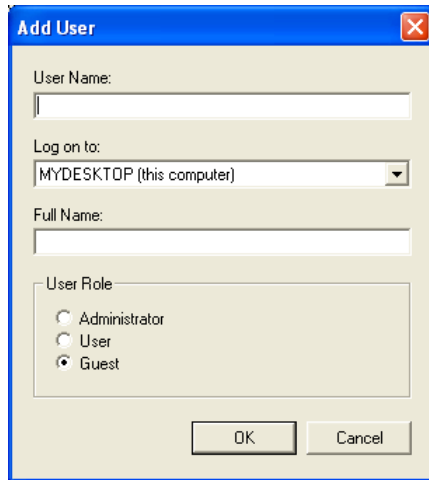

- 4 In the Add User window, enter the User Name and Full Name. This User Name must be a valid User Name on the PC or domain.
- 5 Select the user's role, which will determine the user's ability to make changes in vLog. Choose one of the following roles:
 - **Administrator** - full use of vLog; can change administrator options
 - **User** - full use of vLog; cannot change administrator options
 - **Guest** - read-only access to vLog
- 6 Click **OK** to save, or click **Cancel** to close without saving.

To change a vLog user account:

- 1 In vLog, choose **Tools | Administrator Options**.
- 2 Click the **Security** tab.
- 3 To view and/or change a user's properties, highlight the user name and click **Properties**. The User Properties window opens. You can verify or modify the fields.
- 4 Click **OK** to save, or **Cancel** to close without saving.

To remove a user account from vLog:

- 1 In vLog, choose **Tools | Administrator Options**.
- 2 Click the **Security** tab.
- 3 To remove a user, highlight the name and click **Remove**.
- 4 Click **OK** to save, or **Cancel** to close without saving.

Setting up PC Communication Ports

Before using a Vaisala Veriteq logger with vLog, you must select the correct PC serial communication (COM) port. vLog can support up to 4096 COM ports, numbered COM1 to COM4096.


To select a COM port:

- 1 Using Windows Device Manager, determine which COM port has been allocated to your logger:
 - Windows XP: **Start | Settings | Control Panel | System | Hardware | Device Manager**
 - Windows Server® 2003: **Start | Control Panel | System | Hardware | Device Manager**
 - Windows Server® 2008 or Windows Vista®: **Start | Control Panel | System | Device Manager**
 - Windows® 7: **Start | Control Panel | Hardware and sound | Device Manager**

Under Ports, make note of the COM port number allocated to the logger.

- 2 In vLog, choose **Tools | Options**, then choose the **General** tab.

- 3 From the **COM port** drop-down list, select the COM port used by the logger.

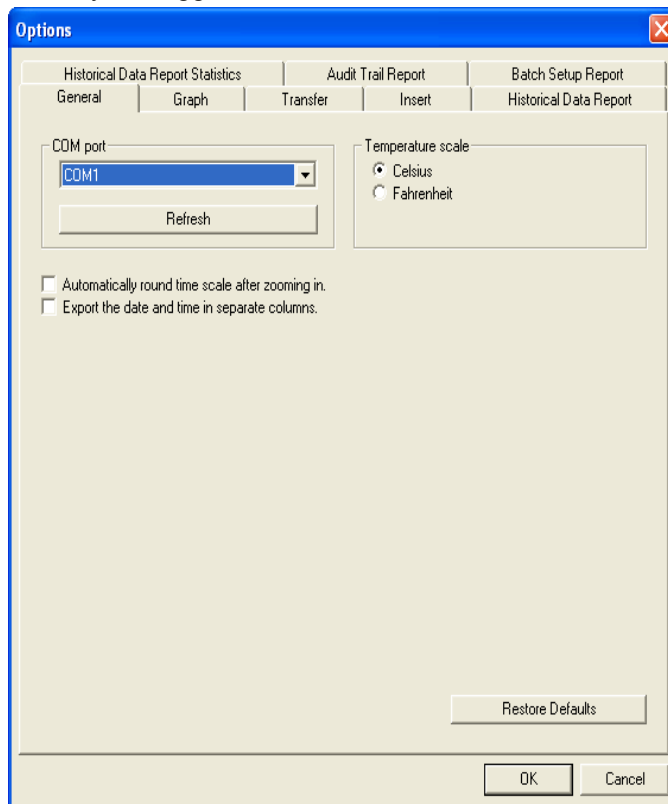

Note: Click **Refresh** to view a list of currently available COM ports.

- 4 Click **OK**.

Once you have selected a COM port, you do not need to change the setting unless you change the communication port to which your loggers are connected.

Configuring Loggers


Setting up a logger varies slightly depending on the logger you are using and the configuration options you prefer. You can configure loggers one at a time (as follows), or do a

Batch Setup of multiple loggers (see “Configuring Loggers: Batch Setup (vLog VL)” on page 46).

Note: If you use vLog VL and have high security enabled, you will need to confirm your password to finalize all logger configuration changes. For more on security, see “Setting Security Levels (vLog VL)” on page 22.

To configure loggers individually:

- 1 If you have not already done so, connect the logger to the PC.
- 2 Select **Logger | Setup (X)**.
- 3 If communication with the logger is established, the Logger Setup window opens.

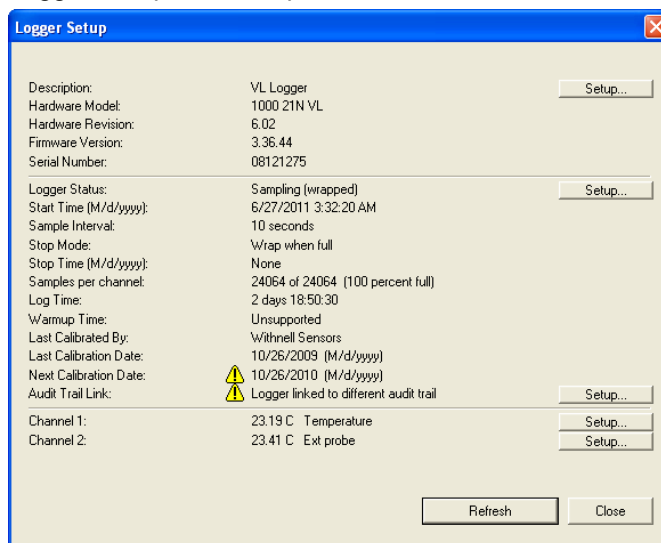

The information shown in the Logger Setup window relates to the logger connected to the COM port you have specified. Some fields shown will vary depending on the model of logger in use.

Note: If your logger has factory-disabled channels, they may not display at the bottom of the window in the list of Channels, or they will display as 'Disabled' (no **Setup** button will be available).

You can use this window to:

- Edit the logger description (see page 31).
- Set the sample timing (see page 33).
- Enable or disable logger channels (see page 39).
- vLog VL: Link loggers to audit trail file (see page 31).

Any changes will only affect the logger currently connected to the COM port.

Understanding the Logger Setup Window

The following table describes the information displayed on the Logger Setup window:

Table 3: Logger Setup window information

Description	Logger description (For information on configuring the logger description, see “Configuring Logger Description” on page 31.)
Hardware Model	Hardware model number of the currently-connected logger
Hardware Revision	Hardware version number associated with the currently connected logger
Firmware Version	Firmware version associated with the currently connected logger
Serial Number	Unique eight-character product serial number of the currently-connected logger. This number should match the serial number label on the back of the logger.
Logger Status	Indicates the sampling status of the currently-connected logger. The possible modes are: • Sampling: The logger is currently taking samples and the memory is not full. • Sampling (wrapped): The logger is currently taking samples with the memory full. The oldest reading is discarded and replaced with the newest.

Table 3: Logger Setup window information

Logger Status	• Stopped at specified time: The logger is not sampling and has stopped because the user-set stop time has been reached. • Stopped when full: The logger is not sampling and has stopped because it is configured to stop when the memory is full. • Start at specified time: The logger is not sampling, but is set to begin sampling as soon as the Start time is reached. No enabled channels: The logger is not sampling because no channels are enabled.
	The Setup button to the right of the Logger Status enables you to clear the logger's memory and change the Sample Interval , Start Mode , and Stop Mode settings.
Start Time	Indicates either the time that the currently-connected logger began taking samples, or the time (in the future) that it is set to begin taking readings Future Start times are indicated by a ⚠ symbol. For information on changing the Start Time, see "Setting the Start Mode" on page 34.
Sample Interval	Indicates the current sample interval setting (how frequently the logger is programmed to take readings). For information on setting sample intervals, see "Setting the Sample Interval" on page 37. • Note that the most frequent sampling rate possible with VL-1000-VLT, VL-1016, VL-1200 and SP-1200, SP-1016, VL-1416 and SP-1416 loggers is once per minute. For other loggers, the most frequent sampling possible is once every 10 seconds; however, sampling of this frequency will have a negative effect on battery life.
Stop Mode	Indicates how the currently-connected logger is set to stop while recording samples in its memory. There are three choices: Wrap when full, Stop when full, and Stop at specified time. For more information on setting Stop modes, see "Setting the Stop Mode" on page 35.

Table 3: Logger Setup window information

Stop Time	Indicates the programmed Stop Time , if set. If a Stop Time has not been set, it reads None . For information on setting the Stop Time , see “Setting the Stop Mode” on page 35.
Samples per channel	Identifies how many samples have been stored for each enabled channel on the currently connected logger, the capacity available (in samples) for each channel, and the percentage of memory used. Click the Refresh button to update values.
Log Time	Identifies the length of time the logger has been recording. The value is based on the number of samples in the memory multiplied by the sampling interval. Click the Refresh button to update values.
Warmup Time	This function applies to VL-4000 and SP-4000 Data Loggers. See “Configuring Warm Up Time” on page 42.
Last Calibrated By	vLog VL: Indicates the name of the organization that last calibrated the currently-connected logger, for example, Vaisala Canada Inc.
Last Calibration Date	vLog VL: Indicates the date of the last calibration.
Next Calibration Date	vLog VL: Indicates the expiry date of the previous calibration; it is recommended you get your logger calibrated before this date.
Audit Trail Link	vLog VL: If a high level of traceability is important to you, link each logger to a specific audit trail file. Do this using the Audit Trail Link setup button on the Setup window.
Channel 1	Temperature or Relative Humidity or other measured condition.
Channel 2 (there can be more than 2 channels)	Temperature or Relative Humidity or other measured condition. Vaisala Veriteq 4000-Series Data Loggers have several channels, each which can be configured for use with various transducers. See Configuring Channels to Work with Transducers on page 39.

Table 3: Logger Setup window information

	Click the Refresh button to update the currently displayed values. Click the Setup button to enable and disable channels. If a channel is disabled and no Setup button is available, the channel was disabled at the factory and cannot be enabled by the user.
--	---

Linking Loggers to an Audit Trail File (vLog VL)

If you use vLog VL and have enabled security, it is very important to link loggers to a specific audit trail file so you can see the exact configuration changes that have occurred over time. This linkage only needs to be done once per logger. You will know a logger is not linked to an audit trail when you open the Logger Setup window and see the alert, "Logger is not linked to audit trail."

To link a logger to an audit trail file:

- 1 From **Logger | Setup**, click the **Setup** button in the row for Audit Trail Link.
- 2 You are prompted to link this logger to the audit trail file. Click **Yes**.


Now that the logger is linked to an audit trail file, changes made to the logger settings are tracked in the audit trail. To make changes to a logger linked to an audit trail, a user must first log on to a computer connected to the same audit trail.

Configuring Logger Description

The Logger Description helps identify a logger every time you communicate with it. Using a description that refers to the application or location of the logger, and checking the logger description prior to each use, is a quick way to ensure you are connected to the correct logger.

To configure the Description for the currently-connected logger:

1 Select Logger | Description.

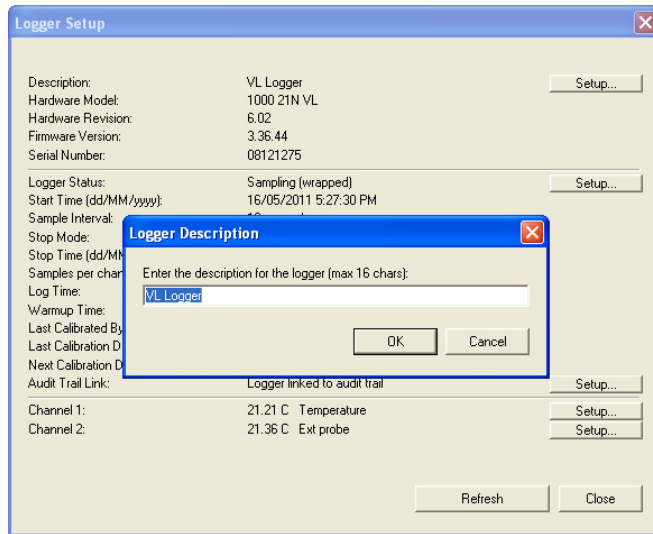

2 In the box provided, type a description of the logger, using up to sixteen alphanumeric characters (the space available for the logger description may be less than 16 characters, depending on your regional computer settings).

3 Click **OK** to save or **Cancel** to close without saving.

The text from the **Description** field is displayed on all on-screen graphs, tabular displays and printouts. The text can also be used in the creation of a default file name when transferring logger data to a PC if you choose to include it.

Note: You can also modify channel descriptions which display on logger files and graph files (see page 56).

Setting Logger Sample Timing

There are four main steps in setting the Logger Sample Timing:

- a Set the **Clear Mode**.
- b Set the **Start Mode** (see page 34).
- c Set the **Stop Mode** (see page 35).
- d Set the **Sample Interval** (see page 37).

Setting the Clear Mode

To configure the Clear mode:

- 1 Select **Logger | Sample Timing**.

The logger is sampling (wrapped).

Clear Mode
☐ Do not clear logger ☒ Clear logger

Start Mode
☐ Start now ☒ Start at specified time

Stop Mode
☐ Wrap when full ☐ Stop when full ☒ Stop at specified time

Start Time
 Date: 19/05/2011 Hours: 12 Minutes: 24 ☐ AM ☒ PM

Sample Interval
 Days: 0 Hours: 00 Minutes: 00 Seconds: 10

Log Time
 Days: 2 Hours: 18 Minutes: 50

Stop Time
 Date: 22/05/2011 Hours: 07 Minutes: 14 ☒ AM ☐ PM

- 2 Under **Clear Mode**, select one of the following options:
 - **Do not clear logger** - This is the default setting if the logger memory is not yet full. Use this setting when the logger memory is not yet full and all you want to do is change Stop mode settings *without* clearing the existing samples. If the logger memory is full, or if Start mode or Sample Interval changes are required, you

must use the Clear Logger setting to make the necessary changes.

- **Clear Logger** - This is the default setting if the logger memory is full. The Clear Logger setting allows you to modify both Start and Stop mode parameters and change Sample Interval values. Choosing this setting results in the logger's memory being cleared following completion of the setup changes. You will lose data unless you transfer it first. If you are in doubt as to whether the information on the logger is still valuable, first make a copy of the data using **Logger | Transfer**, and then make the necessary logger setup changes.

When you choose **Clear logger**, the **Start Mode** and **Stop Mode** radio buttons become available.

Setting the Start Mode

The **Start Mode** enables you to choose when the logger starts taking samples.

To set the Start Mode:

- 1 Under **Clear Mode**, select the **Clear logger** option.

Logger Sample Timing

The logger is sampling (wrapped).

Clear Mode

☐ Do not clear logger ☒ Clear logger

Start Mode

☐ Start now ☒ Start at specified time

Stop Mode

☒ Wrap when full ☐ Stop when full ☐ Stop at specified time

Start Time

Date: 19/05/2011 Hours: 12 Minutes: 29 ☐ AM ☒ PM

Sample Interval

Days: 0 Hours: 00 Minutes: 00 Seconds: 10

Log Time

Days: 2 Hours: 18 Minutes: 50

Stop Time

Date: Hours: Minutes: ☐ AM ☒ PM

- 2 Under **Start Mode**, select one of the following options:
 - **Start now** - This option sets a connected logger to start taking samples immediately (within one sample interval) after the changes in the Logger Sample Timing window have been applied to the logger's memory.
 - **Start at specified time** - This option sets a connected logger to start taking samples at a pre-defined date and time in the future. Until then, the logger remains dormant. This delayed start feature helps to save logger memory and minimizes the collection of unwanted data. You can also use the **Start Time** option to ensure that multiple loggers are synchronized to start logging at the same time.
- 3 When you choose the **Start at specified time** option, the **Start Time** section of the Logger Sample Timing window becomes available. Use the **Date**, **Hours**, and **Minutes** drop-down lists and **AM/PM** radio buttons to set the time you want the logger to begin recording samples.

Note: You cannot choose a date or time earlier than the current date and time values.

If desired, use the **Preset** button to automatically set the start time hours and minutes to the earliest possible setting for the date selected. For any day other than the current day, this value is 12:00AM.

You can use the **Start at specified time** option in conjunction with the **Stop at specified time** option to create a pre-defined logger monitoring period.

Setting the Stop Mode

The **Stop Mode** enables you to choose how the connected logger operates once the memory is full, or when a specified time is reached.

To set the Stop Mode:

- Under **Stop Mode**, select one of the following options:

- **Wrap when full** - This is the default method of collecting samples on the logger. When you choose this setting, samples are recorded in a first in, first out manner. When the memory is full, the logger continues to take samples by discarding the oldest one, and replacing it with the newest one. This setting ensures that you always have the latest data when you transfer from the logger. The length of time that the memory fills and wraps around is a function of the **Sample Interval** (see “Setting the Sample Interval” on page 37) you have set on the logger. vLog automatically calculates this time and displays it in the **Log Time** section of the Logger Sample Timing window.

Note: If the memory wraps around too quickly, you can choose a longer sample interval.

- **Stop when full** - When you choose this method, the logger stops taking readings once the memory is full. This option allows the maximum amount of samples to be taken without overwriting data. It is particularly important in any application where there is an

excessive time delay between the gathering of the data and the transferring of that data to the PC.

- **Stop at specified time** - This method enables you to set a precise time for a logger to stop recording data. When you choose this setting, both the **Log Time** and **Stop Time** sections of the Logger Sample Timing window are available.

The **Log Time** drop-down lists enable you to set an elapsed time for the logger to collect samples. For example, you can set the logger to collect samples for 24 hours, 48 hours, 7 days, and so on.

The **Stop Time** drop-down lists enable you to set absolute date and time values for when the logger is to stop taking samples.

Note: The Log Time and Stop Time options are interrelated. Setting the Log Time to a specified period automatically adjusts the Stop Time values. Adjusting the Stop Time Date and Time values, in turn, automatically adjusts the Log Time values.

Setting the Sample Interval

You can change the sample interval on the logger to suit the requirements of the application. You can also set up more frequent sampling intervals to capture rapidly changing data, or spread out the readings to conserve memory usage and cover extended monitoring periods.

Note: You can only change sample intervals if you set the Clear Mode to the Clear logger setting ("Setting the Clear Mode" on page 33).

To change the sample interval on the currently-connected logger:

1 Select **Logger | Sample Timing**.

- 2 Under **Sample Interval**, use the drop-down lists to set the **Days**, **Hours**, **Minutes**, and **Seconds** to the desired sample interval. If you want to use a sample interval of once per minute, click **Default**. You can choose intervals ranging from a minimum of ten seconds to once every 24 hours (with the exception of the VL-1000-VLT, VL-1200 and models SP-1016, SP-1200, SP-1416, VL-1016, VL-1416, which have minimum sample intervals of one minute). The total period covered by a logger before the memory is full depends on how many channels are enabled and the sample interval setting. This information is automatically calculated in the **Log Time** section of the Logger Sample Timing window.

Note: Running the logger with a sample interval shorter than one minute for a prolonged time shortens the battery life.


Enabling and Disabling Logger Channels

The **Logger Channels** option permits you to enable or disable the channels on the currently-connected logger.

Note: During this process the logger's memory is cleared. If the information on the logger is still of value, transfer the data to a PC before enabling or disabling channels.

To enable or disable logger channels:

- 1 Select **Logger | Channels**.
- 2 In the window that appears, select or deselect the checkbox to enable or disable a channel.

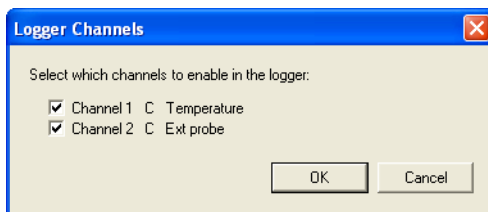

Click **OK**. This initiates the process of clearing the logger's memory. To conserve logger memory it is a good idea to disable unused channels.

Note: If you are using a temperature and humidity logger, remember that the humidity sensor requires temperature values; therefore, you cannot disable the temperature sensor channel if you are using the humidity channel.

Configuring Channels to Work with Transducers

This section applies only to those using Vaisala Veriteq VL-4000 and SP-4000 Series Data Loggers. Used with transducers, you can configure VL-4000 and SP-4000 Series loggers to display measurements in units other than milliAmps or volts.

To configure VL-4000 and SP-4000 Series Data Loggers to work with transducers:

- 1 In vLog, choose **Logger | Setup (⌘S)**.
- 2 From the Setup window, click **Setup** next to the channel to configure. Different loggers offer different channel types, each with different setup windows and slightly

different steps. In general, configure the logger's typical units and range to reflect the units and range you want to display.

3 To set up a DC Current channel:

Setup DC Current Channel

☒ Channel Enabled

Description (max 12 chars): Palm505

Input Range
20 mA

Input Signal
High (mA): 20.000
Low (mA): 0

Display Units
High: 20.000 Units: mA
Low: 0

OK Cancel Warmup Time...

- a** Ensure the Channel is enabled by selecting the **Channel Enabled** check box.
- b** If applicable, enter a new channel description.
- c** From **Input Signal**, specify the range for your transducer output.
- d** Under **Display Units**, specify the range and units to display instead of the input signal range. For example, psi instead of mA.
- e** To set warmup time and save transducer battery life, see “Configuring Warm Up Time” on page 42.

4 To set up a Voltage channel:

- a Ensure the Channel is enabled by selecting the **Channel Enabled** check box.
- b If applicable, enter a new channel description.
- c If applicable, for **Input Range**, choose whether this is a 1 or 10 volt range channel.
- d From **Input Signal**, specify the range for your transducer output.
- e Under **Display Units**, specify the range and units to display instead of the input signal range. For example, 'psi' instead of 'v'.

vLog also supports a Boolean DC Voltage Channel type on some SP- and VL-4000 loggers. For this type, specify the Threshold Voltage and the resulting display units. For example, if you want voltage above the threshold value to show as 1, indicate 1 for Above Threshold and specify

the units). Boolean choices are typically 1 (above threshold) and 0 (below threshold).

- 5 Click **OK** to save your changes.

Configuring Warm Up Time

Note: This section only applies to Vaisala Veriteq VL-4000 and SP-4000 Series Data Loggers.


Some VL-4000 and SP-4000 loggers have the ability to switch power from an external battery onto a transducer using the Warmup Time feature, also known as the XPS excitation control switch. This is useful for conserving power when an external transducer is being powered from an external battery.

Set transducer warm up time to a number of minutes or seconds before data is recorded by the logger. To set when the logger will record data, see “Setting Logger Sample Timing” on page 33.

To configure transducer warmup time:

- 1 Choose **Logger | Setup (X)**
- 2 From the row for the channel you want to set up the warm up time for, click **Setup**.
- 3 From the Setup Channel window, click **Warmup Time**.
- 4 From the Warmup Time window, select the **Timer** radio button. You can also set up warmup time from the Logger Setup window.

- 5 For Timer Value, set the number of minutes or seconds you want the transducer to begin warming up before taking a sample. The default is 10 seconds.

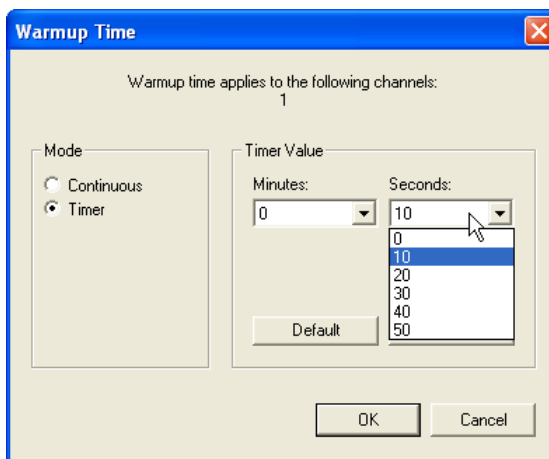

- 6 Click **OK**.


Configuring Logger Channels to Work with Thermocouples

This section applies only to those using Vaisala Veriteq VL-1700 and SP-1700 Series Data Loggers. These loggers work with various types of thermocouples (E, J, K, R, S, and T). You need to configure your channel to specify the type of thermocouple you are using.

To configure VL-1700 and SP-1700 Series Data Loggers to work with thermocouples:

- 1 In vLog, choose **Logger | Setup** (⌘S).
- 2 From the Setup window, click **Setup** next to the channel to configure.

- 3 From the Setup Thermocouple Channel screen, select **Channel Enabled**.

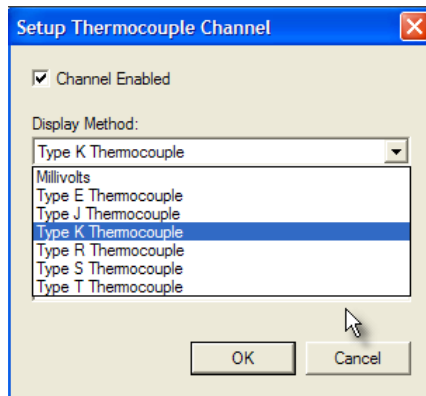

- 4 From the Display Method list, choose your thermocouple type.
- 5 Enter a description.
- 6 Advanced: If your installation includes very long thermocouple cables, Vaisala Canada Inc may recommend you specify a temperature offset. Enter it here. All values will be adjusted by this offset.
- 7 Click **OK** to save or **Cancel** to close without saving.

Clearing Logger Samples

Clearing a logger deletes all recorded samples in a logger's memory.


You can clear the logger without making any setup changes, or as an essential part of the following logger setup functions:

- enabling or disabling channels
- changing **Start Mode** settings
- changing logger sample intervals

Note: You do not have to clear the logger samples before you start a new logging session, but doing so minimizes the storage of unneeded information.

To clear the logger without changing the logger setup:


- 1 Select **Logger | Clear**.

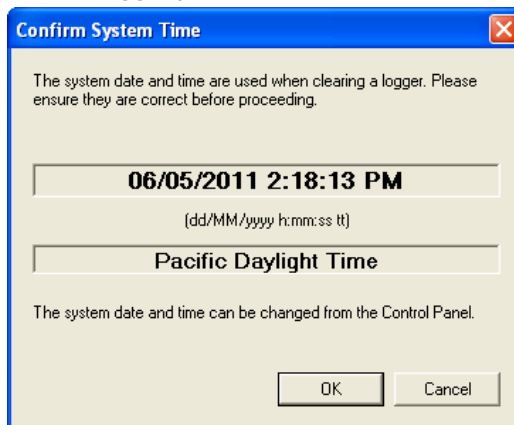

- 2 In the Confirm System Time window, check the system time.
 - a If the time shown is correct, click **OK**.
 - b If the time shown is incorrect, click **Cancel**, and then use the Windows Control Panel to make the appropriate changes. See "Setting System Date and Time" on page 64.

Note: The PC's date and time must be accurate when clearing a logger's memory or the date and time values in future logger files may be incorrect.

After you click **OK**, a confirmation window opens:

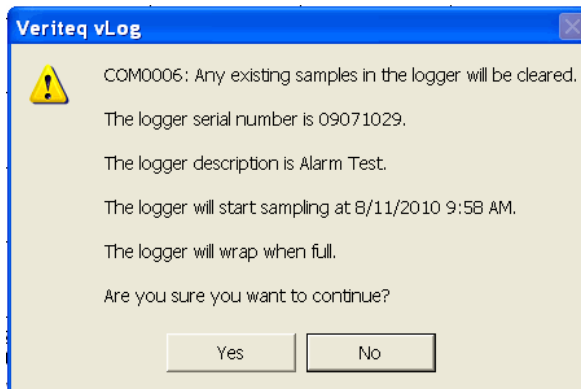

- 3 Click **OK**.

Configuring Loggers: Batch Setup (vLog VL)

Before using this section, read “Configuring Loggers” on page 26, as you need to have an understanding of the settings and choices you will make before using Batch Setup.

Batch Setup, available in vLog VL, allows you to quickly configure several loggers with very similar settings. For example, using Batch Setup, you can quickly apply identical logger sample settings to all loggers in your environment.

To set up the first logger in a batch:

- 1 Ensure a logger is connected to your PC.
- 2 Choose **Logger | Batch Setup**.

Logger Batch Setup

Communication Settings

Select the first COM port: COM1 Refresh

Select the number of COM ports: 6

Logger Sample Timing Settings

Start Mode: Start Now Modify...

Sample Interval: 1 minute

Stop Mode: Wrap when full

Logger Description Settings

Modify Description (max 16 chars):

☒ Logger

Channel Settings

	Enable Channel	Modify Description	Description (max 12 chars):
Channel 1:	☒	☒	Temperature
Channel 2:	☒	☒	Humidity
Channel 3:	☐	☐	Unmodified
Channel 4:	☐	☐	Unmodified
Channel 5:	☐	☐	Unmodified

View Report Setup Logger Close

- 3 Optional: Click **Refresh** to update the COM port list with available ports.
- 4 From the Logger Batch Setup window, choose the COM port your logger is communicating with and specify the number of COM ports in the range of your batch setup. (For example: If setting up loggers on COM ports 2, 4, and 6, the number of COM ports selected should read 5, not 3.)

- 5 To set logger sample timing settings, click **Modify** in the Logger Sample Timing Settings area. In the window that appears, set start mode, stop mode, sample interval, and so on. Click **OK**.
- 6 To set logger description, select the **Modify** checkbox and enter a new description (depending on your regional font settings, the space available for the description may be less than 16 characters).
- 7 If desired, enable or disable channels.
- 8 To modify channel descriptions, select the **Modify Description** checkbox for the channel to modify and enter a new description.
- 9 Click **Setup Logger** (you can click **Close** to close without saving).
- 10 If security is enabled, enter your password, then click **OK**.
- 11 When the Confirm System Time dialog appears confirm the exact time before clearing the logger. Click **OK**.
- 12 Confirm any changes to the description(s). Click **OK**.
- 13 A message appears showing the changes to be made, asking if you want to continue. Click **Yes**.

Second and subsequent loggers if performing Batch Setup using only one COM Port:

- 1 With the Logger Batch Setup window still open, disconnect the current logger from the logger cable and attach the next logger to be configured.
- 2 On the Logger Batch Setup window, change the logger or channel descriptions if desired.
- 3 Click **Setup Logger**.
- 4 Confirm the changes to the logger and descriptions. Click **OK**.
- 5 A message appears showing the changes to be made, asking if you want to continue. Click **Yes**.
- 6 The settings as specified will be applied to this logger.
- 7 Repeat steps 1 through 6 for all remaining loggers.
- 8 When all loggers have been setup, click **View Report** if desired, otherwise, click **Close**.

Second and subsequent loggers if performing Batch Setup using multiple COM ports (vLog VL):

- 1** The Confirm Description dialog box opens for the next COM port.
- 2** If this COM port is assigned to a logger you wish to setup with the same Sample Timing settings as the first logger, skip to step **5**. If this COM port is not assigned to a logger you wish to setup, click **Cancel**.
- 3** A message appears, asking if you would like to abort the batch setup. Click **No**.
- 4** Repeat from step 1 until you reach a COM port which is assigned to a logger you wish to setup.
- 5** In the Confirm Descriptions dialog box, change the Logger and Channel Descriptions as desired. Click **OK**.
- 6** A message appears showing the changes to be made, asking if you want to continue. Click **Yes**.
- 7** Repeat from step 1 for all remaining loggers.
- 8** When all loggers have been setup, click **View Report** if desired, otherwise, click **Close**.

Mounting Loggers

Because they are small and light, you can easily mount Vaisala Veriteq loggers in the field using a variety of techniques including Velcro strips, double-sided tape, duct tape, and magnetic strips.

In harsh environments, place the logger in a standard industrially-rated enclosure, and run external sensors and probes from the enclosure using special fittings in the openings.

Some tips for placing loggers:

- **Ventilation** - Do not obstruct the air passage holes in the logger. There should be a clear path for the air to get to the internal sensors (if they are used or available).
- **Vibration** - Excessive vibration can knock the unit off its mounting. If the location is subject to excessive jarring or vibration, mount the logger more securely.
- **Temperature** - Excessive temperature (low or high) can affect adhesives.

- **Moisture** - Moisture can also affect the adhesive durability.
- **Permanence** - Make sure the mounting location is permanent for the duration of the recording period.
- **Public access** - If the logger location is too visible, people may interfere with the readings. Make sure it is not too accessible.

Using External Temperature Probes

Some Vaisala Veriteq loggers have external temperature channel(s) for use with an external temperature probe. Using the external temperature probe enables you to collect data from hard-to-get-at locations, such as in duct work or under pipe insulation.

You can use an external temperature probe with these loggers, provided the probe was calibrated with the logger.



To use an external temperature probe:

- 1 Enable the desired channel on the logger (see “Enabling and Disabling Logger Channels” on page 39). If you do not need to monitor from the internal temperature sensor, disable it to conserve memory.

Note: You can disable most unused channels - external probe or not - to conserve logger memory. One exception is the internal temperature channel on VL-1700 and SP-1700 Data Loggers.

- 2 Attach the temperature probe to the connector on the logger.



- 3 Use the **Logger | Setup** window to test the logger operation (refresh the readings as required).
- 4 Place the logger in the desired location and start recording.

Caring for Relative Humidity Sensors

Each Vaisala Veriteq SP- and VL-2000 Data Logger has a relative humidity sensor. Although the sensor has been designed to resist contamination and condensation, avoid exposing the logger to contaminants, if possible.

For information on specific environments, contact Vaisala Canada Inc.

Chapter 3: Transferring Logger Files

This section provides details about:

- transferring logger data and setting transfer preferences
- generating a logger file creation report and ID number (vLog VL)
- setting the system date and time
- changing the time zone
- transferring logger data in batches to your PC (vLog VL)

Transferring Data

When you transfer data from a logger's memory, vLog creates a copy of the data and stores it on your PC as a logger file.

Logger transfers can be performed either individually (logger by logger) or in a batch (by selecting a range of COM ports using vLog VL).

Data from secure loggers (VL loggers) is stored in secure logger files. Data from non-secure loggers (SP loggers) is stored in non-secure logger files. All logger files have the .spl extension, logger files from SP loggers are non-secure and can be opened in vLog SP or vLog VL. Files from VL loggers are secure and can only be opened in vLog VL.

After inserting a logger file into vLog, you can view the information as a graph or as a tabular report.

VL Data Loggers (vLog VL)

During the data transfer process, vLog generates a Logger File Creation Report and a unique Logger File ID Number for VL-Series loggers. The Logger File ID Number also appears in the Logger file and on the graph and provides a convenient way to track the data from the logger to the final Graph file. The transfer process does not affect or delete the data in the logger.

Before transferring data, you may want to set the transfer preferences, and ensure that you have write permission to the folder on the PC where you are going to save your transferred data as logger files.

Setting Transfer Preferences

There are a variety of logger configuration preferences available to you:

- Configure the default file name format for logger files to help manage file storage.
- Configure vLog VL to display a Logger File Creation Report after the transfer. This report displays all the details of the transfer. It can be printed and it provides a secure and verifiable link between logger files and printed

records. This report is only produced for individual logger transfers, not for batch transfers.

- Configure vLog VL to interrupt batch transfer operations when an error is detected. If this option is not selected, vLog will skip over any errors that are detected and complete the batch transfer.

Note: In vLog 4.10, the Logger File Creation Report feature was added to eliminate the need for using usernames and passwords to verify logger files. If you want to continue to configure vLog to prompt you for a username and password when you transfer data, see page 106.

Configuring the Default File Name Format

To configure the default file name format:

- 1 Select **Tools | Options**, then choose the **Transfer** tab.

The screenshot shows the 'Options' dialog box with the 'Transfer' tab selected. The dialog has a title bar with a close button. Inside, there are three tabs: 'Historical Data Report Statistics', 'Audit Trail Report', and 'Batch Setup Report'. The 'Transfer' tab is active, showing sub-tabs 'General', 'Graph', 'Transfer', 'Insert', and 'Historical Data Report'. The 'Transfer' sub-tab is selected. The main area contains the text: 'When transferring samples, use the following format for the default name of the logger file:'. Below this are three columns for 'Field 1', 'Field 2', and 'Field 3'. Each column has radio buttons for 'None', 'Serial Number', 'Description', 'Date and time', 'Date', and 'Time'. 'Serial Number' is selected for Field 1, 'Description' for Field 2, and 'None' for Field 3. A '+' sign is between Field 1 and Field 2, and another '+' is between Field 2 and Field 3. Below the fields is a 'Sample:' label and a text box containing '99999999_Description.spl'. At the bottom, there are three checkboxes: 'Prompt for a username and password to store in the logger file when transferring samples. Note that this feature is optional. Vaisala now uses Logger File ID Numbers to provide a secure and verifiable link between logger files and printed records.' (unchecked), 'Display a Logger File Creation Report after transferring the samples from a logger.' (checked), and 'Interrupt batch transfer operations when an error is detected.' (unchecked). There is a 'Restore Defaults' button and 'OK' and 'Cancel' buttons at the bottom right.

- 2** In each of the three Field columns, select the information you want to include in the file name. The options for each field are:

Option	Format
None	This field is not used in the Logger file name format
Serial number	Eight digit logger serial number
Description	Description of logger
Date and time	yyyy-MM-dd_hh-mm, where: yyyy is the year MM is the month dd is the day hh is the hour mm is the minutes
Date	yyyy-MM-dd, where: yyyy is the year MM is the month dd is the day
Time	hh-mm, where: hh is the hour mm is the minutes


- 3** Click **OK**.

Transferring Logger Data to the PC

To transfer the information stored in one data logger to the PC:

- 1** Connect the data logger to the PC (see page 16).
- 2** Ensure that the COM port is configured (see page 25).
- 3** Ensure that you have write permission to the folder on the PC where you are going to save the transferred data as a logger file.

- 4 Select **Logger | Transfer** ( or F12). This starts the communications process between the logger and the PC and opens the Transfer Wizard window.

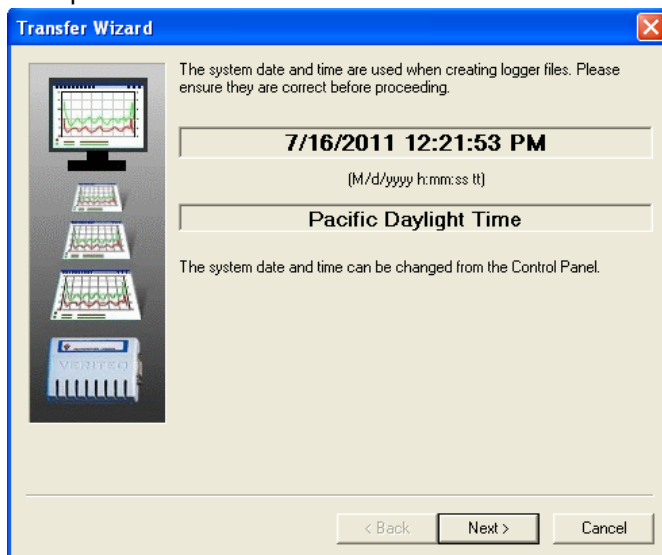

- 5 Confirm that the current date and time shown in the window are correct. If the date and time are not correct, use the Windows Control Panel to update the date and time. See “Setting System Date and Time” on page 64.

Note: The PC's date and time must be correct when setting up a data logger, or transferring logger data. Otherwise, the timestamps in the Logger file may be incorrect.


6 Click **Next**.

- 7 Edit or confirm the description of the logger using a maximum of 64 characters. The default entry for this field uses the existing **Description** in the logger's memory (see "Configuring Loggers" on page 26). You can accept the default logger description, or use the extra character spaces to further describe the application or location of the data logger.


Note: Changing the description here does not modify the description in the logger's memory.

- 8 Edit or confirm the description of listed logger channels.
9 Click **Next**.

- 10 Accept the default name for the Logger file or click **Browse** to save the Logger file under a different name or folder location.

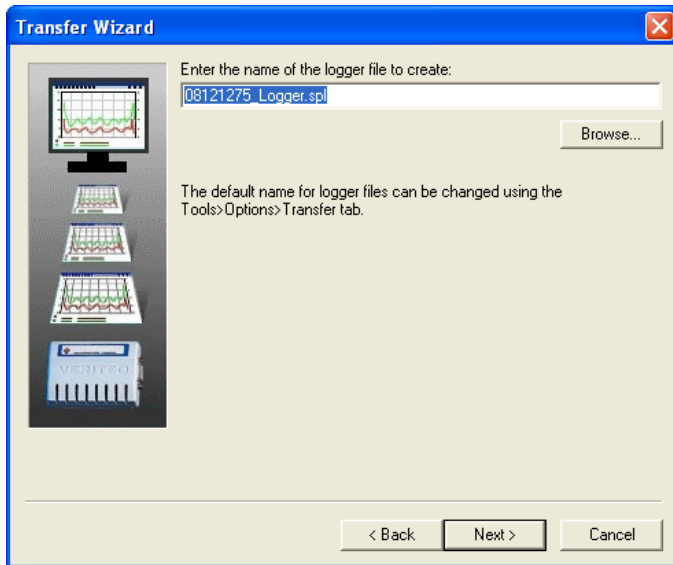


- 11 If you use vLog SP, continue to step 14.
- 12 If vLog VL is configured to prompt for a username and password (Tools | Options | Transfer) when transferring files, click **Next** and the following window opens.

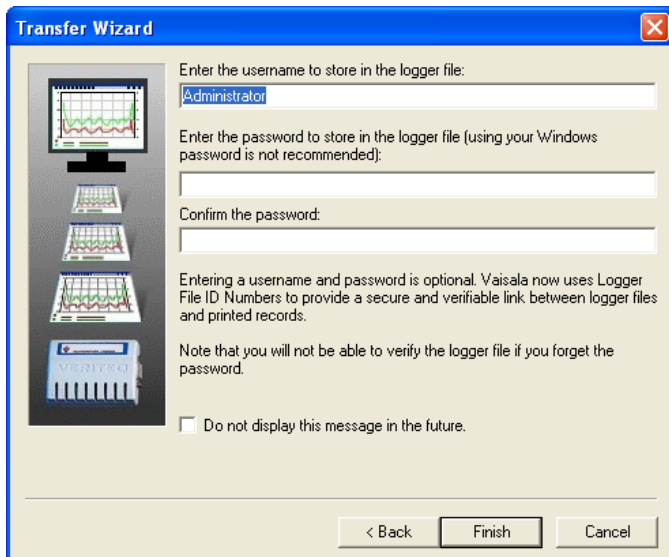


- 13 Type a username and password. Retype the password to confirm.

CAUTION: If you use this feature, it is vital to remember the username and password you enter here, as you will need to re-enter the password when you verify the graph (see page 82). **This information cannot be recovered.**

- 14 Click **Finish**. The logger transfers its data to the PC while displaying the progress of the transfer in the Transferring Samples window.

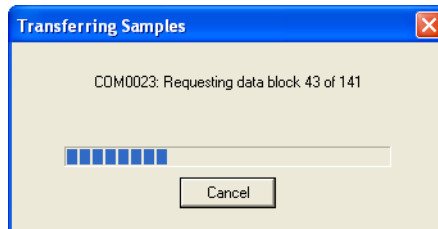

- 15 If you use vLog VL, once all samples are transferred, the Logger File Creation Report opens.

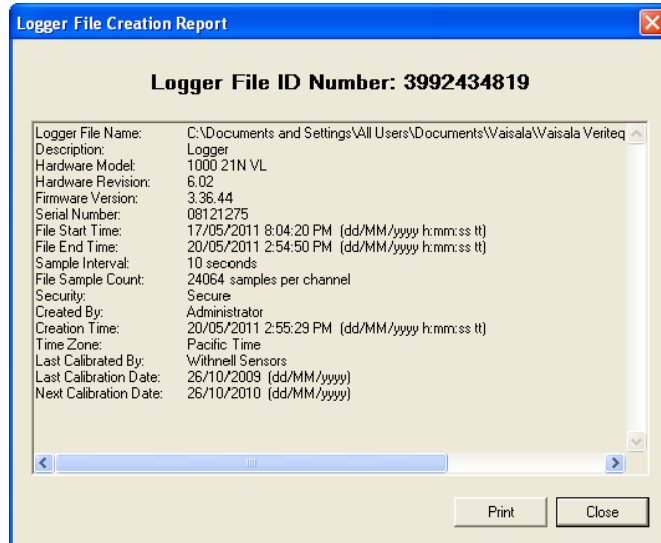


If you need a copy of this report, you can print it now (the information in this report is available in the Logger file).

To print a copy of the report, click **Print**. To close this window and return to the vLog main window, click **Close**.

Note: After the transfer, the data remains in the logger. If the logger memory was not full, the data logger continues to store additional data in its internal memory without interruption. If you want to clear the data logger, see “Clearing Logger Samples” on page 44.

- 16 Insert the newly transferred logger file into a new graph or the active (current) graph, if one is active.
- 17 Click **OK**.
- 18 From the Select Channels window that appears, confirm the channels you want to display in the graph.

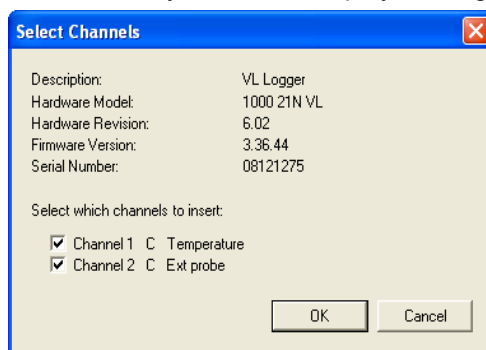


- 19 Click **OK**.
- 20 If one of the channels you want to display records relative humidity, the following window opens. Confirm the display method for the relative humidity channel.



- 21 Click **OK**. The newly transferred file will appear in the active graph or a new graph based on your choice in step 16.

- 22** If one of the channels you want to display is a thermocouple channel, the following window opens.

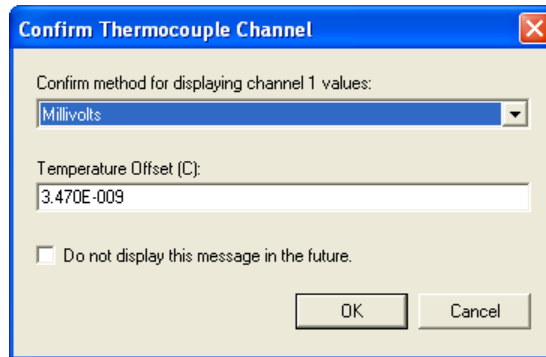

- 23** Confirm the display method and temperature offset. If your installation includes very long thermocouple cables, Vaisala Canada Inc may recommend you specify a temperature offset (all values will be adjusted by this offset).
- 24** Click **OK**. The newly transferred file will appear in the active graph or a new graph based on your choice in step **16**.


Logger File Creation Report (vLog VL)

vLog VL software can be configured to generate and display the optional Logger File Creation Report after a VL logger's data is transferred to a PC through an individual transfer. This report is not created for loggers transferred in a batch transfer.

This report provides a detailed record of each Logger file you create. Each report contains a unique, automatically-generated Logger File ID Number. This number provides a means of positively identifying individual logger files when reports or graphs are printed for review and approval.

In addition to the Logger File ID Number, the Logger File Creation Report contains the following information necessary as part of the validation process:

- Logger file name
- Logger description
- Hardware model of the data logger
- Hardware revision for the data logger
- Firmware version number for the data logger
- Serial number of the data logger
- File start time
- File stop time
- Sample interval
- File sample count
- Whether or not the file is secure
- Username of the person who created the file (if that option is selected)
- Creation time
- Time zone
- Name of the organization who last calibrated the data logger
- Date of the last calibration
- Recommended date for the next calibration

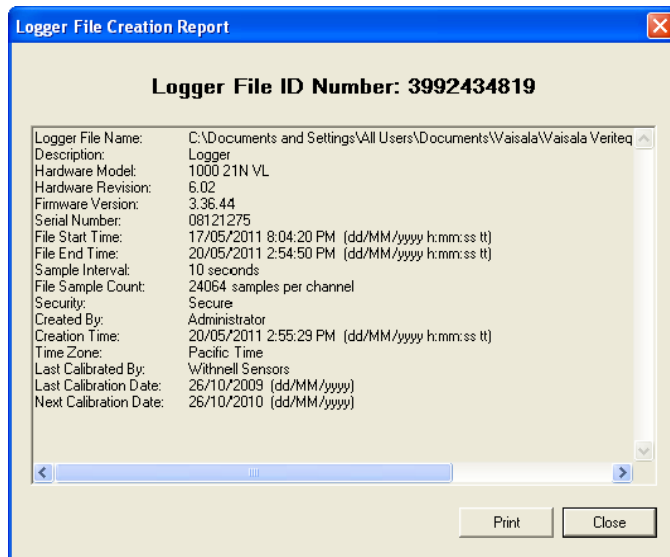

To configure vLog to generate a Logger File Creation Report after data transfers:

- 1 Select **Tools | Options**, then choose the **Transfer** tab.
- 2 Select the option, **Display a Logger File Creation Report after transferring the samples from a data logger**.

Options

Historical Data Report Statistics | Audit Trail Report | Batch Setup Report

General | Graph | **Transfer** | Insert | Historical Data Report

When transferring samples, use the following format for the default name of the logger file:

Field 1: ☐ None ☒ Serial number ☐ Description ☐ Date and time ☐ Date ☐ Time

Field 2: ☐ None ☐ Serial number ☒ Description ☐ Date and time ☐ Date ☐ Time

Field 3: ☒ None ☐ Serial number ☐ Description ☐ Date and time ☐ Date ☐ Time

Sample: 99999999_Description.spl

☒ Prompt for a username and password to store in the logger file when transferring samples. Note that this feature is optional. Veriteq now uses Logger File ID Numbers to provide a secure and verifiable link between logger files and printed records.

☒ Display a Logger File Creation Report after transferring the samples from a data logger.

☐ Interrupt batch transfer operations when an error is detected.

Restore Defaults

OK Cancel

3 Click **OK**.

Logger File ID Number (vLog VL)

The Logger File ID Number is a unique checksum value that is automatically generated every time a VL-logger file is created. By recording the Logger File ID Number as shown in the Logger File Creation Report, you can ensure this ID number appears on all reports or graphs you have to sign as part of regulatory processes.

The Logger File ID Numbers are included in the following files and documents:

- Logger File Creation Reports
- Logger files
- Graph files that contain inserted VL logger file data.
- Reports

Setting System Date and Time


To set the PC's date and time using the Windows operating system:

- 1 Select **Start | Control Panel**, and then double-click **Date/Time**.
- 2 Reset the time to the correct values.
- 3 Click **OK**.

Batch Transfer of Logger Data to the PC (vLog VL)

To transfer the information stored in multiple loggers to the PC in a batch:

- 1 Connect the data loggers to the PC (see page 16).
- 2 Ensure that the COM ports are configured (see page 25).
- 3 Ensure that you have write permission to the folder on the PC where you are going to save the transferred files.
- 4 Select **Logger | Batch Transfer**. The Batch Transfer Wizard window opens.

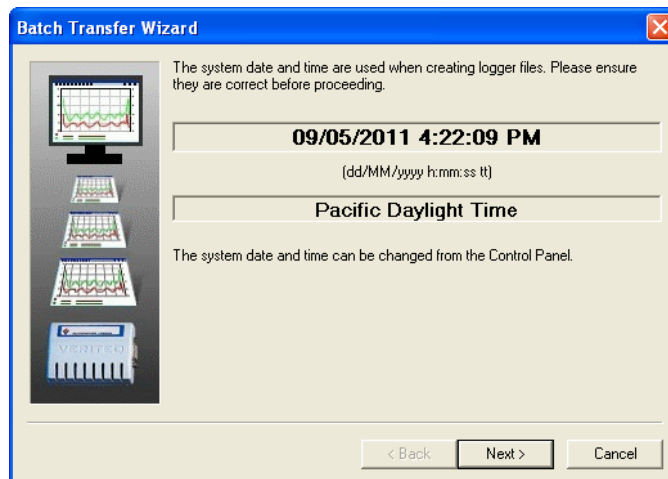


- 5 Confirm that the current date and time shown in the window are correct. If the date and time are not correct, use the Windows Control Panel to update the date and time (see page 64).

Note: The PC's date and time must be correct when setting up a data logger, or transferring logger data. Otherwise, the timestamps in the logger files may be incorrect.


6 Click **Next**.

7 Select the first COM port in your batch, then select the number of COM ports in the range of your batch transfer. (Example: If transferring data from loggers on COM ports 2, 4, and 6, the number of COM ports selected should read 5, not 3.)

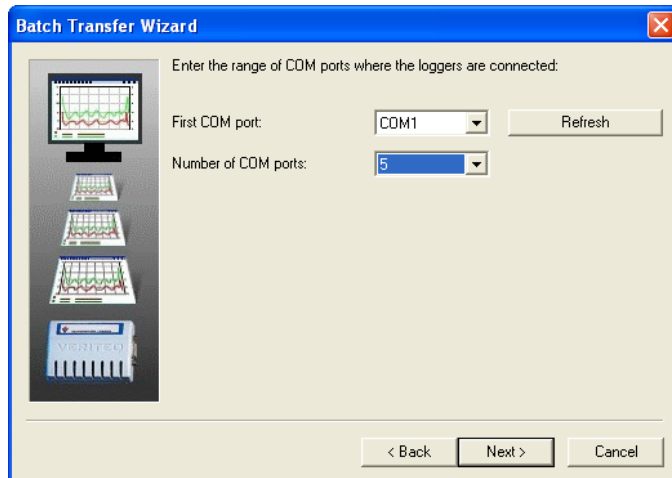
The screenshot shows a Windows-style dialog box titled "Batch Transfer Wizard" with a red close button in the top right corner. On the left side, there is a vertical stack of four icons: a computer monitor displaying a line graph, a data logger device, a multi-channel terminal block, and a Vaisala Veriteq data logger unit. The main area of the dialog has a light beige background and contains the text "Enter the range of COM ports where the loggers are connected:". Below this text are two input fields: "First COM port:" with a dropdown menu showing "COM1" and a "Refresh" button to its right; and "Number of COM ports:" with a dropdown menu showing "5". At the bottom of the dialog, there are three buttons: "< Back", "Next >", and "Cancel".

8 Click **Next**.

- 9 Accept the default folder where the Logger files will be created, or click **Browse** to select a different location.

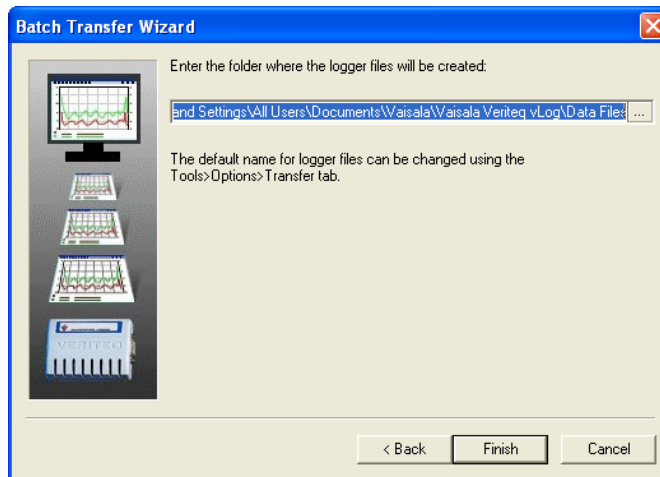

- 10 Click **Finish**.
- 11 vLog transfers the logger data to the PC while displaying the progress of the transfer. After all data is transferred, a dialog box appears prompting you to insert the new logger files into a new graph or the active (current) graph. Choose the appropriate radio button.

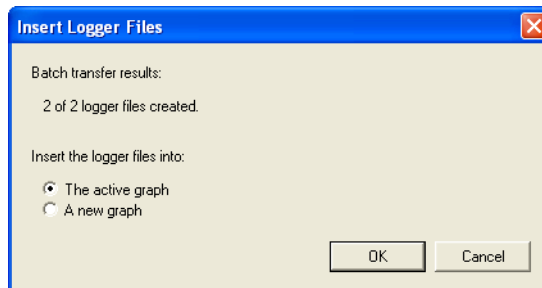

- 12 Click **OK**.

Chapter 4: Working with Graphs

Graphs are a key component of vLog, allowing you to display and analyze data.

This section provides details about:

- displaying existing and creating new graph files
- formatting graphs: adding graph titles, modifying fonts, formatting time zones and measurement units
- working in graphs: zooming in or out, changing the temperature scale, sample smoothing, displaying line numbers and horizontal lines
- moving graph data to new locations: copying and pasting graphs, exporting graph data and saving and printing graphs

Displaying and Creating Graph Files

To open an existing Graph file (.spg):

- 1 Select **File | Open** (📂 or Ctrl+O).
- 2 Navigate to the Graph file you want to open.
- 3 Click **Open**.

The graph opens in the vLog window.

To create a new graph:

- ▶ Select **File | New** (📄 or Ctrl+N).

An empty graph window opens, ready for you to insert one or more Logger files. The new Graph has a default name (for example, Graph1). You can rename the graph when you save it.


Inserting Logger Files

To insert a Logger file into a graph:

- 1 Select **File | Insert Logger Files** (📁).

Note: If you have logger files created using sGo that you want to insert, select **File | Insert Logger Files from sGo Folder** to automatically open the correct folder (Veriteq sGo software for Palm OS® allows users to setup, read, and download SP-Series Data Loggers using a Palm™ hand-held device, instead of a laptop computer).

- 2 Select the Logger file you want to insert.
- 3 Click **Open**. The Select Channels window opens.

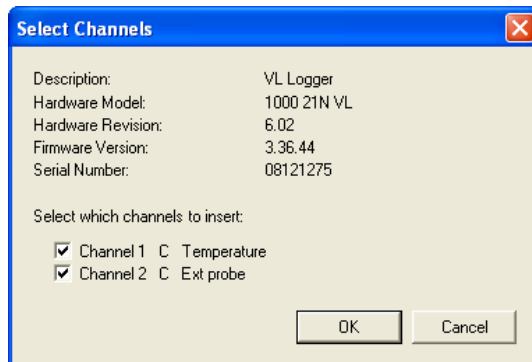


- 4 Select the channels you want to display in the graph.
- 5 Click **OK**.
- 6 If you are inserting a file with temperature and humidity, and the option for confirming display method for relative humidity has been set to yes, the Confirm Relative Humidity Channel window opens, prompting you to choose a display method for relative humidity.



Confirm the method for displaying values, then click **OK**.

If you do not want to see this window every time you work with humidity logger files, do one of the following:

- Select the **Do not display this message again** check box before clicking **OK**, or,
 - Select **Tools | Options**, choose the **Insert** tab, and clear the **Confirm display method for relative humidity channels** check box.
- 7 If you are inserting a file with a thermocouple channel, confirm the method to display the channel, and, if required, the temperature offset:

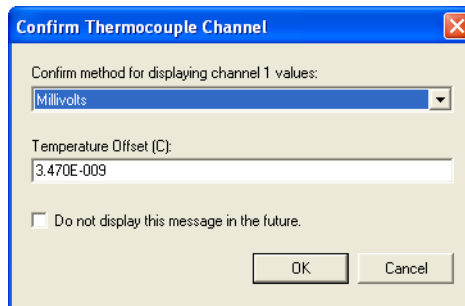

The logger file is displayed in the graph. If the logger file has more than one channel, vLog automatically selects distinct colors for each graph line.

Note: VL loggers: Check the Logger File ID Number displayed on the graph to ensure that the correct file is displayed.

Inserting Multiple Logger Files

You can insert multiple logger files in one graph. This works best if the range and dates measured in the logger files are similar. This enables you to easily compare data, and is especially useful for before-and-after studies or to review data collected from different loggers.

To insert additional Logger files into a Graph, follow the instructions in “Inserting Logger Files” on page 68 for each logger file you want to insert.

Note: If you are comparing SP-series and VL-series logger files, the graph will indicate Not Secure in the top left corner. This means that at least one of the logger files was generated from a logger that uses a non-secure format. To evaluate the security issue you can generate a Graph Security Problems Report (see page 100).

Removing Channels from a Graph


After viewing a Graph that contains one or more Logger files -- or a single Logger file with multiple channels -- you may want to remove a channel and its associated graph lines for a cleaner view.

To remove a channel from the active graph:

- 1 From the Channel Information table, highlight the channel to remove and do one of the following:
 - From the menu, select **Edit | Remove Channel**.
 - Right-click the channel in the Channel Information table, and select **Remove Channel**.

LN	Serial #	Description	CH	Value	Maximum	Av
1	08121275	VL Logger			21.27	
2	08121275	VL Logger			21.50	
3	11012038	Logger	1		21.76	
4	11012038	Logger	2		21.76	

- 2 A confirmation window opens. Click **Yes** to confirm.

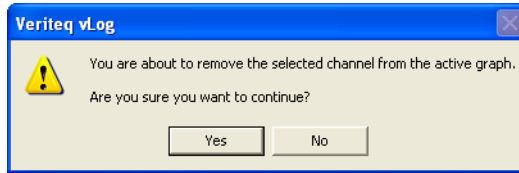

The channel and its associated graph line are removed from the active Graph.

Scrolling Through Graphs

Scrolling allows you to shift the view of an on-screen graph to the right or to the left. Each time you scroll, the graph view shifts by one graph division.

To scroll the on-screen graph to the left:

- ▶ Select **Format | Scroll Backward** (◀, or Alt + left arrow).


To scroll the graph to the right:

- ▶ Select **Format | Scroll Forward** (▶, or Alt + right arrow).

Pinpointing Graph Values and Times

To pin-point exact values and times on the currently-active graph:

- 1 Select **Format | Position Cursor** (⌘ or F7).
- 2 Position the cursor line along the X axis using the mouse or the arrow keys. The values at the bottom left portion of the window indicate the exact date and time where the cursor line is positioned as well as the exact values that correspond to that point.

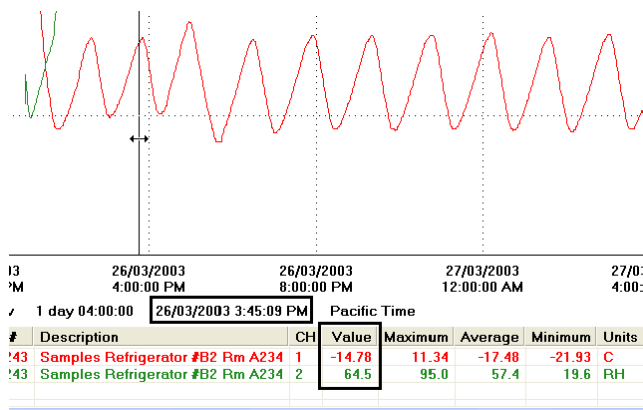

- 3 To lock the cursor in place, click the mouse or press **Enter**.

To remove the vertical cursor bar:

- Select **Format | Remove Cursor** (🗑️ or F6).

Changing the Time Zone

If you are sharing information with others around the globe, you may wish to view a graph according to a specific time zone.

To change the time zone for a graph:

- 1 Select **Format | Time Zone**.
- 2 In the Format Time Zone window, use the drop-down list to select the time zone for which you want to display your data.

Adding Graph Titles

To add a descriptive title to the active graph:

- 1 Select **Format | Title** The Title window opens.

Format Title

Enter the title for the graph (max 50 chars):

Samples Refrigerator B2 Rm A34 Temp /RH

OK Cancel


- 2 Type the desired title (up to 50 characters long).
- 3 Click **OK**.

Selecting a Graph Font

You can select any font that is on the PC to use for all displayed and printed graphs.

To change the graph font:

- 1 Select **Tools | Options**, then choose the **Graph** tab.
- 2 Click **Graph Font**.

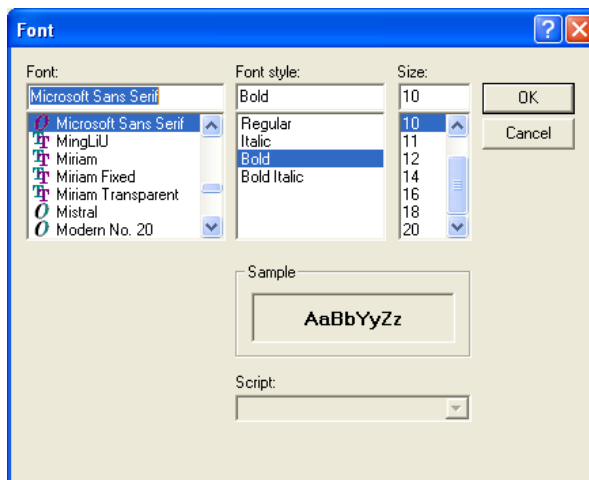

- 3 Select the desired **Font**, **Font style**, and **Size**.
- 4 Click **OK**.
- 5 If the font style and size are inappropriate, choose **Restore Defaults**.

Zooming In

Zooming allows you to look closer at a portion of the currently-displayed graph. It also enables you to eliminate extraneous data and produce a graph and a tabular printout that cover a set time period, such as 24 hours.


You can zoom-in two ways:

- Use the zoom selection tool. This is the quickest and easiest way to zoom.


- Use the **Format | Time Scale** option. This gives you precise control over the X-axis time scale of the zoom. You can zoom-in by entering exact start and end times and dates.

To zoom in on the graph using the zoom selection box:

- 1 Select **Format | Zoom-in** ( or F9). A vertical line appears on the graph.

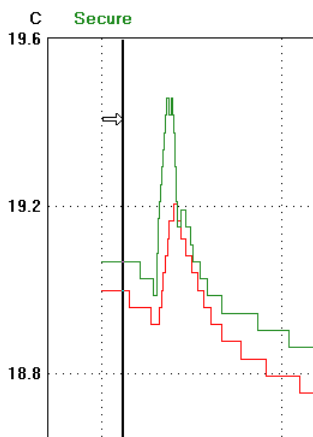

- 2 Use the mouse or arrow keys to move the vertical line to the beginning of the section you want to zoom-in on..

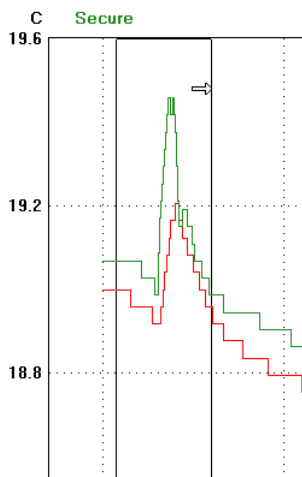


3 Either:

- a Press and hold the left mouse button, and move the second cursor line to the end of the section you want to zoom in on.
- b Release the mouse button and the zoomed-in area is displayed.

Or:

- a Press **Enter** to anchor the vertical line.
- b Use the arrow keys to stretch the zoom selection box.
- c Press **Enter** and the zoomed in area is displayed.

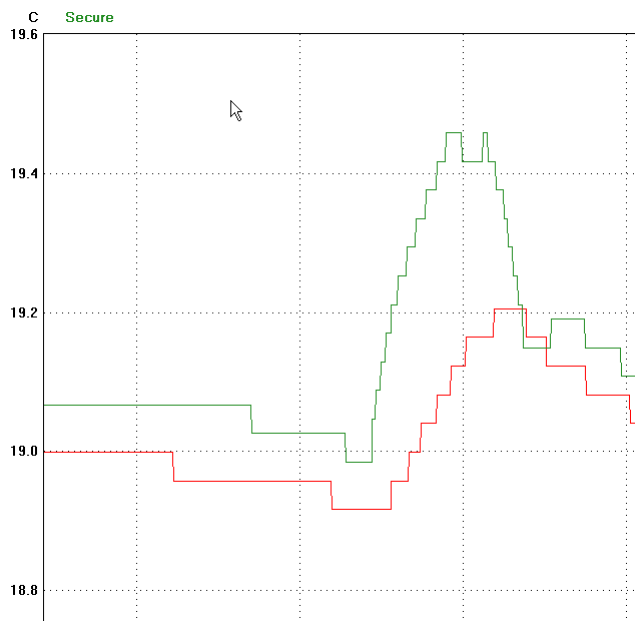


Note: If you have difficulty stretching the Zoom box accurately over the data, try turning off the Time Scale rounding setting. Select **Tools | Options**, choose the **General** tab, then clear the **Automatically round time scale after zooming in** check box.

To zoom in by formatting the time scale:

- 1 Select **Format | Time Scale**, or place the cursor over the date and time at the bottom of the graph, and right-click to select **Format | Time Scale**.

The Format Time Scale window opens.

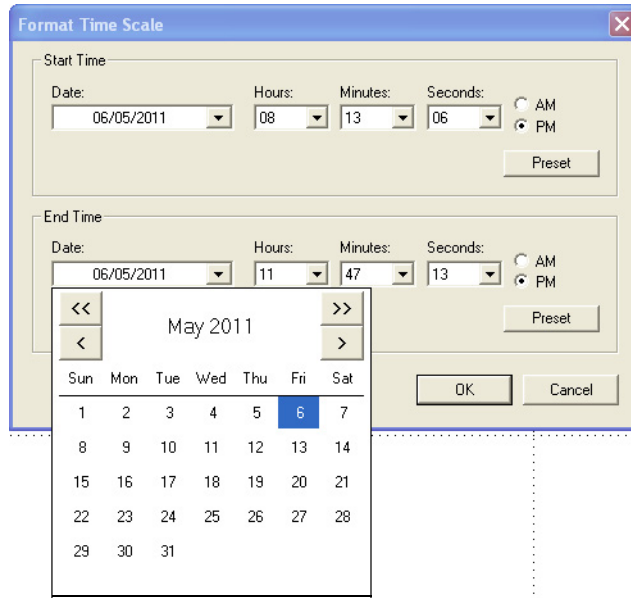

- 2 Under **Start Time**, use the drop-down lists to select the start date and time that correspond to the beginning of the section of the graph you want to view in greater detail.
- 3 Under **End Time**, use the drop-down lists to select the end date and time that correspond to the end of the section of the graph you want to view in greater detail.
- 4 You can also use the **Preset** buttons to automatically set the time to the beginning of the Calendar day chosen (12:00:00 AM).
- 5 Click **OK**.

Zooming Out

To zoom-out:

- ▶ Select **Format | Zoom-Out** ( or F10).

You can keep zooming out until you have returned the graph to its original condition.


Formatting Measurement Units

vLog automatically determines Y-axis scales appropriate for any displayed graph and sets a default number of decimal places. If desired, you can override these settings.

Specifying the Y-axis Scale

To change a Y-axis scale:

- 1 With the Graph you want to modify open, select **Format | Units**. The Format Units window opens.

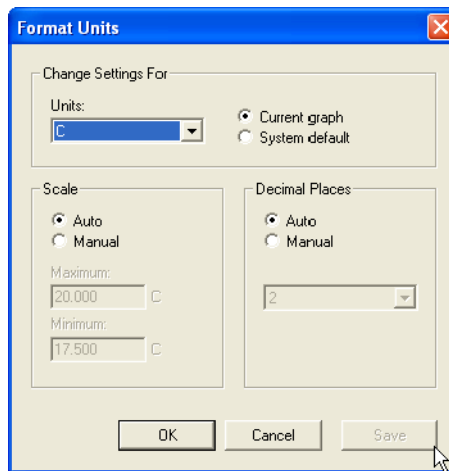


- 2 Under **Units**, use the drop-down list to select the desired units of measurement.
- 3 Under **Change Settings For**, use the drop-down list to select where you want the changes to apply—the **Current Graph** or **System Default**. If you select **System Default**, the settings entered here automatically apply to other new graphs that use the same units.
- 4 Under **Scale** select **Manual**, then enter the maximum and minimum values for the axis that you want displayed.
- 5 Click **Save** to save changes, and **OK** to close.

Changing the Decimal Place Precision

To change decimal place precision:

- 1 With the Graph file you want to modify open, select **Format | Units**. The Format Units window opens.

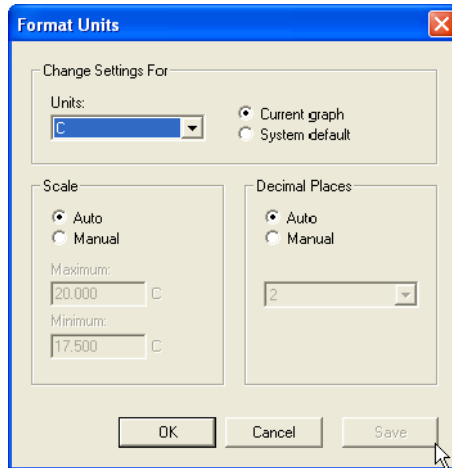

- 2 Under **Units**, use the drop-down list to select the measurement units you want to change.
- 3 Under **Change Settings For**, use the drop-down list to select where you want the changes to apply—the **Current Graph** or **System Default**. If you select **System Default**, the settings entered here automatically apply to other new graphs that use the same units.
- 4 Under **Decimal Places** select the **Manual** option and use the drop-down list to select the desired number of decimal places.
- 5 Click **Save** to save changes, and **OK** to close.


Changing the Temperature Scale

You can display temperature in Fahrenheit or Celsius.

To change the temperature scale:

- 1 Select **Tools | Options**, then choose the **General** tab.

2 Under Temperature scale, select **Fahrenheit** or **Celsius**.

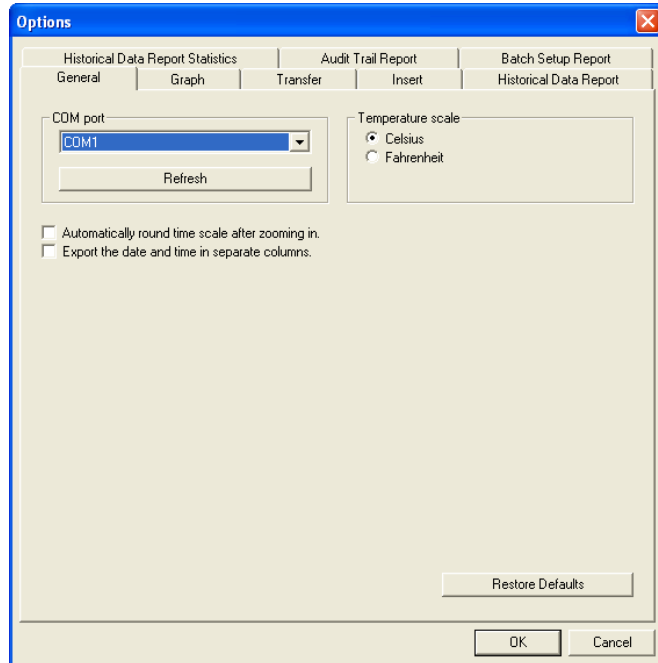

This setting affects the data in all graphs, displays, and printed reports.

Setting Sample Smoothing

In vLog, sample smoothing controls how samples in a logger file are interpreted when they are inserted into a graph, allowing you to dampen the effects of 'noise' in Logger files.

Sample smoothing is the minimum change in input signal required to cause a change in displayed channel values, and is defined in raw analog-to-digital converter counts.


vLog allows you to set two different sample smoothing values: one for thermocouple channels, the other for other channel types.

Note: Sample smoothing only affects how samples are interpreted and presented in graphs and reports; the data in the loggers and the logger files is not

affected by any settings related to sample smoothing.

To set sample smoothing:

- 1 Select **Tools | Options**, then choose the **Insert** tab.

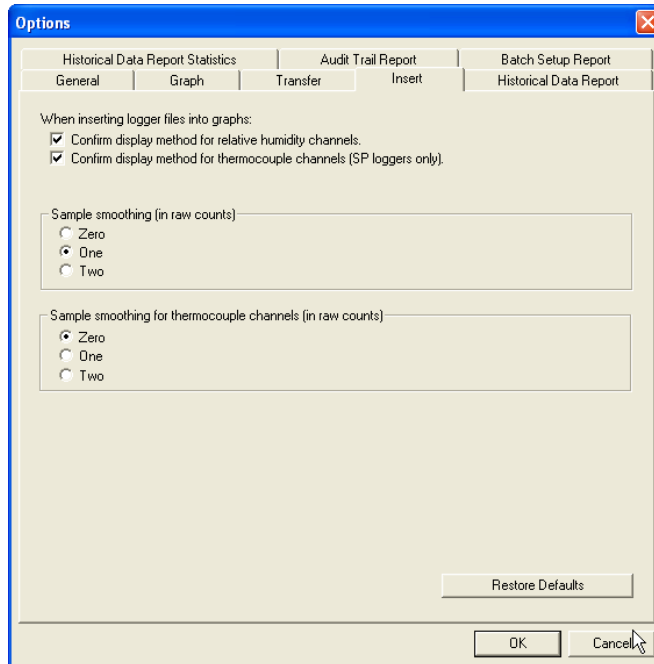

- 2 Under Sample smoothing (in raw counts), select one of the following options.
 - **Zero** - produces the most responsive graph, but makes the recorded readings susceptible to noise
 - **One (default)** - suitable for most applications
 - **Two** - decreases input noise and produces a more stable graph
- 3 Click **OK**.

Adding Horizontal Lines to Graphs

Horizontal lines add value to vLog graphs by showing maximum, minimum, and average values, as well as indicators as chosen by you.

Maximum, minimum and average value lines can be added automatically, or you can customize the lines to add.

To add lines automatically, your graph must only be showing one unit type, such as temperature or humidity. Maximum, minimum and average values will be marked, including the time and date the maximum and minimum were achieved.

You can also add lines manually, such as adding lines to show when a particular threshold has been reached, or to add lines to graphs where more than one unit type is displayed.

To add lines automatically:

- From a valid graph, choose **Format | Auto Horizontal Lines**.

The Graph displays the maximum, average and minimum values in color-coded horizontal lines. Maximum is shown in red, average is green and minimum is blue.

If you zoom in on the Graph, the horizontal lines change to show the maximum, average and minimum values over the selected time period.

To add lines manually:

- 1 From a valid graph file, choose **Format | Manual Horizontal Lines**.
- 2 In the **Format Manual Horizontal Lines** window, enter the value, label text, color, and display units for the horizontal lines to add (up to four colored lines).

Enable:	Value:	Units:	Label:	Color:
☒	0.0	C	Orange Alert Threshold	Orange
☒	0.0	C	Low Temperature Alarm Threshold	Green
☐	0.0	C		Black
☐	0.0	C		Black

OK Cancel


Lines appear on the Graph based on the colors, labels, units, and values you selected.

Displaying Line Numbers on Graphs

vLog displays lines on graphs in different colors. If desired, you can add numbers to the lines in the graph to accommodate monochrome printing.

To display numbers on the graph lines:

- 1 Select **Tools | Options**, then choose the **Graph** tab.

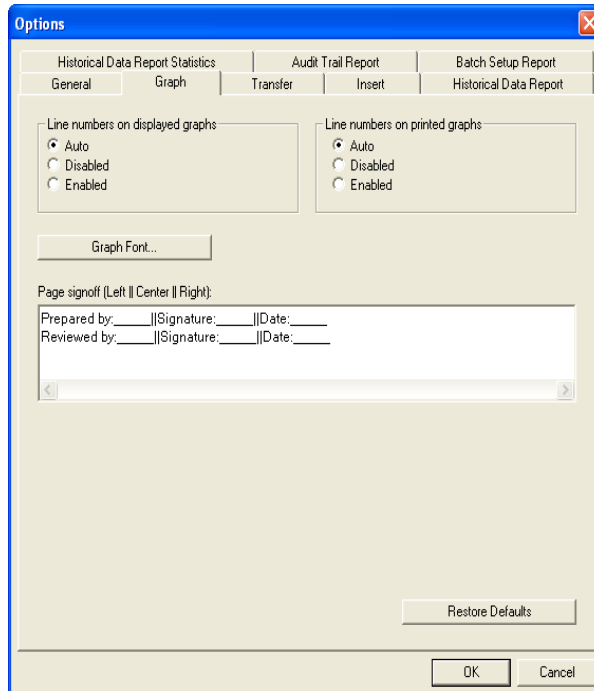

- 2 Under Line numbers on displayed graphs, select **Enabled**.
- 3 Click **OK**.

Verifying a Graph (vLog VL)

If you transferred files using the username and password feature enabled in vLog VL, your graph will display **Not Veri-**

fied in the upper right corner of the graph. Refer to the instructions on verifying graphs on page 106.

Note: This verification method is no longer recommended as vLog now allows you to track Logger File ID numbers (see “Setting Transfer Preferences” on page 52).

Copying and Pasting Graphs

vLog allows you to copy graphs into the Clipboard for pasting into other applications, such as a Microsoft Word file or PowerPoint presentation file.

To copy graph data:

- 1 From vLog, choose **Edit | Copy** (Ctrl+C).
- 2 Go to the application you want to paste the graph image into. In the selected application, choose **Edit | Paste**.


Exporting Graph Data

vLog allows you to export graph data (either .txt or .csv format) for import into spreadsheet and database programs for further analysis and manipulation. You can also specify whether to have date and time saved within one column or separated. As an alternative, you can view and save the Historical Data Report; it contains identical data, plus much more (see “Historical Data Reports” on page 88).

To export graph data to .txt or .csv files:

- 1 By default, date and time are exported in a single column. To have date and time exported into separate columns:
 - a Select **Tools | Options**, then choose the **General** tab.
 - b Select the option, **Export the date and time in separate columns**.
 - c Click **OK**.

- 2 To export the file, choose **File | Export**. The File Save As window opens.

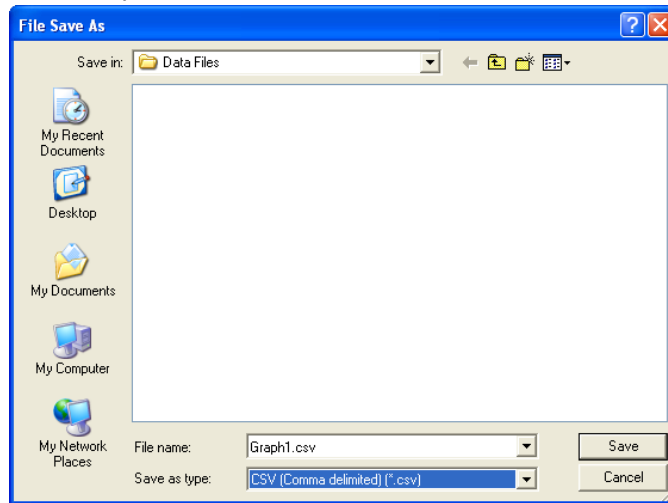

- 3 Navigate to the location where you want to save the file.
- 4 In the File name box, type the desired file name.
- 5 In the Save as type box, accept the default **CSV (comma-separated values)(.csv)** file type, or select **Text (Comma delimited)(*.txt)**.
- 6 Click **Save**.

Saving Graphs

To save changes to a graph file (.spg) for later use:

- 1 Choose **File | Save** or **File | Save As**.
- 2 Navigate to the location to save your file, enter a file name, and click **Save**.

To export saved graph data, see “Exporting Graph Data” on page 83.


Printing Graphs

Printing Line Numbers on Graphs

vLog displays lines on graphs in different colors as visual indicators. For monochromatic printers, vLog automatically adds numbers to the graph lines when the graph is printed.

To change line numbers on graphs:

- 1 Select **Tools | Options**, then choose the **Graph** tab.

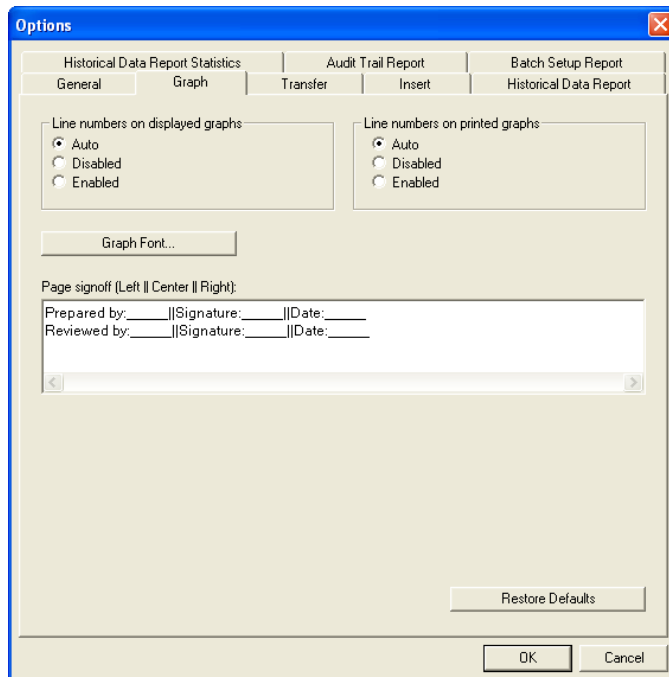

- 2 Under **Line numbers on printed graphs**, select one of the following options.
 - **Auto** - Lines are numbered if the graph is printed on a monochromatic printer and not numbered if the graph is printed on a color printer.
 - **Disabled** - Lines are not numbered on the printed graph, even if it is printed on a monochromatic printer.
 - **Enabled** - Lines are numbered on the printed graph, even if it is printed on a color printer.
- 3 Click **OK**.

Printing Footers on Graphs

To print footers on graphs:

- 1 Select **Tools | Options**, then choose the **Graph** tab.
- 2 Edit the footer content in the Page signoff box. To justify items, use pipe “|” characters, as in left-justified content
|| center justified content || right-justified content.
- 3 Click **OK**.

Print Preview

To view a print preview of the current graph:

- 1 Select **File | Print Preview**. The Print Preview window appears.



- 2 If you want to view the graph in more detail, click **Zoom In** (or **Zoom Out**).
- 3 When you are ready to print the file, click **Print**.

Printing Graphs

To print a currently-active graph:

- 1 Select **File | Print (Ctrl+P)**. The Print window opens.
- 2 Configure the printer settings.
- 3 Click **Print**.

Chapter 5: Viewing Reports and Statistics

Using vLog you can generate a variety of reports to assist you with the analysis of monitoring changes in data logger readings. These reports include:

- Historical Data Reports
- Audit Trail Reports (vLog VL)
- Batch Setup Reports (vLog VL)
- Graph Security Problems report (vLog VL)

You can also copy and paste report data for use in other applications.


Historical Data Reports

Historical data reports display logger file information in a tabular format. You must have at least one logger file inserted into an active vLog graph to view a report. To insert a file, see “Inserting Logger Files” on page 68.

Printed reports can be customized to include specific text in headers and footers, and both print and on-screen reports can be customized to display or hide various statistics and column headings.

To view an historical data report for the currently-active graph:

- 1 Select **View | Historical Data Report**, press F6, or double-click anywhere on the graph.

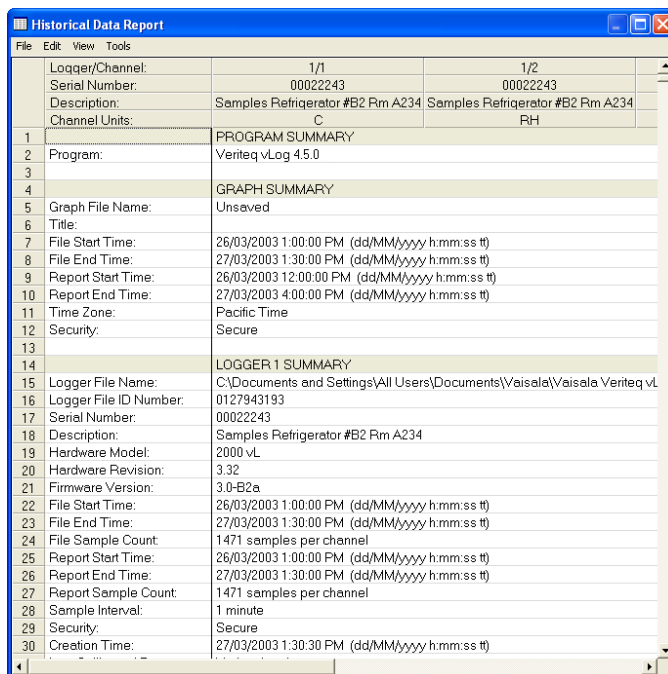

Historical Data Report			
File	Edit	View	Tools
Logger/Channel:	1/1	1/2	
Serial Number:	00022243	00022243	
Description:	Samples Refrigerator #B2 Rm A234	Samples Refrigerator #B2 Rm A234	
Channel Units:	C	RH	
1	PROGRAM SUMMARY		
2	Program:	Veriteq vLog 4.5.0	
3			
4	GRAPH SUMMARY		
5	Graph File Name:	Unsaved	
6	Title:		
7	File Start Time:	26/03/2003 1:00:00 PM (dd/MM/yyyy h:mm:ss tt)	
8	File End Time:	27/03/2003 1:30:00 PM (dd/MM/yyyy h:mm:ss tt)	
9	Report Start Time:	26/03/2003 12:00:00 PM (dd/MM/yyyy h:mm:ss tt)	
10	Report End Time:	27/03/2003 4:00:00 PM (dd/MM/yyyy h:mm:ss tt)	
11	Time Zone:	Pacific Time	
12	Security:	Secure	
13			
14	LOGGER 1 SUMMARY		
15	Logger File Name:	C:\Documents and Settings\All Users\Documents\Vaisala\Vaisala Veriteq vL	
16	Logger File ID Number:	0127943193	
17	Serial Number:	00022243	
18	Description:	Samples Refrigerator #B2 Rm A234	
19	Hardware Model:	2000 vL	
20	Hardware Revision:	3.32	
21	Firmware Version:	3.0-B2a	
22	File Start Time:	26/03/2003 1:00:00 PM (dd/MM/yyyy h:mm:ss tt)	
23	File End Time:	27/03/2003 1:30:00 PM (dd/MM/yyyy h:mm:ss tt)	
24	File Sample Count:	1471 samples per channel	
25	Report Start Time:	26/03/2003 1:00:00 PM (dd/MM/yyyy h:mm:ss tt)	
26	Report End Time:	27/03/2003 1:30:00 PM (dd/MM/yyyy h:mm:ss tt)	
27	Report Sample Count:	1471 samples per channel	
28	Sample Interval:	1 minute	
29	Security:	Secure	
30	Creation Time:	27/03/2003 1:30:30 PM (dd/MM/yyyy h:mm:ss tt)	

- 2 The Historical Data Report window summarizes the information in the Graph file. You can customize which headings and related statistics appear in reports. For more, see “Adding Statistics to Historical Data Reports” on page 89.

- 3 To print the report, from the Historical Data Report window, select **File | Print**. To customize headers and footers on printed reports, see “Customizing Headers, Footers and Font Sizes in Historical Data Reports” on page 91.

Adding Statistics to Historical Data Reports

As you view or print a historical data report, you can also include statistical information derived from the logger file, including maximum, average and minimum values for various intervals within the data in the currently-active graph. Standard deviation, mean kinetic temperature, lethality, sample count and/or group statistics can also be included.

To view statistical summaries in historical data reports:

- 1 With the desired Graph file active, do one of the following:
 - Select **Tools | Options**, then choose the **Historical Data Report Statistics** tab.
 - From the Historical Data Report window, select **Tools | Options**, then choose the **Historical Data Report Statistics** tab.
 - From the Historical Data Report window, right-click and choose **Options | Historical Data Report Statistics**.

- 2 From the Historical Data Report Statistics tab, choose the content to include, such as samples, statistics, or combinations of the two (interlaced or separated).

The screenshot shows the 'Options' dialog box with the 'Historical Data Report Statistics' tab selected. The settings are as follows:

- Report content:**
 - ☐ Samples only
 - ☐ Statistics only
 - ☒ Interlaced samples and statistics
 - ☐ Separate samples and statistics
- Statistics time interval:**
 - Days: 0
 - Hours: 12
 - Minutes: 00
- Statistics time offset:**
 - Days: 0
 - Hours: 08
 - Minutes: 00
- Mean kinetic temperature (MKT):**
 - Activation energy (in kJ/mol): 83.144
- Include:**
 - ☒ Maximum value
 - ☒ Average value
 - ☒ Minimum value
 - ☒ Standard deviation
 - ☐ Mean kinetic temperature (MKT)
 - ☐ Lethality (Fo)
 - ☒ Sample count
 - ☐ Group statistics
- Lethality (Fo):**
 - Reference temperature (in C): 121.11
 - Z value (in C): 10.00
 - D value (in minutes): 1.00

Buttons at the bottom: Restore Defaults, OK, Cancel.

- 3 From the same window, choose your statistics time intervals and offset. For example, to show statistical summaries throughout the report every 12 hours at 8:00am and 8:00pm, set the interval to 12 hours, and the offset to 8 hours.
- 4 If including mean kinetic temperature in your report, specify the activation energy for the item you're monitoring temperature for. This value will be used in the mean kinetic temperature calculation. Enter the activation energy in kJ/mol.
- 5 Choose the statistics to include. Choose from maximum values, average values, minimum values, standard deviation, mean kinetic temperature, lethality, sample count and/or group statistics. Group statistics is for showing data from channels of the same type (for example, temperature or humidity) together in one column.

- 6 If you're including lethality in your report, specify the following values:
 - reference temperature in degrees
 - Z value in degrees (the Z value is the temperature required to change the D value by a factor of 10)
 - D value in minutes (the D value is the time required to destroy 1 log cycle (90%) of the microbial population at the reference temperature)
- 7 Click **OK**. The report is displayed. Scroll through the report to see your statistical information.
- 8 To print your report, from the Historical Data Report window, choose **File | Print**. To customize your report's appearance, see "Customizing Headers, Footers and Font Sizes in Historical Data Reports" on page 91.
- 9 To perform a more in-depth statistical analysis, you can save the report using **File | Save As** from the Historical Data Report window, then analyze the report contents in a spreadsheet program such as Microsoft Excel. See "Saving Reports" on page 102.

Customizing Headers, Footers and Font Sizes in Historical Data Reports

Historical data reports can be customized to include certain items in column headings, certain items and statistics within the report itself, as well as specific headers, footers, and font size in a printed report.

You can also choose to include or not include logger serial numbers, descriptions, and channel descriptions in report column headings for both printed and on-screen reports.

To customize column headings:

- 1 From vLog, select **Tools | Options** then choose the **Historical Data Report** tab.

- 2 From the Historical Data Report tab, select or deselect the items that you want to appear in column headings.

- 3 Click **OK**.

You can also change the font size for printed Historical Data reports. Changing this setting does not affect the on-screen appearance of the report.

To change the font size in printed reports:

- 1 From vLog, select **Tools | Options**, then choose the **Historical Data Report** tab.
- 2 From the Historical Data Report tab, on the **Printer font size** box, select an option from the drop-down list, or enter the desired font size.
- 3 Click **OK**.

To help keep track of your report contents or add instructions, you can add custom headers and footers to printed Historical Data Reports. For example, you can add sign-off lines in footers for signatures and dates.

When you specify a header, what you specify is added to every page in the printed report.

When you specify a footer, you can choose to have what you specify displayed on:

- no pages
- first page
- last page
- first and last pages
- all pages

To edit headers or footers:

- 1 From vLog, select **Tools | Options**, then choose the **Historical Data Report** tab.
- 2 Edit the text in the **Page header** or **Page footer** areas to contain what you want to appear.
 - a To justify items, use the keyboard pipe “||” characters, as shown in the example on-screen. For example, to center-justify a header with nothing on the left or the right, enter, “ || Your Header Text || “. For headers with items to appear on the left and right, enter, “ Left Header Info || || Right Header Info “.
 - b For long underlined spaces (for signatures, for example), type a few “_” or underline characters. vLog is designed to recognize the underline character in headers and footers and stretch the underlines to better fit the page.
 - c For footers, select a **Display page footer on** option to specify the pages on which you want the footer to appear.

Page header (Left || Center || Right):

Page footer (Left || Center || Right):

Prepared by: _____ | Signature: _____ | Date: _____

Display page footer on:

☐ No pages.

☒ First page.

☐ Last page.

☐ First and last page.

☐ All pages.

Printer font size: 7

Restore Defaults

- 3 Click **OK** to save your changes. When printed, the Historical Data report will display the headers and/or footers you configured.

Audit Trail Reports (vLog VL)

The Audit Trail Report, available in vLog VL, contains information from the audit trail file. Audit trail files are similar to event logs, tracking operations that users perform in vLog so auditors can monitor changes to the vLog system. For example, the Audit Trail Report includes the time of an event, event type, severity, description, serial number and description of the logger involved, and computer and user name involved.

The Audit Trail Report is displayed in tabular format and can be saved as a file in multiple formats (see “Saving Reports” on page 102).

vLog also allows you to filter the Audit Trail Report, enabling you to show only events from certain dates, times, performed by certain users, or performed on certain loggers.

Viewing and Filtering Audit Trail Reports

To view and filter the Audit Trail Report:

- 1 Select **View | Audit Trail Report**. The Audit Trail Report Filter window appears. From this window, you can specify the dates and times to show in the Audit Trail Report, as well as certain types of events or events associated with a certain logger, user, or computer.
- 2 To filter by time (optional):
 - a Select **Enable start time filter** checkbox and specify the beginning of the date and time range you want to show in your report. Buttons like **1 month ago**, **1 week ago**, and others make time selection faster.

- b** Select **Enable end time filter** and then specify the end time of the last data you want to appear in your report.

Audit Trail Report Filter

Audit trail folder: C:\Documents and Settings\All Users\Documents\Vaisala\Vaisala Veriteq vLog\Audit Trail\

This folder contains 1 file(s) with 58 event(s) that occurred between 02/05/2011 2:11:52 PM and 19/05/2011 9:41:49 PM.

Start time

☒ Enable start time filter

Date: 19/05/2011 Hours: 12 Minutes: 00 Seconds: 00 AM PM

1 month ago 1 week ago Yesterday Today

End time

☒ Enable end time filter

Date: 20/05/2011 Hours: 12 Minutes: 30 Seconds: 40 AM PM

Today Now

Filter list

Filter	Status
Event Severity	Unfiltered
Event Type	Unfiltered
Event Description	Unfiltered
Logger Serial Number	Unfiltered
Logger Description	Unfiltered
Computer Name	Unfiltered
User Name	Unfiltered

Filter details

☐ Enable event severity filter

Include selected items

☒ Info

Select All Deselect All

Additional filtering

OK Cancel

- 3** To filter by user name, highlight the **User Name** row in the Filter list.
- a** From Filter details, select the **Enable user name** filter box.
- b** From the Include selected items list, select or deselect the user names to include in your report.
- c** If you want to add an additional filter, for example, to display only names beginning with 'S', enter the condition with the wildcard character (*) in the **Additional filter** field (enter **S***).

- 4 To filter by computer name, highlight the Computer Name row in the Filter list.

- a From Filter details, select the **Enable computer name filter** box.

Filter details

☒ Enable computer name filter

Include selected items

☒ CC-LAB-EGDESKTOP

☒ MYDESKTOP

Select All Deselect All

Additional filtering:

*cc

- b From the Include selected items list, select or deselect the computer names to include in your report.
 - c If you want to add an additional filter, for example, display only computer names beginning with 'CC', enter the condition with the wildcard character (*) in the **Additional filter** field (enter **CC***).
- 5 **Other filter options:** You can also filter by event severity, event type, event description, logger serial number, and logger description; following the steps outlined in steps 3 and 4.
- 6 Click **OK**. The Audit Trail Report displays your filtered data.

	Time	Severity	Event Type	Event Description	Serial
1	06/05/2011 1:17:27 PM	Info	System	User account added, User: Guest (Guest Access). ...	
2	06/05/2011 1:17:27 PM	Info	System	User account added, User: User (User Access). Us...	
3	06/05/2011 2:38:58 PM	Info	System	User logoff, User: admin (Administrator), User role: ...	
4	09/05/2011 2:50:40 PM	Info	Logger Transfer	Logger transfer started, File: C:\Documents and Sett...	08
5	09/05/2011 2:51:08 PM	Info	Logger Transfer	Logger transfer succeeded, File: C:\Documents and...	08
6	10/05/2011 3:21:30 PM	Info	System	User logoff, User: admin (Administrator), User role: ...	
7	11/05/2011 3:42:43 PM	Info	Logger Setup	Logger setup started	11
8	11/05/2011 3:42:44 PM	Info	Logger Setup	Logger setup succeeded, Logger cleared, Old CH2 ...	11
9	11/05/2011 3:45:09 PM	Info	System	User logoff, User: admin (Administrator), User role: ...	
10	11/05/2011 3:45:53 PM	Info	Logger Transfer	Logger transfer started, File: C:\Documents and Sett...	11
11	11/05/2011 3:45:54 PM	Info	Logger Transfer	Logger transfer succeeded, File: C:\Documents and...	11
12	11/05/2011 3:46:57 PM	Info	Logger Setup	Logger setup started	11

Configuring Date and Time Columns in Audit Trail Reports

In vLog VL, you have the option of showing date and time in one or two columns in the Audit Trail Report.

To show date and time in two (or one) columns (vLog VL):

- 1 From the Audit Trail Report window, select **Tools | Options**.
- 2 Select or deselect the **Display the date and time in separate columns** check box.
- 3 Click **OK**.

Customizing Headers, Footers and Font Sizes in Printed Audit Trail Reports

To customize the printed Audit Trail Report (vLog VL):

- 1 From the Audit Trail Report window, select **Tools | Options**.
- 2 Edit the text in the **Page header** or **Page footer** areas to contain what you want to appear.
 - a To justify items, use the keyboard pipe “||” characters, as shown in the example on-screen. For example, to center-justify a header with nothing on the left or the right, enter, “ || Your Header Text || “. For headers with items to appear on the left and right, enter, “ Left Header Info || || Right Header Info “.
 - b For long underlined spaces (for signatures, for example), type a few “_” or underline characters. vLog is designed to recognize the underline character in headers and footers and stretch the underlines to better fit the page.

- c For footers, select a **Display page footer on** option to specify the pages on which you want the footer to appear.

- 3 To restore defaults, click **Restore Defaults**.
- 4 Click **OK** to save your changes.

Batch Setup Reports (vLog VL)


The Batch Setup report is a list of loggers set up using the Batch Setup feature. Details include: serial number, description, model number, sample interval, channel unit and description. You can view, edit headers and footers, save, and print the Batch Setup Report.

For more on Batch Setup and configuring data loggers, see “Configuring Loggers: Batch Setup (vLog VL)” on page 46.

To view Batch Setup Reports:

- 1 Select **View | Batch Setup Report**.
- 2 Select the Batch Setup Report (*.bsf) to open. Batch Setup files can be created when you perform batch setups. The Batch Setup Report window opens. The Batch Setup Report window summarizes which, and how, loggers were set up using Batch Setup, and includes the logger serial number, description, logger

model number, channels enabled and their description, and the start time, sample interval, and stop time.

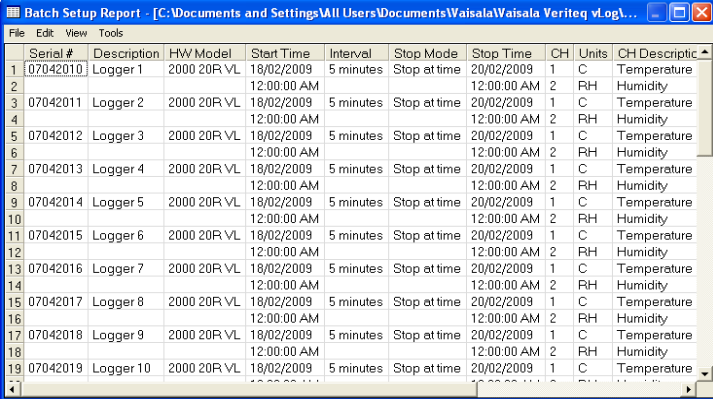

	Serial #	Description	HW Model	Start Time	Interval	Stop Mode	Stop Time	CH	Units	CH Descriptio
1	07042010	Logger 1	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
2				12:00:00 AM			12:00:00 AM	2	RH	Humidity
3	07042011	Logger 2	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
4				12:00:00 AM			12:00:00 AM	2	RH	Humidity
5	07042012	Logger 3	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
6				12:00:00 AM			12:00:00 AM	2	RH	Humidity
7	07042013	Logger 4	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
8				12:00:00 AM			12:00:00 AM	2	RH	Humidity
9	07042014	Logger 5	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
10				12:00:00 AM			12:00:00 AM	2	RH	Humidity
11	07042015	Logger 6	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
12				12:00:00 AM			12:00:00 AM	2	RH	Humidity
13	07042016	Logger 7	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
14				12:00:00 AM			12:00:00 AM	2	RH	Humidity
15	07042017	Logger 8	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
16				12:00:00 AM			12:00:00 AM	2	RH	Humidity
17	07042018	Logger 9	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
18				12:00:00 AM			12:00:00 AM	2	RH	Humidity
19	07042019	Logger 10	2000 20R VL	18/02/2009	5 minutes	Stop at time	20/02/2009	1	C	Temperature
20				12:00:00 AM			12:00:00 AM	2	RH	Humidity


Customizing Headers and Footers in Printed Batch Setup Reports

You can modify the existing header and footer on a currently open Batch Setup report, or modify the format for batch setup reports that will be printed at a later time.

To customize the header and footer in a Batch Setup Report:

- 1 You can modify the header and/or footer in two ways:
 - In a currently open Batch Setup Report, from the Batch Setup Report window, select **View | Page Header and Footer**.
 - To modify the header and/or footer before displaying a Batch Setup Report, from vLog select **Tools | Options**, then select the **Batch Setup Report** tab.
- 2 Add or remove items from the **Page header** and **Page footer** areas.
 - a To justify items, use pipe “|” characters as shown in the example on-screen. For example, for a center-justified header with nothing on the left or the right, enter “ || Your Header Text || ”.

- b To display information on the left or right of a header, enter, “ Left Header Text || || Right Header Text ”.

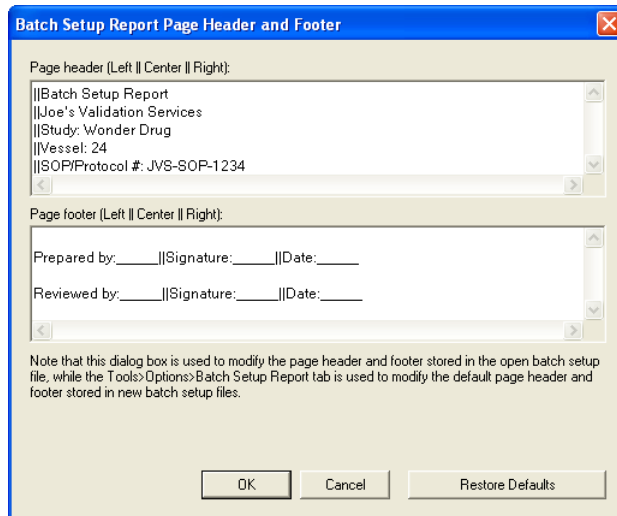

- c For long underlined spaces (for signatures, for example), type a few “_” or underline characters. vLog is designed to recognize the underline character in headers and footers and stretch the underlines to better fit the page.
- 3 Click **OK** to save your changes. When printed, the Batch Setup Report will include the header and footer elements you configured.

Graph Security Problems Reports (vLog VL)


If you see a status indicator appearing in the upper left corner of a graph - Not Secure or Tampered - you can generate a report to review the details.

To generate a Graph Security Problems Report:

- 1 From vLog, select **View | Graph Security Problems Report** (or right-click on the status indicator)..



The Graph Security Problems Report appears:

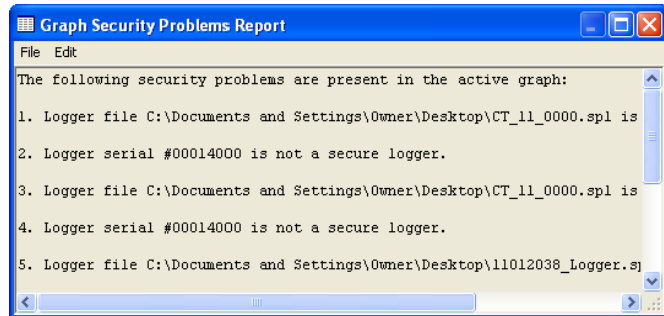

- 2 To print the report, select **File | Print**.
- 3 To save the report as a .txt file, select **File | Save As** and select a file storage location.

Copying and Pasting Report Data

You can copy report data from all vLog reports for use in other applications, such as Microsoft Word or Microsoft Excel.


To copy report data:

- 1 From vLog, choose a report to view from the **View** menu.
- 2 From the report window, choose **Edit | Select All** (Ctrl+A).
- 3 With all contents highlighted, choose **Edit | Copy** (Ctrl+C).
- 4 Go to the application into which you want to paste the report data. In this application, choose **Edit | Paste**.

Saving Reports

To save report data:

- 1 From the report window, choose **File | Save As**.

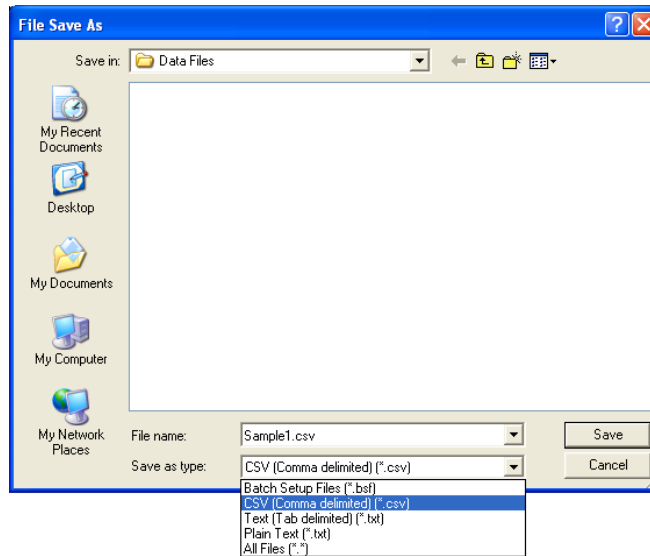

- 2 Navigate to the location where you want to save the file.
- 3 In the **File name** box, enter the desired file name.
- 4 In the **Save as type** box, choose one of the following file types for your file:
 - BSF (Batch Setup file)(.bsf)
 - CSV (Comma delimited)(.csv)
 - Text (Tab delimited)(.txt)
 - Plain text (.txt)

You can readily import these file formats into most popular spreadsheet and database programs for further analysis and manipulation.

- 5 Click **Save**.

Printing Reports

To print reports, from the report window of your selected report, choose **File | Print**.

Appendix: FAQs and Troubleshooting

Software Error Messages

Unable to allocate COM port because it is busy

vLog was unable to allocate the COM port because it was already allocated by a different program. You may have to wait until the COM port becomes available.

Unable to allocate COM port because it is unavailable

vLog was unable to allocate the COM port because it did not exist in the operating system. Check your available COM ports.

Unable to receive sync byte from COM port

vLog sent a request message to the data logger, but did not receive any response from the data logger. This could mean the data logger is unplugged or the cable is not connected properly.

Unable to receive all of the bytes from COM port

vLog sent a request message to the data logger, but received an incomplete response message from the data logger. This may mean there is a problem with the COM port on this PC, a problem with the cable, or a problem with the logger itself. Check these possibilities.

Unable to communicate with logger on COM port

You may receive this error message for miscellaneous error conditions, other than those listed previously. This message could mean that there is a problem with the COM port on this PC, a problem with the cable, or a problem with the logger.

Things to Try:

- Does the logger and interface cable work on a different PC? If so, the problem likely resides in the serial port configurations on the PC.

- Are there other devices that use the serial port in question without problems? If so, that port is likely not the problem.
- Can you transfer data using a **different** interface cable? If so, the cable may be damaged.
- Can you transfer data using a different logger (with the same interface cable)? If so, the data logger may be damaged.

Invalid hardware model

vLog did not recognize the hardware model of the data logger, and as a result, did not know how to interpret the information stored in the data logger. Ensure you are using a data logger supported by vLog.

Unable to create logger file

This error message is generated when vLog attempts to create a logger file on a disk and an error occurs. Some possible reasons include:

- The disk is full.
- The user does not have write access to the disk. Ensure the user has write permission to the folder where the logger file is being saved.
- The file is already open by a different program. Close it.
- The file name is invalid. Try using a more simple .spl file name.

Frequently Asked Questions

Q: I have recently upgraded Vaisala Veriteq vLog software to the latest version. Will this new software work with my existing loggers and files?

A: Yes. Each version of Vaisala Veriteq vLog software is designed with backward compatibility. This means that it works with all previous versions of loggers or files created from those loggers.

Q: I have recently upgraded to Vaisala Veriteq vLog software from an older version of Spectrum software. Will this new software work with my existing Spectrum loggers and files?

A: Yes. Vaisala Veriteq vLog (vLog VL or vLog SP) software is designed with backward compatibility and supports Spectrum

loggers and Spectrum software files. This means it works with all previous versions of loggers or files created from those loggers. The key difference is that Spectrum loggers use a non-secure file format.

Q: I have recently purchased a new logger. Will it work with my older version of vLog software?

A: Not always. In our pursuit of continual product improvement, we periodically upgrade logger hardware or release new logger models. Unfortunately, these changes may not be recognized by an older version of Vaisala Veriteq vLog software and can result in difficulties. For that reason, each logger that you purchase includes a note indicating which software version it must be used with. The safest choice is to always upgrade to the latest version of software. Contact Vaisala Canada Inc if you require a software update.

Q: What is the difference between Vaisala Veriteq VL-Series Loggers and Vaisala Veriteq SP-Series Data Loggers (formerly Spectrum)?

A: The primary difference is the VL-Series models provide you with the ability to obtain additional information and utilize security features. Vaisala Veriteq VL-Series Loggers allocate a portion of their internal memory to store important calibration-related parameters.

The stored parameters include the following:

- The company or person who last calibrated the data logger
- The last calibration date
- The next calibration date
- Security information (to prevent unauthorized calibration or manipulation of the data in the logger).

Q: Does the software account for daylight savings time changes?

A: Yes. Daylight savings time changes are automatically accounted for in all graph and tabular displays.

Q: How many logger files can I overlay onto the same graph?

A: You can overlay many channels of data onto the same graph; in testing, we inserted more than 100 channels successfully.

Q: Do all graphs have to be created from one data logger?

A: No. You can create graphs using files from different data loggers and different data logger models.

Q: I have a version of vLog earlier than version 4.5. What happened to that Verification process? And why, in Version 4.10 and later, has this been replaced by the Logger File ID Number?

In Vaisala Veriteq vLog, logger data and graph files can be generated at different times by different people. Consequently, there is a need for verification that the logger files used to create graph files are the right ones. This has been handled in vLog versions prior to version 4.10 by the verification process that required a transfer username and password to be entered in order to verify graphs.

The Logger File ID Number is a more flexible approach that allows users to compare the printed records they have against the logger files they are using to create graphs.

However, if you are an existing customer with VL-Series loggers you can still use usernames and passwords for verification.

Note: Do not use your Windows username and password, as you may be required to share the username and password with others. Also, the username and password option cannot be used for Batch Transfers in vLog VL.

To configure the username and password prompting option:

- 1 Select **Tools | Options**, then choose the **Transfer** tab.
- 2 Select the **Prompt for username and password when transferring samples** checkbox.
- 3 Click **OK**. See page 107 for instructions on how to verify a graph with a username and password.

Q: How can I verify a graph on the screen if I have forgotten my password?


A: You can't. Vaisala Veriteq vLog VL software requires the original password when transferring the data logger in order to verify the status of the file. This is a security feature for validation to ensure that the person responsible for transferring the data is also the person who creates the graph. For vLog 4.10 and later, the username and password feature is optional and can be

turned off. It has been replaced with a new Logger File ID Number (see, “Transferring Data” on page 52).

If you transferred files using the username and password feature enabled in vLog VL, you are required to complete these next steps to verify the graph.

To verify a graph:

- 1 Select **Tools | Verify Passwords in Logger Files** or right-click on the words **Not Verified**, in the upper right corner of the graph.

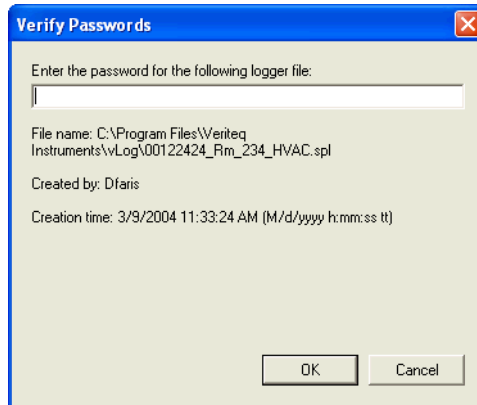

- 2 Type the password you used when you originally created the file.
- 3 Click **OK**.

Glossary

#

21 CFR Part 11 (vLog VL)

A Food and Drug Administration (FDA) rule relating to the use of Electronic Records and Electronic Signatures. The FDA's view is that the risks of falsification, misinterpretation, and change without leaving evidence are higher with electronic records than paper records, and, consequently, 21 CFR Part 11 has been developed to provide adequate control over these risks.

A–B

Ambient temperature

The average or mean temperature of the surrounding air that comes in contact with the equipment or instrument.

Axis

A line that borders one side of the graph area, providing a frame of reference for measurement or comparison. In vLog software, measurement values are plotted on the vertical Y-axis while date and time values are plotted on the horizontal X-axis.

C–E

Calibrate

To compare and correlate readings against a standard measurement

Calibration certificate

Documented evidence that a product has been calibrated against a standard measurement reference and that it falls within the specified requirements

CFR

The acronym for Code of Federal Regulations; this is a book of federal rules and regulations that are followed by government agencies

Channel

The path through which an instrument receives the input signals it measures

Dewpoint

Dewpoint is defined as the temperature at which the amount of moisture present in the air is the MAXIMUM amount the air can hold (at a constant pressure and water vapor content). This means that when the temperature falls to the dew point, the relative humidity is 100% and the air is considered saturated.

Dewpoint does not vary with temperature, as does relative humidity (RH). Rather, dewpoint is a DIRECT measure of the actual moisture content in the air and thus is preferred as the unit of measurement in many scientific and industrial applications—particularly in tightly controlled environments where precise indications of water content in the air are required.

The dewpoint is never greater than the air temperature. Once the dewpoint is reached and the air temperature continues to drop, water vapor is forced to leave the air through condensation. This condensation may take the form of water droplets, dew, or may remain suspended in the air as clouds or fog.

To better understand how dewpoint and relative humidity are related to each other, try inserting a logger file with relative humidity (RH) channel values into a graph twice—once, displaying RH using dewpoint and once displaying RH using relative humidity.

in vLog, dewpoint readings below 0 degrees C are recorded as frostpoint.

F–K

FDA

An acronym for the US Food and Drug Administration; this is a government agency responsible for enforcing the Federal Food, Drug and Cosmetic Act

GMP

An acronym for Good Manufacturing Practices; these are the regulations found in the CFR, parts 210 and 211

Graph file

A file that is created to display the data from one or more Logger files. Graph files retain formatting information such as graph title, zoom levels, Y-axis scaling and Logger file data and are identified by their .spg extension (for example, LAB_1.spg)

L–Q**Local time**

A value equal to the current time used at a particular location. A relative indicator of time meaningful only if Time Zone Information is also known, Local Time contrasts with UTC Time, an absolute time reference. The software records Local Time and UTC Time in each Logger file.

Logger file

A file that is created after transferring a logger's data to the PC. The file contains the collected readings (from all active channels) for the transferred logger including the logger's description, model number, serial number, firmware version number, file description, and Local time. Logger files are identified by their .spl extension (for example, PRIMARY_GAS_FLOW.spl). Logger files cannot be viewed directly—you have to insert the Logger file into a Graph file before the data is displayed as a graph

NIST

An acronym for the National Institute of Standards and Technology

R–T**Range**

Normal operating limits, specified by the lowest calibration point to the highest calibration point.

Relative humidity

Relative humidity is the measure of the amount of water vapor in the air compared to how much it can possibly hold at that temperature. A relative humidity of 50% means that the air contains only one-half as much water as it could possibly hold. The capacity of air to hold water is heavily dependent on the temperature of the air. The colder the air, the less water it can hold. Air at 50°F can hold only one-quarter the amount of water as air at 95°F

Sample

A measurement reading that is stored to data logger's memory

Sample smoothing (deadband)

A setting in vLog (**Tools | Options | Insert** tab - Sample smoothing) that controls how the samples in a logger file are interpreted when they are inserted into a graph. Sample smoothing (or sample deadband) is the minimum change in the input signal required to cause a change in the displayed channel values. Sample smoothing is used to dampen the effects of noise, and is defined in terms of raw analog-to-digital converter counts.

Sample interval

The sample interval is the rate at which the logger takes and records readings to memory. The sample interval applies to all enabled channels on the logger. You may want to occasionally change the sample interval on the logger because:

- Faster sample rates may be required by the application
- You need to record readings over a longer duration
- You want to minimize the collection of unnecessary data

Scale

The range of values a graph displays

Secure

The status of an on-screen or printed graph file (.spg) that has been created in vLog using data from a VL-Series Data Logger. Although vLog displays graph files produced from SP data loggers, they are labeled **Not Secure**. Only VL-Series loggers produce secure files.

Serial port

The communications port to which you can attach devices such as a modem, mouse, or serial printer to the computer. The serial communication ports on the computer are commonly referred to as COM1, COM2, COM3, and so on. Serial ports exchange data sequentially, one unit at time.

Stability

The ability of a sensor or measurement instrument to maintain its performance characteristics over a specified period of time

Thermistor

A temperature-sensing element (a *thermally*-sensitive *resistor*) composed of sintered semiconductor material that exhibits a large change in resistance in proportion to a small change in temperature. The relationship between temperature and resistance is approximated for most thermistors by the Steinhart-Hart equation.

Traceable

Capable of being traced back to a recognized measurement standard or reference

Transfer

The term used when the software makes a copy of the data contained in a connected data logger. The transfer process creates a file on the PC but does not affect or delete the data on the logger. Other terms commonly used to describe the transfer process are download, upload, backup, or copy. Transfers can be performed one logger at a time or in batches.

U–Z

USB

An acronym for *Universal Serial Bus*, a standard for connecting external devices such as a mouse, keyboard, scanner, printer. USB offers many benefits over serial and parallel port connections, including thinner and cheaper cables, greater expandability (with the addition of a USB hub, a single USB port can handle over 100 peripheral devices) and greater speed.

UTC

An acronym for Universal Time Coordinated, UTC is a standard time common to every place in the world. Also known as Greenwich Mean Time (GMT) and World Time; UTC is expressed using a 24-hour clock but can be converted into a 12-hour clock (AM and PM). An absolute time reference that does not rely on Time Zone information, UTC contrasts with Local time. UTC and Local Time information is contained in Logger files and is used to form the time base on logger graphs and data.

Validatable

Capable of being validated. Three issues are essential before equipment can be considered validatable:

- There must be documented evidence of performance verification available. As an example of this requirement for documented evidence, a NIST – traceable calibration certificate is provided with every VL-Series Data Logger
- The equipment must match the specific requirements of the application. Each application may involve special accuracy, operating range, performance stability or other requirements that the product must meet
- The equipment must be capable of maintaining a validated state throughout the application process. The product or its output cannot be manipulated, changed or tampered with

Validation

Validation is the process of proving that something does what it is supposed to do. The FDA's definition is “establishing through documented evidence, a high degree of assurance that a specific process consistently produces a product that meets its predetermined specifications and quality characteristics”.

Verified

The status of an on-screen or printed graph file (.spg) that has been “verified” by entry of the same password (using **Tools | Verify Passwords in Logger Files**) used when the file was created (during logger transferring).

Note that the **Verified** status (**Verified** or **Not Verified**) does not appear in vLog graphs if the optional password was not entered during the transferring of the original logger data.

The default verification status of any Graph file that has used a transfer password is **Not Verified**.

Zooming

Enlarging a portion of an on-screen graph

Index

Numerics

21 CFR Part 11 109

A

accounts, user 23

activation energy 90

administrator options 20, 22

audit trail file path 21

security levels 22

user accounts 23

ASCII file

creating 83, 84

audit trail 19, 28

link 30

linking loggers to 31

paths to 21

reports

date & time columns 97

saving reports 102

Audit Trail Reports 94–98

averages 89, 90

axis 109

unit precision 78

B

batch setup 27, 46–48

saving reports 102

Batch Setup Reports 98–100

batch transfers 64

C

cables 16

calibration 15, 30, 109

celsius 78

CFR 109

21 CFR Part 11 109

changing

sample interval 37–38

time zone 72

user accounts 24

channels 39, 110

DC current 40

description 40, 44, 47

disabling or enabling 39

factory-disabled 27

samples per channel 30

thermocouple 60, 69

transducers 39

voltage 41

checklist for setting up vLog 14

checksums 63

clear loggers 34

clear mode 33

clipboard, copying graphs or report

data to 8

COM ports 103

selecting 25

setting up 25

troubleshooting 103

configuring

data loggers 26

data loggers, individually 27

loggers

thermocouples 43

loggers in a batch 46–48

copying report data 101

csv files 83

cursors

positioning 71

D

data

exporting 83

transferring to PCs 54, 64

data logger. See loggers.
 data recorder. See loggers.
 date & time columns in audit trail reports 97
 dates, setting 64
 daylight savings time
 loggers awareness of 105
 DC current channels 40
 deadband 112
 decimal place precision 78
 description 28, 31, 32
 channels 40, 44, 47
 loggers 32, 47
 descriptions 31
 dewpoint 110

F

factory-disabled channels 27
 fahrenheit 78
 file security status 5
 files 83, 102
 graph 111
 inserting into graphs 68
 logger 52
 spg 111
 spl 52, 111

filtering
 audit trail report 94
 firmware revision 28
 fonts
 in graphs 73
 in printed reports 91, 92, 97
 footers, in reports 97, 99

G

Graph Security Problems Report 100
 graphs 111
 adding max, min, avg lines 80
 copying into clipboard 8, 83

creating 68
 displaying 68
 exporting data 83
 files 111
 fonts 73
 horizontal lines 80
 line number display 82
 opening 68
 positioning cursors 71
 print preview 86
 printing 85, 86
 printing with line numbers 85
 saving 84
 scrolling 71
 security problems report 100
 time zone 72
 titles 72
 verification 107
 zooming 73, 115

H

hardware model 28
 hardware revision 28
 headers, in reports 97, 99
 Help
 contacts 9
 high security 20, 22
 Historical Data Reports 88–94
 customizing headers and footers 92
 saving 102
 horizontal lines 80
 adding automatically 81
 adding manually 81
 humidity sensor 50

I

icons on the toolbar 7
 inserting logger files 68
 installation 12–14

intervals, in statistics 90

intervals, sample 37–38

IQOQ documents 15

L

lethality 89, 90

levels of security, normal security 20, 22

line numbers 82, 85

lines, horizontal 80, 81

log time 30, 37

Logger File Creation Report 53, 58, 60–63, 106

logger file ID 60

logger file ID number 63

loggers

- batch setup 27

- clear mode 33

- clearing 34

- clearing samples 44

- configuration 26

- batch setup 46–48

- configuring individually 27

- connecting to PCs 16

- daylight savings 105

- description 28, 31, 32, 47, 61

- disabling or enabling channels 39

- factory-disabled channels 27

- linking to audit trail files 31

- logger file filenames 53

- logger file ID number 63

- models 4, 28

- mounting 48

- mounting tips 48

- operation when full 35

- recording data 36, 37–38

- recording time 111

- sample intervals 37–38

- sampling 28

- serial numbers 28

- setting up 3

- setting up COM ports 25

- setup 26

- status 28

- transferring data 54, 64

- validation 15

- ventilation 48

- VL-series

- calibrating 15

- wrapping data 36

- logging on to vLog 17, 19

- logging on to vLog VL 21

M

maximums 89, 90

mean kinetic temperature 90

menus & keyboard shortcuts 6

minimum requirements, system 12

minimums 89, 90

models, loggers 4

modes

- clear mode 33

- start 34

- stop 29, 34

- stop mode 35

mounting data loggers 48

N

numbers, on graph lines 82

O

offset

- temperature 44, 60, 69

- time interval 90

options

- administrator 22

- security levels 22

P

passwords

- for file verification 6

- logger configuration 27
- logging on 17, 21
- security level 20
- ports
 - serial 113
 - setting up 25
 - troubleshooting 103
- position cursor 71
- printing
 - graphs 85, 86
 - print preview 86
 - reports 102
- probes 49
- protocol documents 15
- R**
- recorders, data. See loggers.
- relative humidity 112
- relative humidity sensor 50
- removing
 - user accounts 25
- reports 88–102
 - adding statistics to 89
 - Audit Trail report 94
 - date & time columns 97
 - headers and footers 97
 - Batch Setup report 98
 - headers and footers 99
 - column headings available 91
 - copying data to the clipboard 8, 101
 - filtering audit trail 94
 - font size in printed 91, 92, 97
 - Graph Security Problems report 100
 - Historical Data report 88
 - customizing 92
 - Logger File creation 106
 - Logger File Creation report 58, 60
 - printing 102
 - saving 102
 - statistics in 90

- requirements 12

S

- sample count 90
- sample intervals 29, 37–38, 112
- sample smoothing 112
- sample timing 38
- samples 30
 - clearing from loggers 44
 - starting taking 34
 - stopping taking 37
- sampling 28
- saving graphs 84
- saving reports 102
- scale
 - x-axis 6, 109
 - y-axis 6, 77, 109
- scrolling, graphs 71
- secure status 113
- security 5, 22
 - report 100
- security key 13
 - vLog-VL 18
- security levels 20, 22
- Security Wizard 19
- selecting COM ports 25
- serial numbers 28
- serial ports 113
- setting up
 - checklist 14
 - COM ports 25
 - connecting loggers to PCs 16
 - loggers in a batch 46–48
 - vLog 3
- Spectrum loggers. See SP-series loggers.
- 68
- spl files 52, 111
- SP-series loggers
 - sample interval 29, 38
 - security 70

- using with vLog 2
 - vs. VL-series loggers 105
- standard deviation 90
- start mode 34
- start time 29
- statistics 89, 90
 - calculation intervals 90
 - lethality 89
- status
 - file security 5
 - logger 28
- stop mode 29, 34, 35
- stop time 30, 37
- Stop when full, defined 36
- system requirements 12

T

- tabular listing 88, 98
- temperature offset 44, 60, 69
- temperature probe 49
- temperatures
 - changing units 78
 - mean kinetic 90
- thermistor 113
- thermocouple channel 60, 69
- thermocouples 43
 - configuring for 43
- time 64, 111, 114
 - log time 30
 - sample 38
 - setting system time 64
 - stop time 30
 - UTC 114
- time zone, changing 72
- titles, graphs 72
- toolbar 7
- transducers 39
 - configuring VL-4000 loggers to work with 40

- transferring logger data to PCs 52, 64, 113
- troubleshooting 103
- txt files 83

U

- user accounts 23
 - adding 24
 - changing 24
 - removing 25
- user names 17, 24

V

- validatable 114
- validation 14, 15, 115
- verification 6, 106, 115
 - Logger File Creation Report 106
- verifying graphs 107
- VL-1700 loggers 4, 43
- VL-4000 loggers 4, 39, 40
- vLog
 - calibrating loggers 15
 - features overview 2–3
 - getting started 14
 - installation 12–14
 - interface overview 5
 - logging on 17
 - menus & keyboard shortcuts 6
 - protocol documents 15
 - security key 13, 18
 - Security Wizard 19
 - setting up 3
 - setup checklist 14
 - toolbar 7
 - vLog VL or vLog SP? 3
- vLog VL
 - and Windows security 21
 - audit trail file path 21
 - audit trail files 19
 - batch setup 46–48

- features 3
- linking loggers to audit trail 31
- logging on after first log on 21
- logging on for the first time 17
- security 22
- user accounts 23
 - changing 24
- VL-series loggers
 - calibrating 15
 - sample interval 29, 38
 - vs. SP-series loggers 105
- VL-Series Validatable Data Loggers 4
- voltage channels 41

W

- warmup time, configuring 42
- Windows security 21
- wrapping data 36

X

- x-axis 6, 109
- XPS excitation control switch 42

Y

- y-axis 6, 77, 109

Z

- zooming in and out 73, 76, 115



Distributed by:
Kenelec Scientific Pty Ltd
1300 73 22 33
sales@kenelec.com.au
www.kenelec.com.au