

# AEROTRAK® PORTABLE PARTICLE COUNTER MODEL 9500



The TSI AeroTrak® 9500 Portable Particle Counter offers best-in-class features and versatility. Perform ISO 14644-1 classifications and complete certifications of EU GMP facilities by sampling 1 m<sup>3</sup> of air in just 10 minutes! Particle count data can be conveniently viewed on screen, downloaded using TrakPro™ Lite Secure software, or printed directly to the integrated printer. Compatible with TSI's accurate ventilation probes that measure air velocity/temperature/RH, the Model 9500 is a complete solution for your certification needs. The combination of an optional electronic filter scanner with up to 26 ft of tubing, makes this the best instrument for in-place filter testing.

This particle counter can generate pass/fail reports for ISO 14644-1, EU GMP Annex 1, and FS209E, via the onboard printer and the supplied TrakPro Lite Secure download software. Multiple particle counter configurations can be conveniently stored and uploaded as needed. It can be used as a stand-alone particle counter or integrated into a facility monitoring system like TSI's FMS 5.


The Model 9500 complies with all the stringent requirements set forth in ISO 21501-4. It is calibrated with NIST traceable PSL spheres using TSI's world-class Classifier and Condensation Particle Counters, the recognized standard for particle measurements. Backed with a standard two-year warranty and TSI's long-standing reputation for high quality, there are no other particle counters like it on the market today.

## Features and Benefits

- + 0.5 to 25 µm size range
- + Up to six channels of simultaneous data
- + 3.53 CFM (100 L/min) flow rate
- + Complies with all requirements of ISO 21501-4
- + Best choice for filter testing in cleanrooms
- + Data storage: 250 zones, 999 locations, and 10,000 samples
- + USB, USB Storage Device and Ethernet data transfer
- + Industry's loudest audible alarm
- + Easy icon-driven touch screen operation
- + Onboard reports for ISO-14644-1, EU GMP Annex 1, and FS209E
- + External certification reports via TrakPro Lite Secure
- + Barcode scanner supports recall of specific Recipe and Zone configurations
- + Displays up to three environmental parameters
- + Stainless steel enclosure
- + Lightweight
- + Stand-alone or integrate into a facility monitoring system
- + Battery life provides for >3 hours of continuous operation



UNDERSTANDING, ACCELERATED

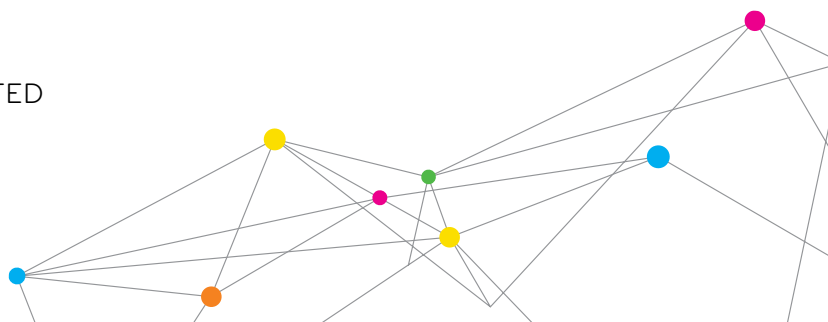

## SPECIFICATIONS

### AEROTRAK PORTABLE PARTICLE COUNTER MODEL 9500

9500			
Size Range	0.5 to 25 µm	Data Storage	250 Zones
Particle Channel Sizes	0.5, 0.7, 1.0, 3.0, 5.0, 10.0 µm		999 Locations
Size Resolution	<15% @ 0.5 µm (per ISO 21501-4)		10,000 sample records including: date, time, six total particle size channels, flow status, and instrument status
Counting Efficiency	50% at 0.5 µm, 100% for particles > 0.75 µm (per ISO 21501-4 and JIS)	Data Transfer	Via USB storage device, connection to the TrakPro Lite Secure software through Ethernet or USB cable, or optional TSI FMS Software
Concentration Limit	425,000 particles/ft <sup>3</sup> (15,000,000/m <sup>3</sup> ) @ 10% coincidence loss	Reports	Provides Pass/Fail on ISO 14644-1, EU GMP, and FS209E reports
Light Source	Long life laser diode	Calibration	NIST traceable using TSI calibration system
Zero Count	<1 count per 5 minutes (per ISO 21501-4 and JIS B9921)	Calibration Frequency	Recommended minimum of once per year
Flow Rate	3.53 CFM (100 L/min) ±5% accuracy (meets ISO 21501-4 and JIS B9921)	Dimension (H x W x D)	9.5 in. x 9.0 in. x 9.2 in. (24.1 cm x 23.0 cm x 23.2 cm) (without isokinetic inlet)
Flow Rate Control	Electronic, automatic closed loop (patented* flow control technology)	Weight	14.2 lbs (6.44 kg) with two batteries
Sample Tube Extension	Up to 26.2 ft (8 m)	Power	110 to 240 VAC universal power supply
Sampling Modes	Manual, automatic, beep; cumulative/differential; count or concentration	Battery	Removable/rechargeable Li-Ion batteries (2 battery slots)
Sampling Time	1 second to 99 hours	Battery Life	>3 hours of continuous operation
Sampling Frequency	1 to 9999 cycles or continuous	Battery Recharge	5 hours (when idle)
Exhaust	Internally HEPA filter	External Surface	Stainless steel
Vacuum Source	Internal pump	Operating Range	36° to 95°F (2° to 35°C), 20% to 95% RH noncondensing
Communication Mode	Modbus® TCP over Ethernet or USB	Storage Range	32° to 122°F (0° to 50°C), up to 98% RH noncondensing
Audible Alarm	Built-in; >85 dB at 1 meter (adjustable)	Standards	ISO 21501-4, CE, JIS B9921
Status Indicators	Low battery, flow, laser	Warranty	Two years
Alarm Output	Dry contacts, closed when alarm is engaged	Included Accessories	Printed QuickStart guide, operating manual on CD, power supply, battery, barbed inlet, isokinetic probe, tubing, purge filter, printer paper, USB cable and TrakPro Lite Secure Software
Alarm Limits	Programmable for all particle channels	Optional Accessories	Electronic filter scanning probe, basic filter scanning probe, TSI velocity probes, Temp/RH probe, additional battery, dual port external battery charger, isokinetic probes, sample tubing, carrying case
Display	VGA 5.7-inch (14.5 cm) touch-screen display		
Languages	English, German, French, Spanish, Japanese, Chinese (simplified), Italian		
Printer	Built-in thermal printer (also available without printer)		
Printer Output	Prints in all available languages		
Environmental Sensor Interface	Supports TSI Probe Models: 960, 962, 964, 966 (air velocity/temp/RH) and Temp/RH probe		
External Chemical Resistance	Isopropyl alcohol, chlorinated solution, hydrogen peroxide		
Software	Compatible with TrakPro Lite Secure and FMS 5 software		
Unit ID	Configurable IP address		
Security	2-level password protection to lock out usage and configuration		

\*Patent Number 6,167,107

Specifications are subject to change without notice.

AeroTrak, TSI and the TSI logo are registered trademarks, and TrakPro is a trademark of TSI Incorporated.

Modbus is a registered trademark of Modicon, Inc.



UNDERSTANDING, ACCELERATED

TSI Incorporated - Visit our website [www.tsi.com](http://www.tsi.com) for more information.

**USA** Tel: +1 800 874 2811  
**UK** Tel: +44 149 4 459200  
**France** Tel: +33 1 41 19 21 99  
**Germany** Tel: +49 241 523030

**India** Tel: +91 80 67877200  
**China** Tel: +86 10 8219 7688  
**Singapore** Tel: +65 6595 6388



Distributed by:

Kenelec Scientific Pty Ltd  
 1300 73 22 33  
[sales@kenelec.com.au](mailto:sales@kenelec.com.au)  
[www.kenelec.com.au](http://www.kenelec.com.au)