

# FITPRO+™ FIT TEST SOFTWARE FOR PORTACOUNT® PRO AND PRO+ RESPIRATOR FIT TESTERS

FitPro+™ Fit Test Software makes fit testing simple and foolproof. Included with every PortaCount Pro and PortaCount Pro+ Respirator Fit Tester, the program automatically leads you through the fit test protocol step by step – you can't make a mistake. Your fit test results are stored in a database for easy report generation, retrieval and regulatory compliance.



## Computer Requirements


- + Microsoft® Windows® 7 (32 and 64 bit), Windows 8 (32 and 64 bit), or Windows 10 (32 and 64 bit) operating system
- + Microsoft .Net version 4 (installation software will install Microsoft .NET automatically, if necessary)  
**Note**—An Internet connection is required to install Microsoft .NET
- + >1 GHz processor
- + >1280x800 screen resolution monitor
- + 500 MB free space on hard drive
- + 1 free USB port
- + CD-ROM reader or internet connection (to install the software)

## Features and Benefits

- + Operate up to four (4) PortaCounts simultaneously using one computer (PortaCount Models 8030/8038 only)
- + Select between new Horizontal or Vertical Fit Test Dashboard
- + Automates fit test process using built-in fit test exercise protocols
- + Microsoft Access® default database (optional support for server-based Microsoft SQL Server®)
- + User Interface in five languages (English, German, Polish, French and Spanish)
- + Comprehensive database management functionality (merging, moving, copying and deleting records; exporting databases; and backup and restoring of data)
- + Maintains separate databases for different departments/locations
- + Easily transfers records between PortaCount Fit Testers using a flash drive
- + Wide-range of reporting and printing capabilities, including individual fit test reports, due or past due fit test reports, custom reports, and printing of individualized fit test cards (using card-form paper or optional card printer)
- + Real-time fit factor display and fit testing emulators for respirator training and familiarization



UNDERSTANDING, ACCELERATED

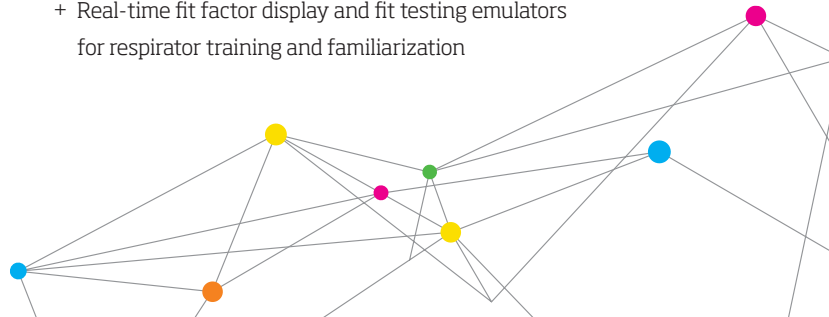

# SPECIFICATIONS

## FITPRO+ FIT TEST SOFTWARE FOR PORTACOUNT PRO AND PRO+ RESPIRATOR FIT TESTERS

### Fit Test Report and Card Printing Options

FitPro+ is designed to improve your respirator program and help protect workers through easy generation of individualized fit test reports inclusive of personal wallet cards. The wallet cards provide documentation of the completed fit test and easy field verification of proper respirator size and model for authorized respirator users.

The reports can be printed by plugging a USB printer directly into the PortaCount Pro/Pro+ Respirator Fit Tester. Exclusive card-form paper is available direct from TSI that provides for the convenient, durable tear-out wallet card.

FitPro+ Fit Test Software also maintains the ability to print fit test results to an external ID card printer.

FIT TEST REPORT			
ID NUMBER	1234	PORTACOUNT S/N	803885303
LAST NAME	BRUCE	N95 COMPANION	N
FIRST NAME	MORTON		
COMPANY			
LOCATION			
TEST DATE	1/22/2014 10:52 AM	PROTOCOL	OSHA 29CFR1910.134
DUE DATE	1/22/2015	PASS LEVEL	500
RESPIRATOR	MSA ULTRATWIN FULL FACE (500)	APPROVAL	EFFICIENCY=99% Fail#
MANUFACTURER	MSA		
MODEL	ULTRATWIN		
MASK STYLE	FULL FACE		
MASK SIZE	F		
EXERCISE	DURATION(S)	FF FACTOR	PASS
NORMAL BREATHING	64	110391	Y
DEEP BREATHING	64	34756	Y
HEAD SIDE TO SIDE	64	219463	Y
HEAD UP AND DOWN	64	82000	Y
TALKING	64	70031	Y
GRIMACE	15	Encl.	Y
BRENDING OVER	64	165212	Y
NORMAL BREATHING	64	64501	Y
OVERALL FF		77582	Y
FITTEST OPERATOR	JDM	DATE	
NAME	MORTON BRUCE	DATE	
NOTE:			

**Respirator Fit Test Card**

Name: MORTON BRUCE      Test Date: 1/22/2014  
 ID: 1234      Next Test Date: 1/22/2015

<b>Respirator</b>	<b>Results</b>
Mfg: MSA	Overall FF: 77582
Model: ULTRATWIN	FF Pass Level: 500
Style: FULL FACE	Pass: Y
Size: F	Operator: JDM
	Protocol: OSHA 29CFR1910.134
	Fit Test Method: QNFT using TSI PortaCount
*** Your company contact information here ***	

Licensed under US Patent No. 5,117,190 and other patents from Commissariat a L'Energie Atomique, France. U.S. Patent No. 6,125,845. Other patents pending.

PortaCount, TSI and the TSI logo are registered trademarks, and FitPro is a trademark of TSI Incorporated.

Access, Microsoft, SQL Server and Windows are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.



UNDERSTANDING, ACCELERATED

TSI Incorporated - Visit our website [www.tsi.com](http://www.tsi.com) for more information.

<b>USA</b>	<b>Tel:</b> +1 800 874 2811	<b>India</b>	<b>Tel:</b> +91 80 67877200
<b>UK</b>	<b>Tel:</b> +44 149 4 459200	<b>China</b>	<b>Tel:</b> +86 10 8219 7688
<b>France</b>	<b>Tel:</b> +33 1 41 19 21 99	<b>Singapore</b>	<b>Tel:</b> +65 6595 6388
<b>Germany</b>	<b>Tel:</b> +49 241 523030		



Distributed by:  
 Kenelec Scientific Pty Ltd  
 1300 73 22 33  
[sales@kenelec.com.au](mailto:sales@kenelec.com.au)  
[www.kenelec.com.au](http://www.kenelec.com.au)